

2024/2025

CATALOG OF FILTERS
AND FILTER SETS

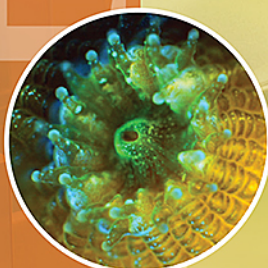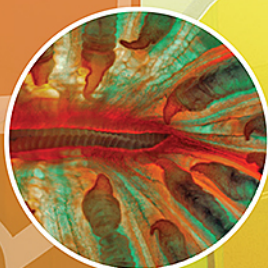Filter sets by Chroma Technology Corp

- 1 Recessed large polyps of *Goniastrea sp.*, commonly known as brain star or honeycomb coral, viewed from the side.
James H. Nicholson, NOAA.
- 2 Close-up of opening mouth of single polyped *Fungia sp.*, commonly known as mushroom coral.
James H. Nicholson, NOAA.
- 3 Recessed large polyps of *Goniastrea sp.*, commonly known as brain star or honeycomb coral, viewed from the side.
James H. Nicholson, NOAA.
- 4 Close-up of a single *Montastraea sp.* (commonly known as boulder star coral) polyp within a colony.
James H. Nicholson, NOAA.

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














Single Bandpass and Single Edge Filters

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





































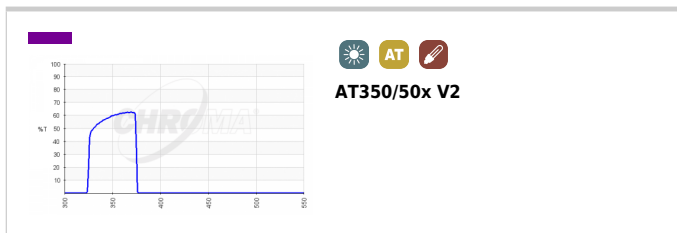
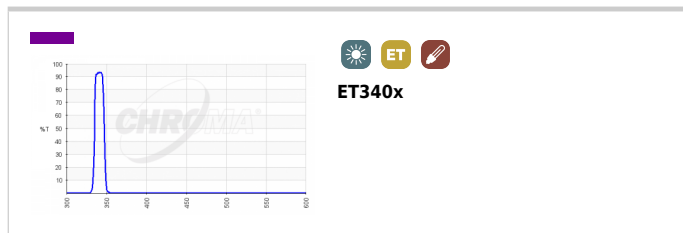
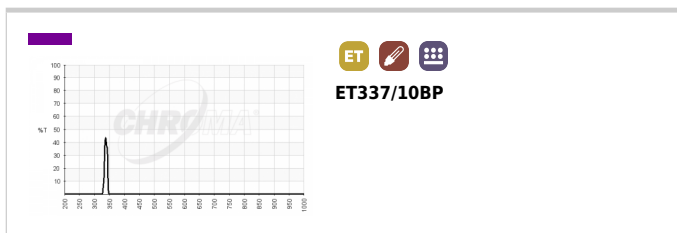
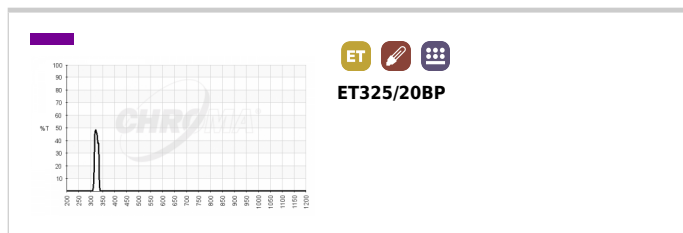
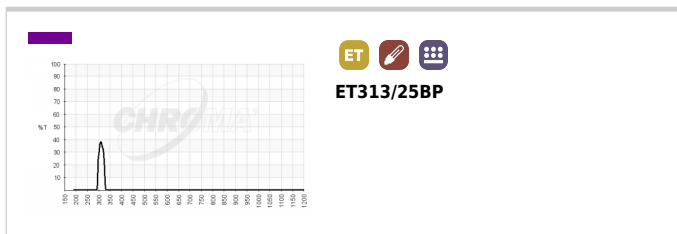
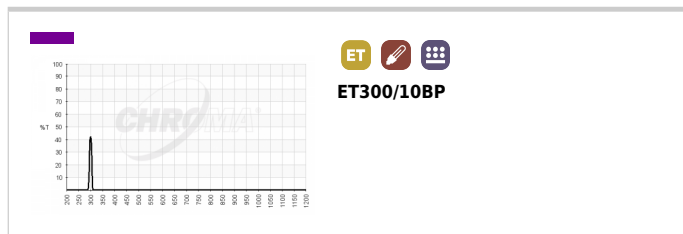
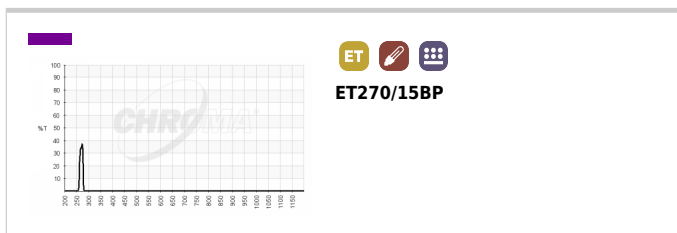
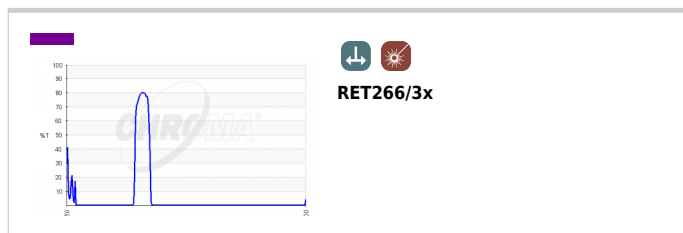
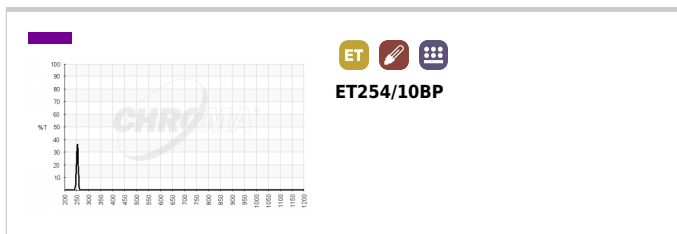
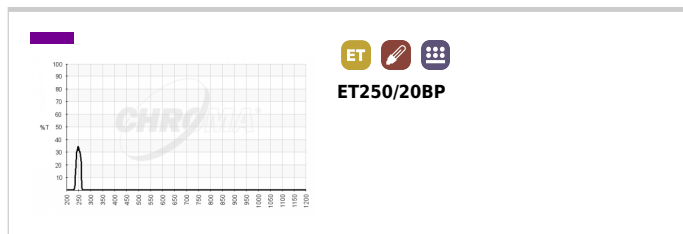
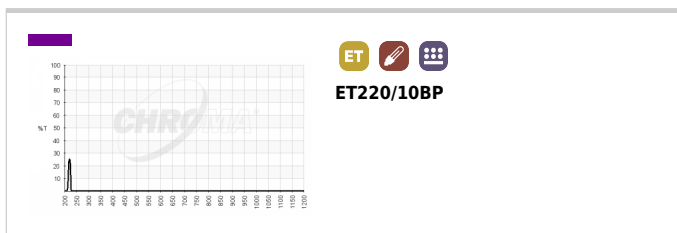
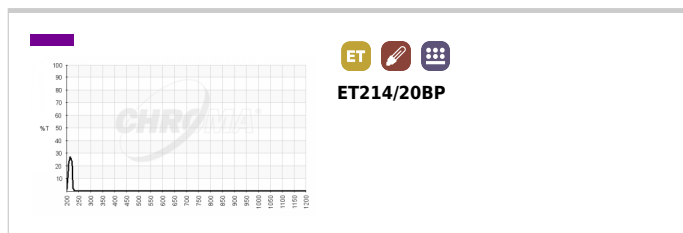
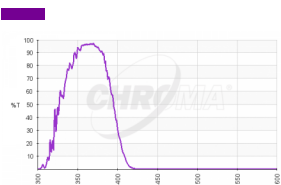
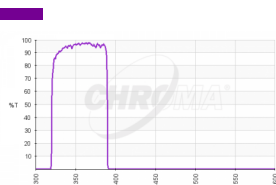
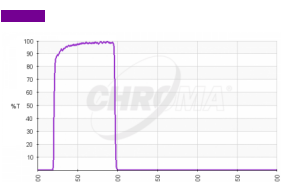
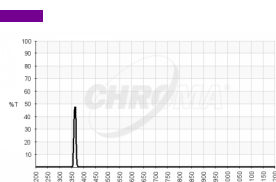
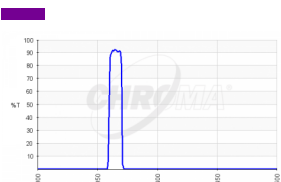
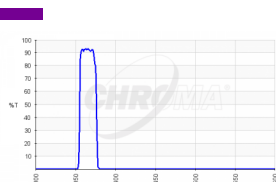
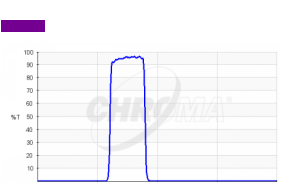
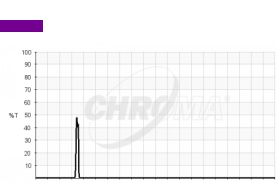
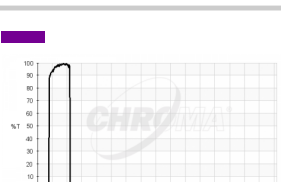
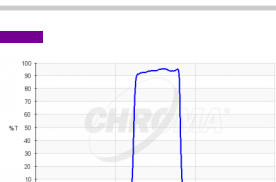
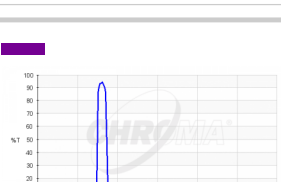
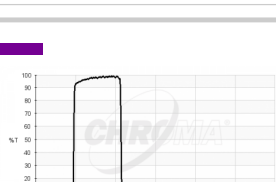
**ET365/10BP****ET365/10x****ZET365/20x****AT375/28x****ET375/10BP****MV375/70 - UV Vision Filter****ET376/30x****ET380x****MV380/60 - UV Vision Filter**

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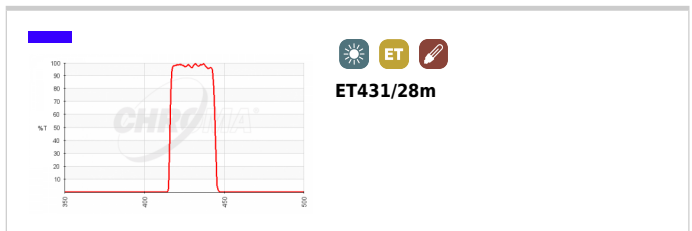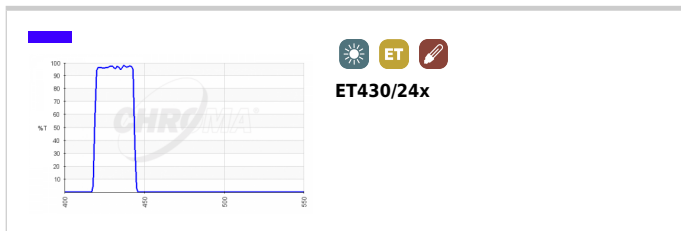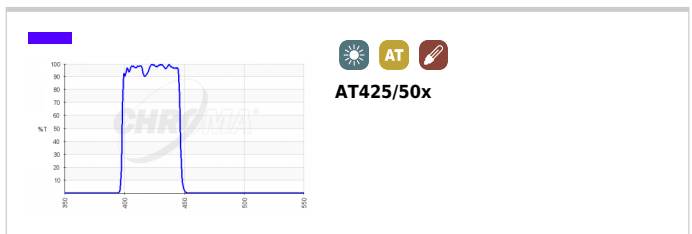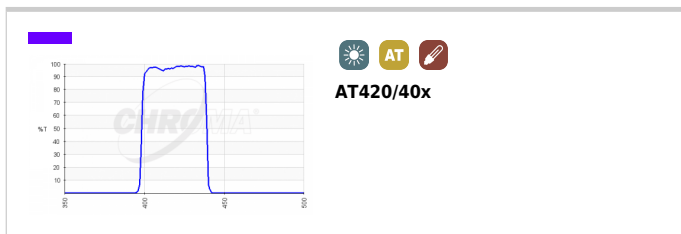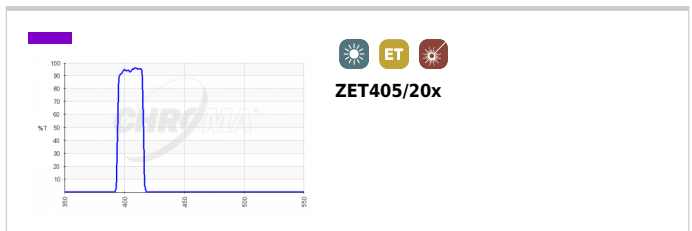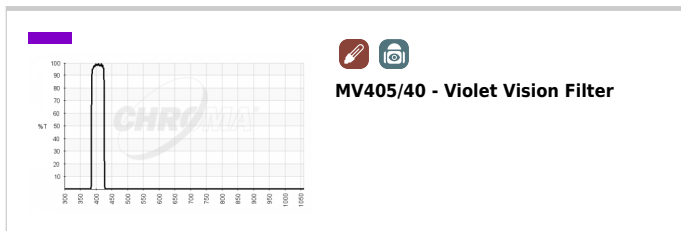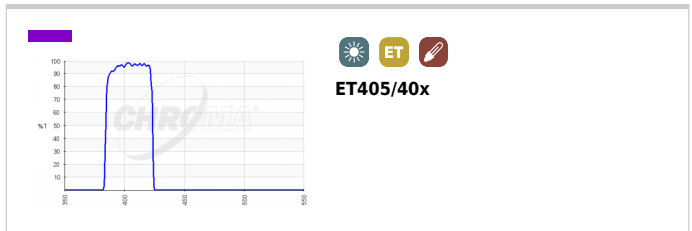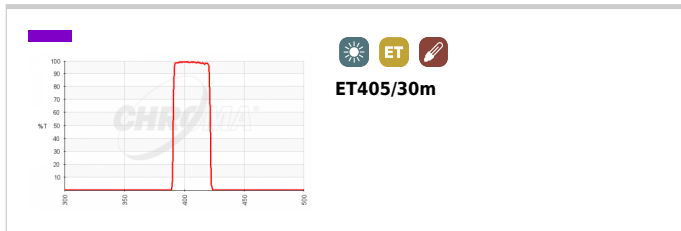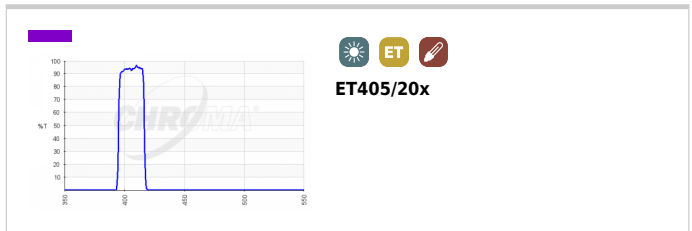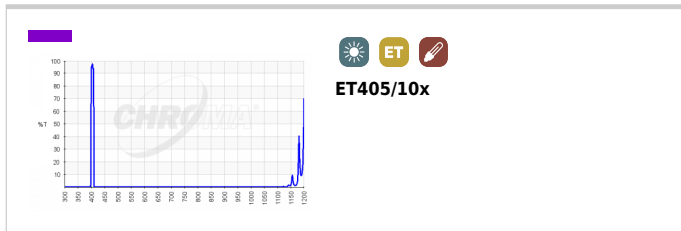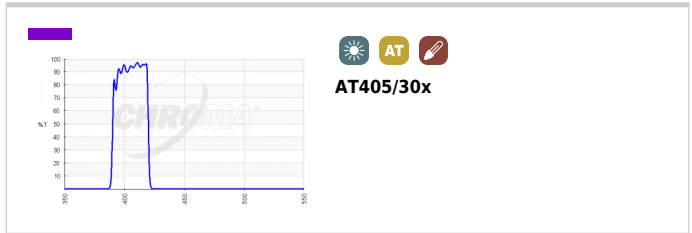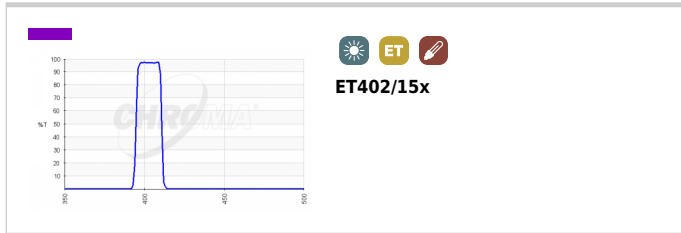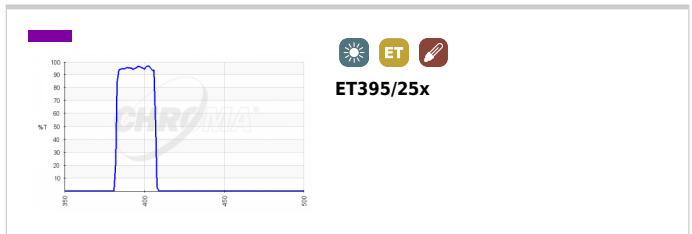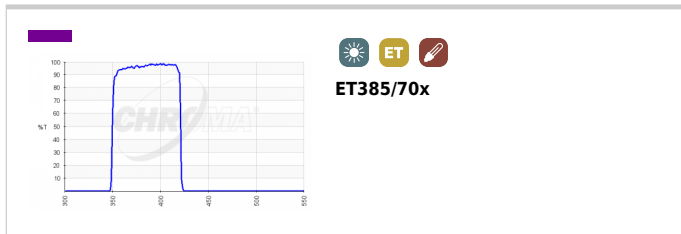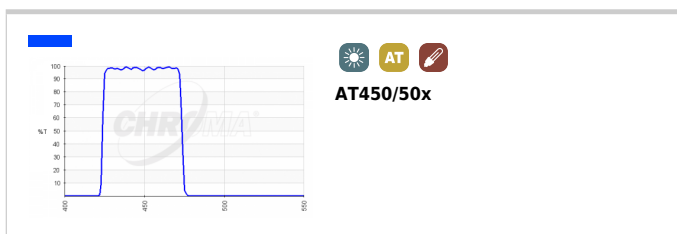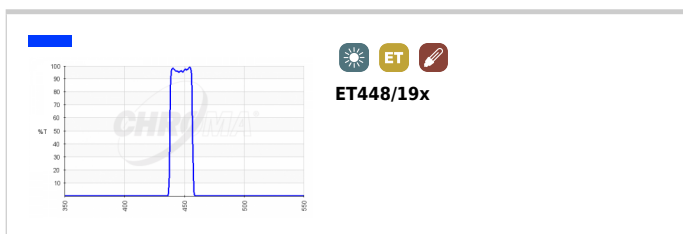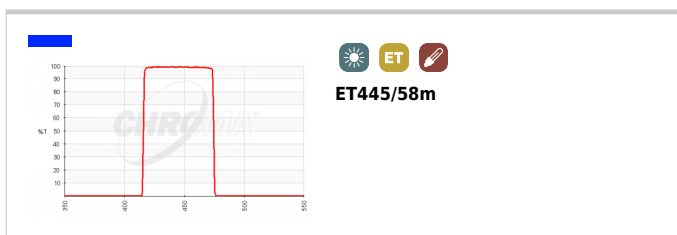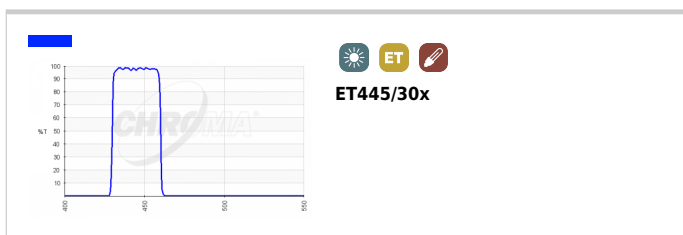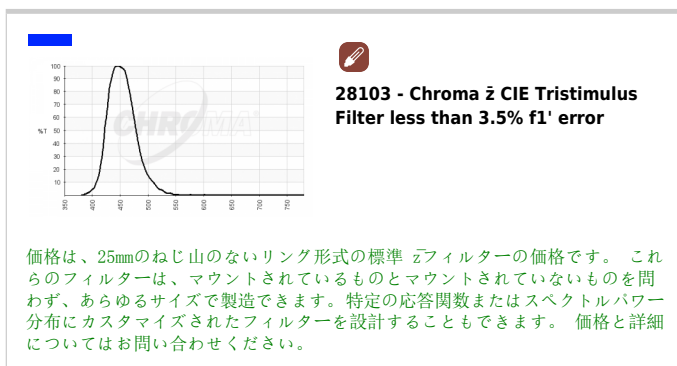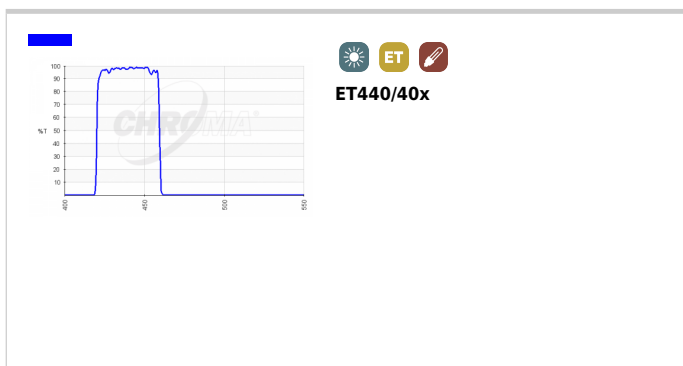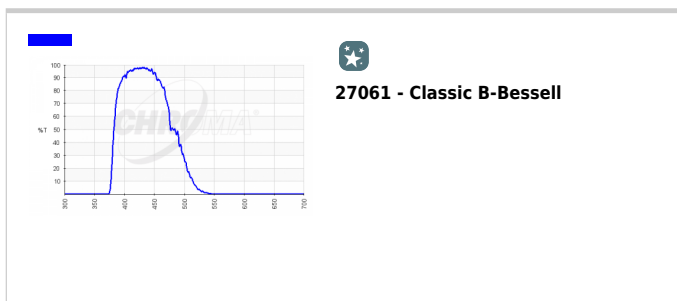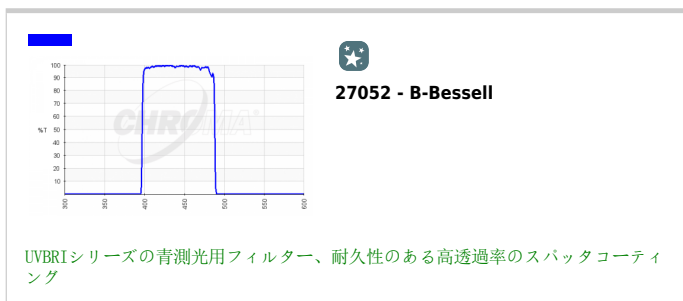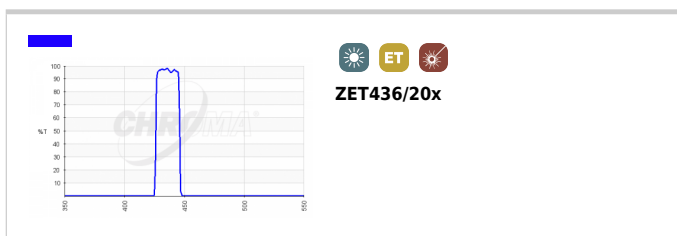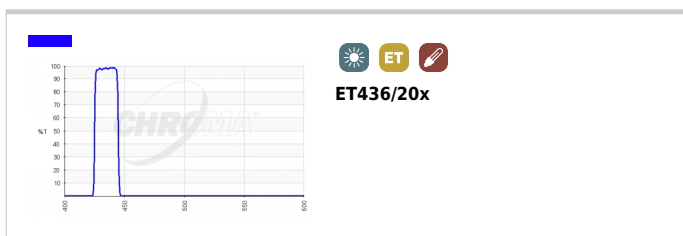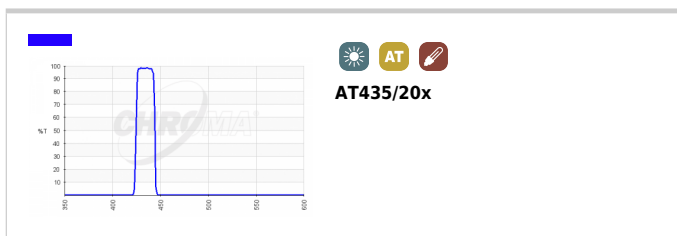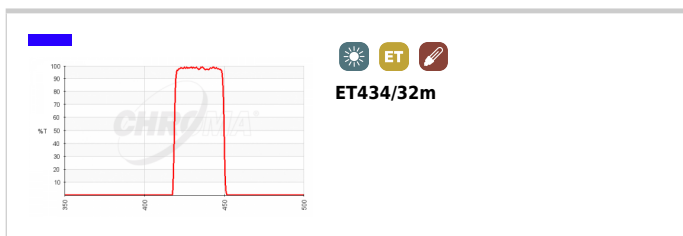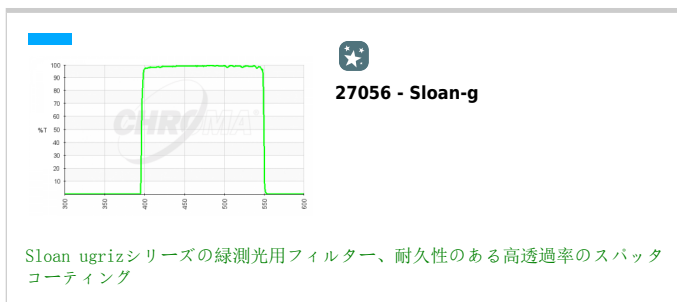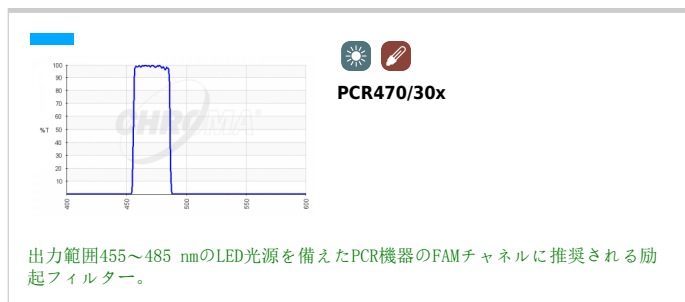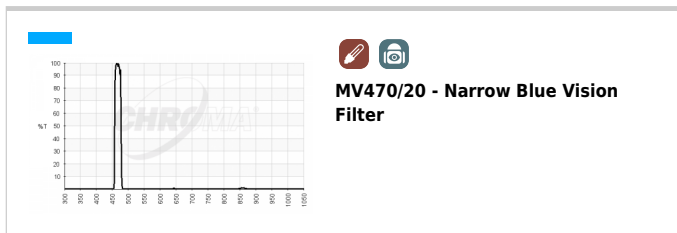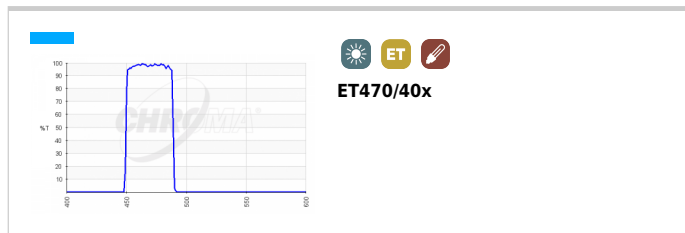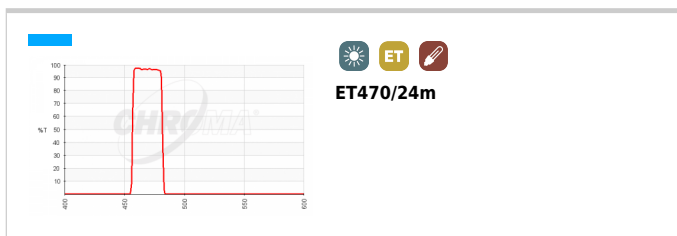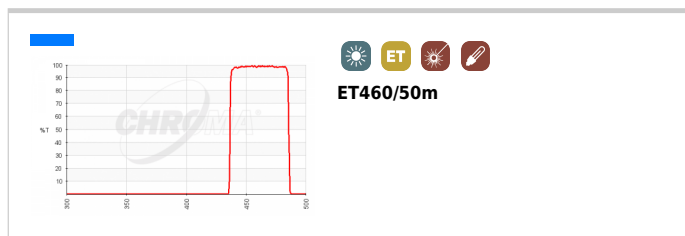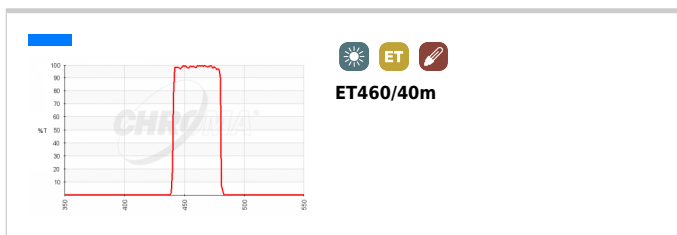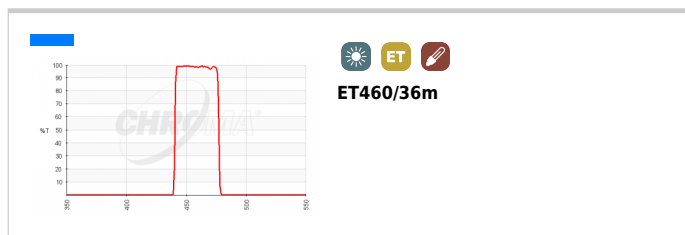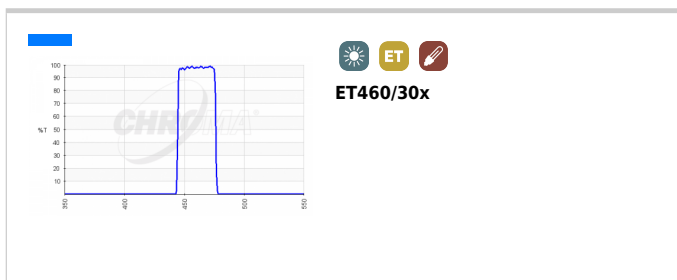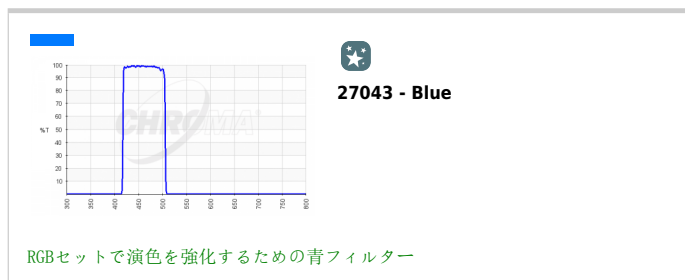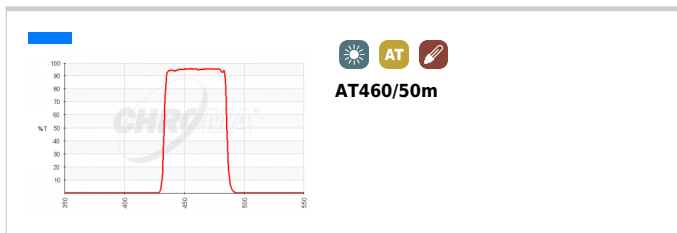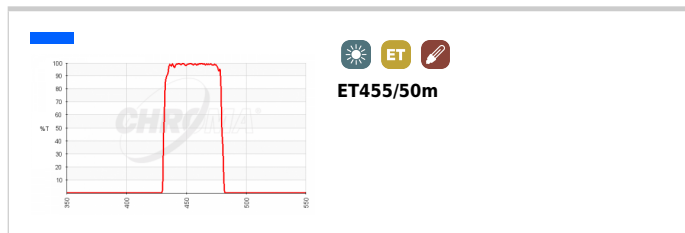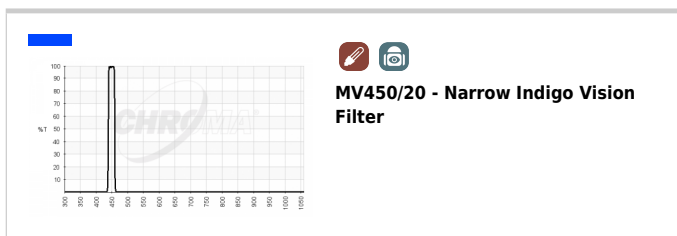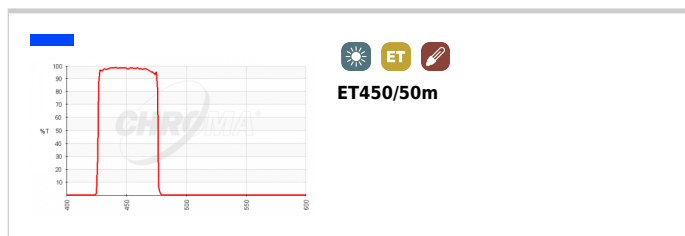
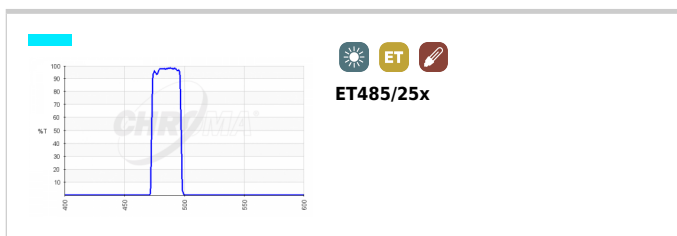
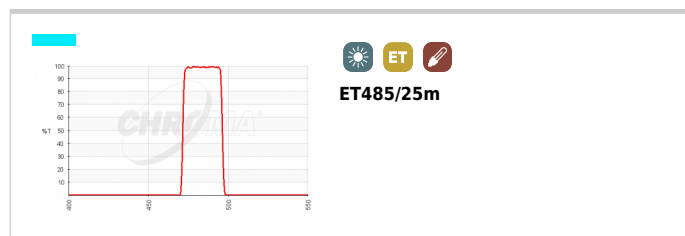
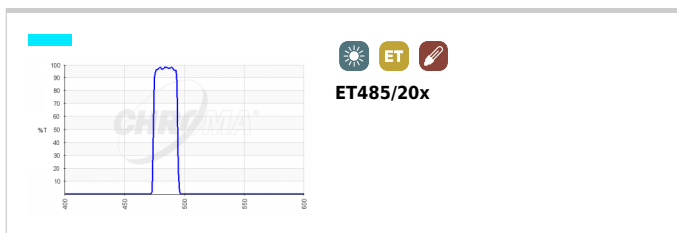
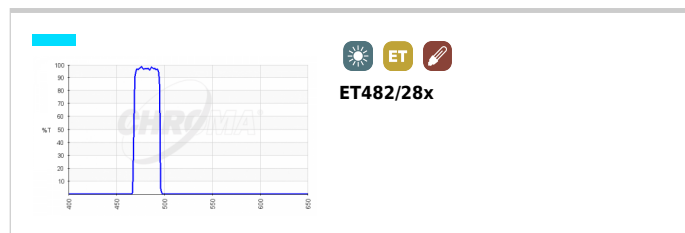
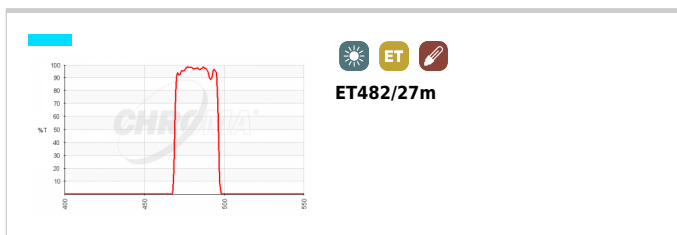
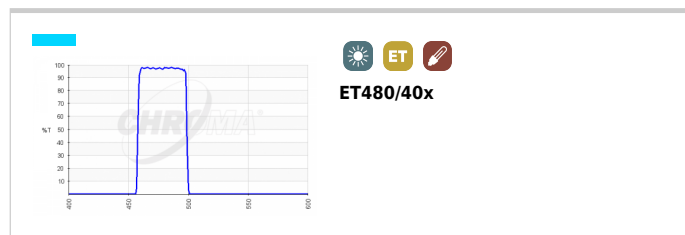
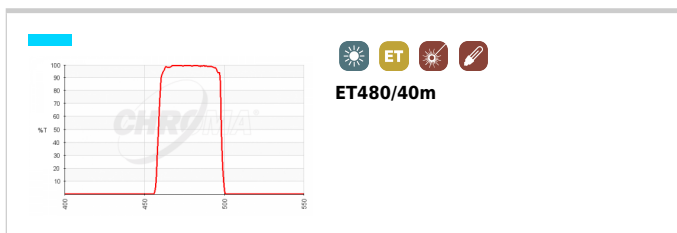
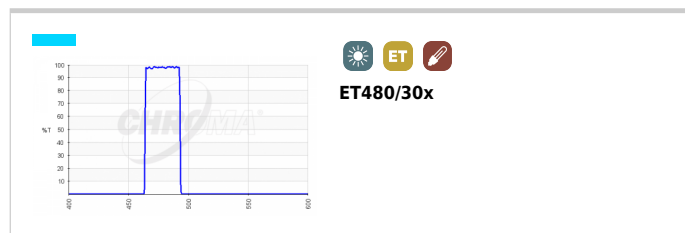
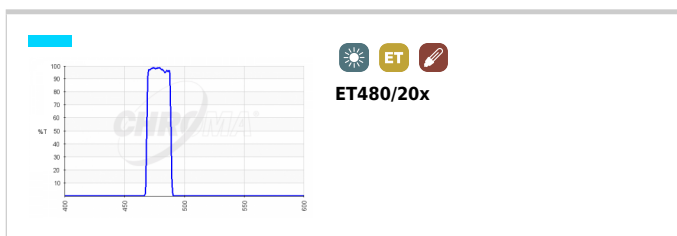
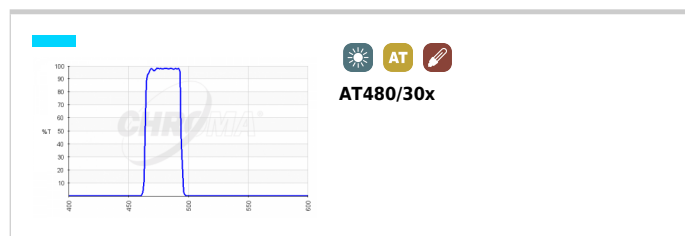
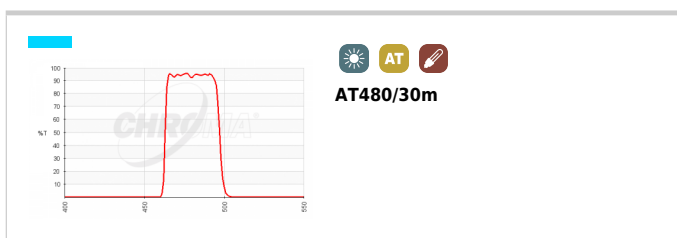
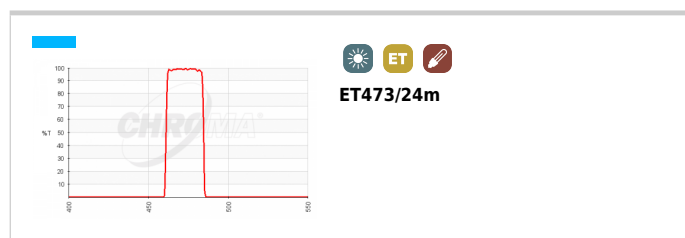
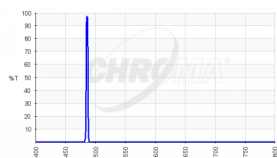
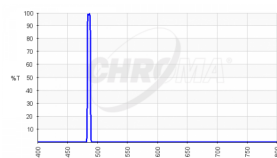
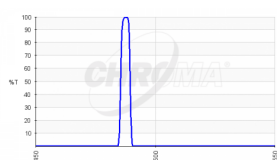
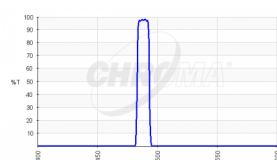
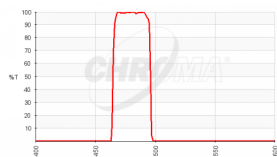
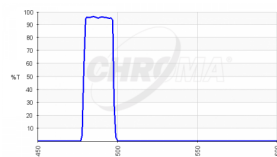
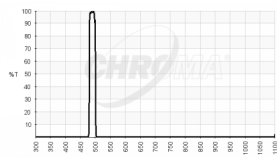
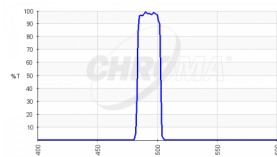
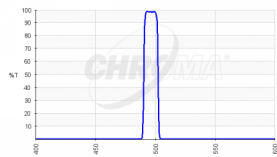
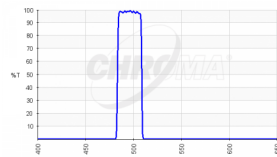
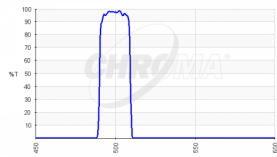
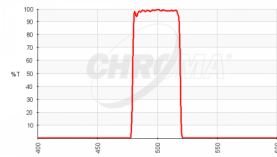
**RET488/4x****ZET488/10x****ET480/30m****ET490/20x****MV490/20 - Narrow Blue Vision Filter****AT495/20x****ET495/10x****ET495/25x****ET500/20x****ET500/40m**

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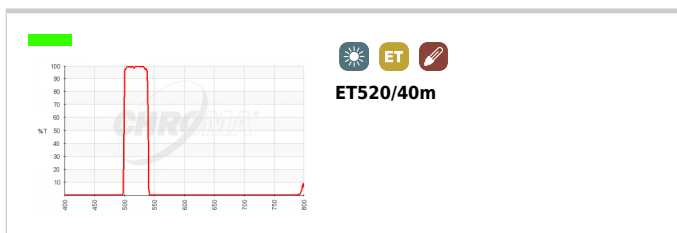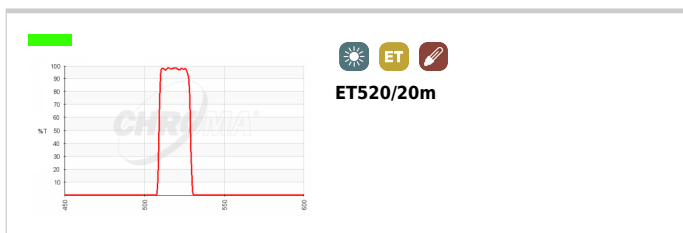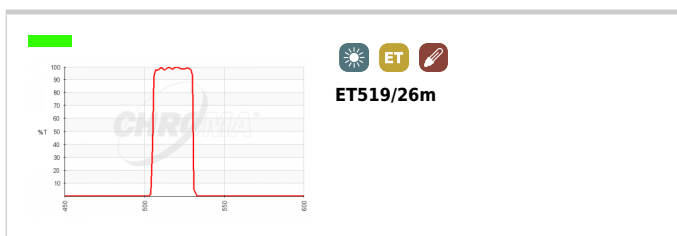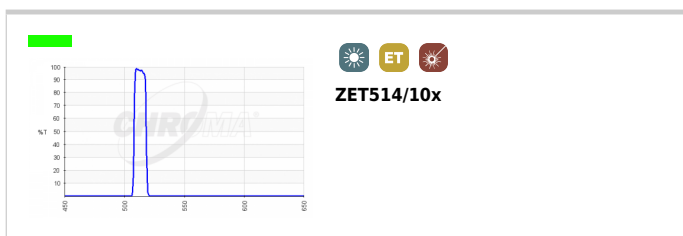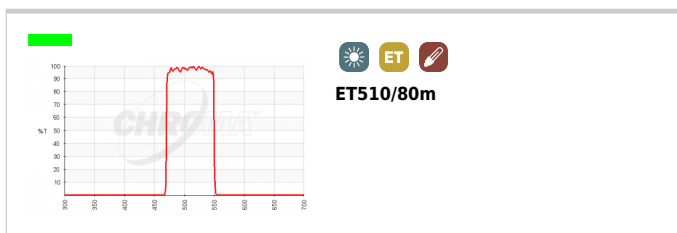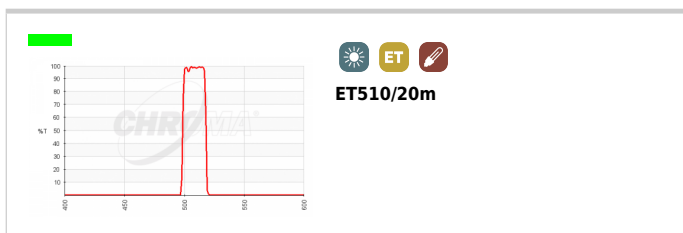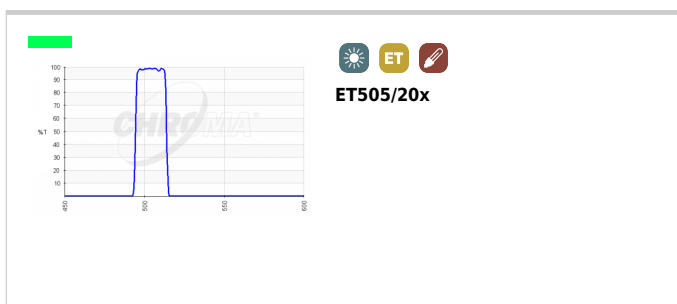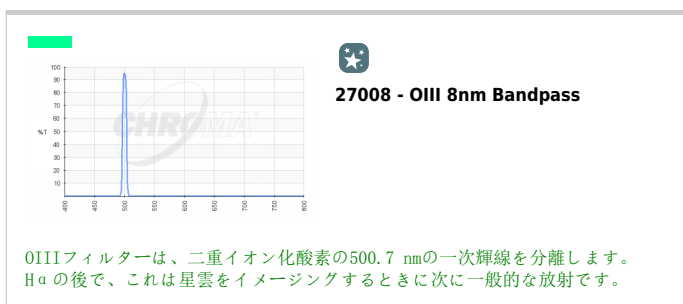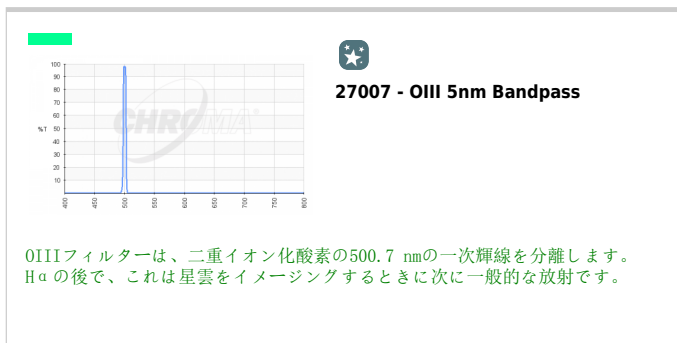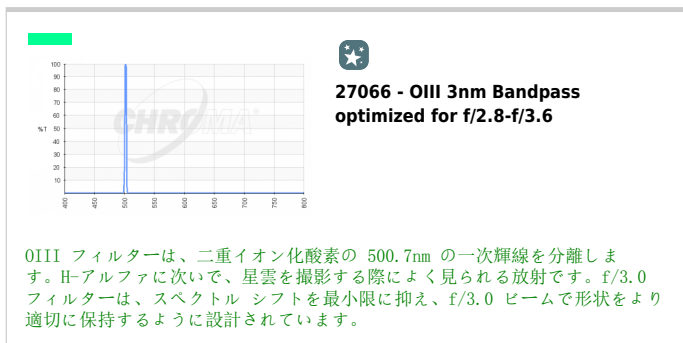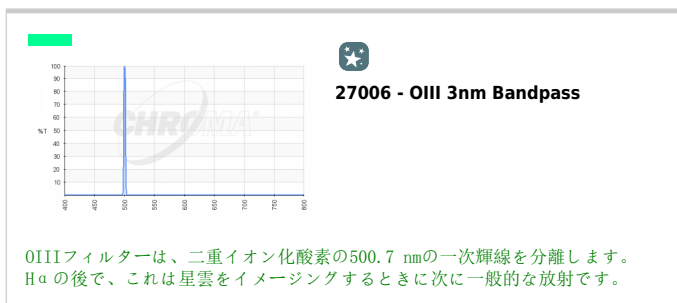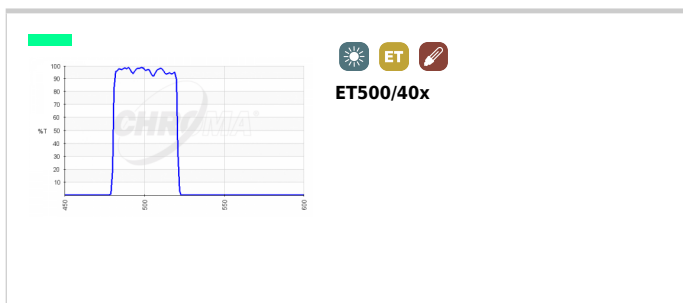
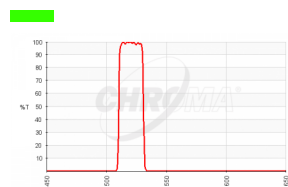
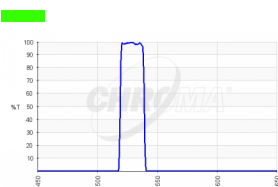
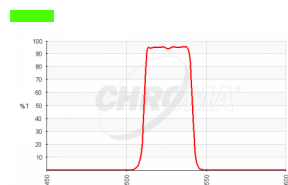
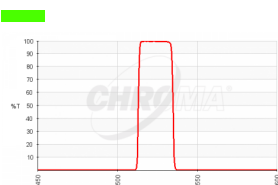
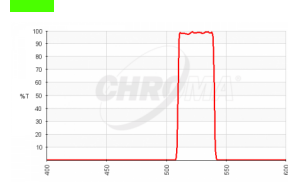
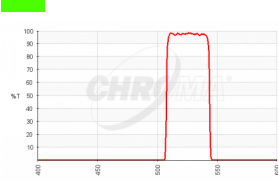
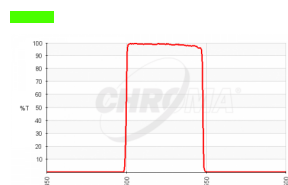
   **ET525/50m**

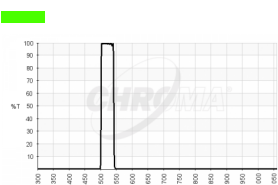
  **MV525/40 - Green Vision Filter**

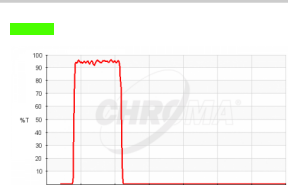
   **ET525/70m-2p - Multiphoton bandpass emission filter**

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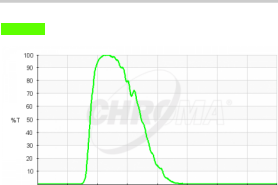
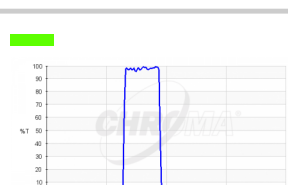
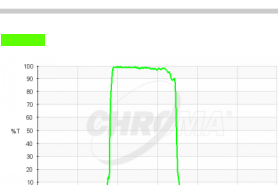
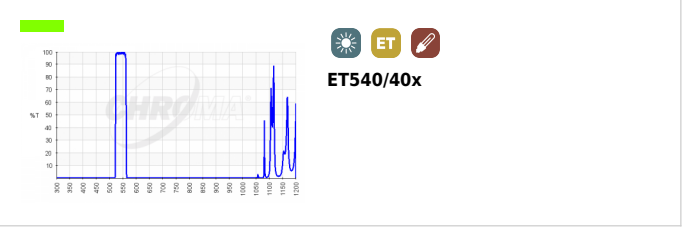
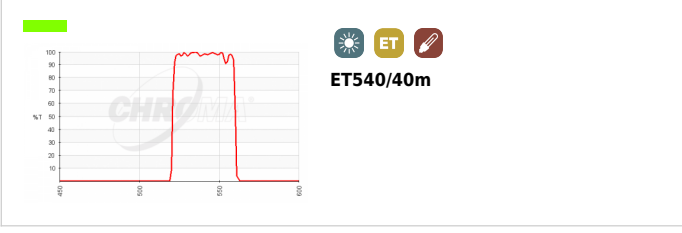
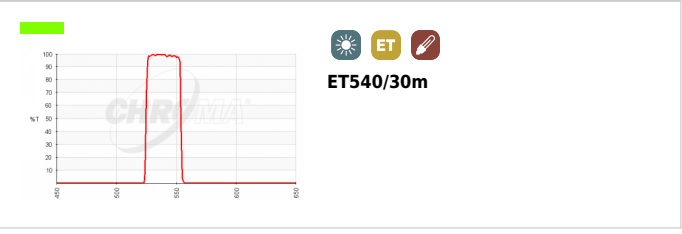
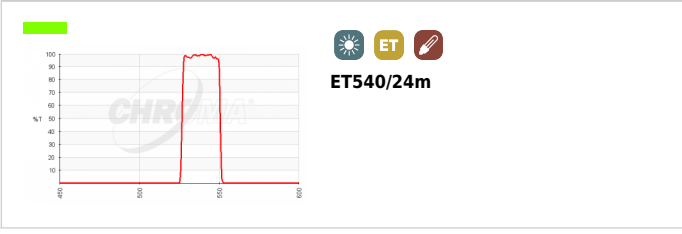
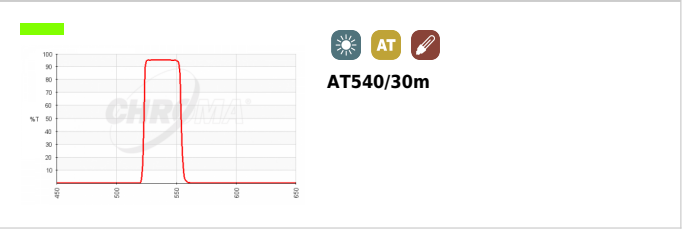
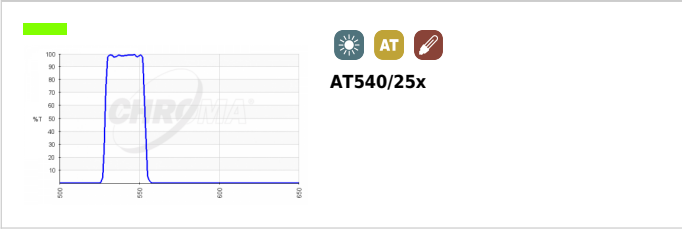
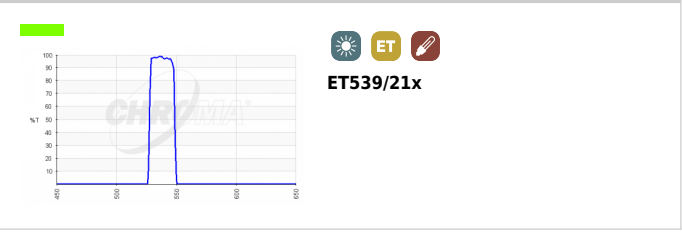
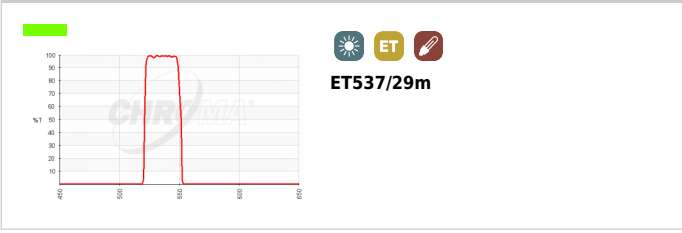
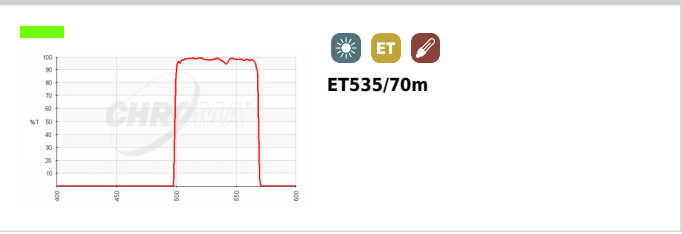
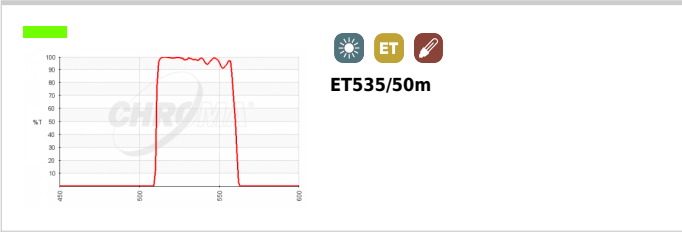
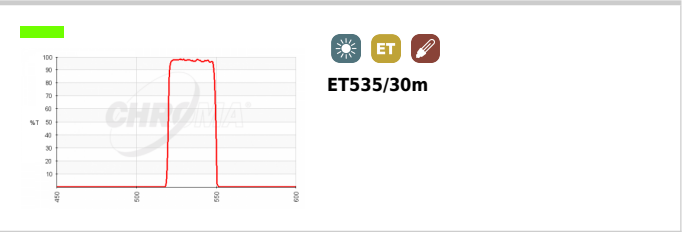
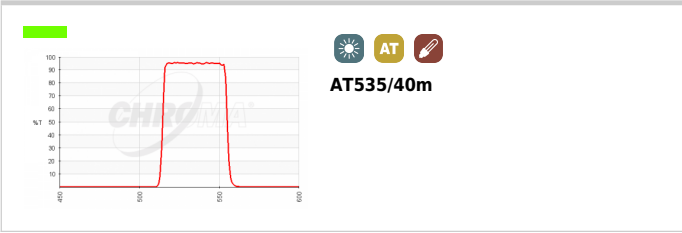
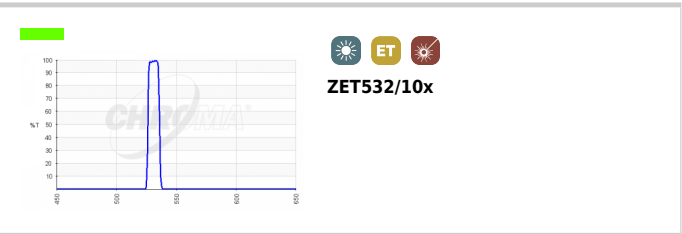
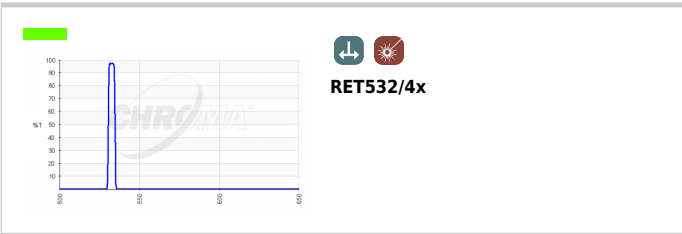
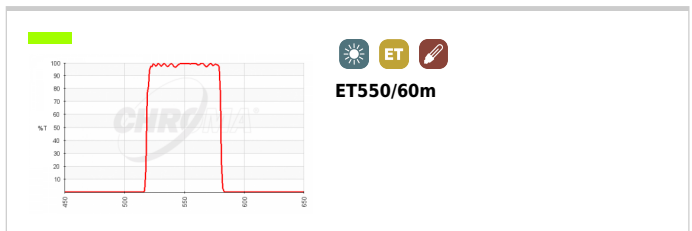
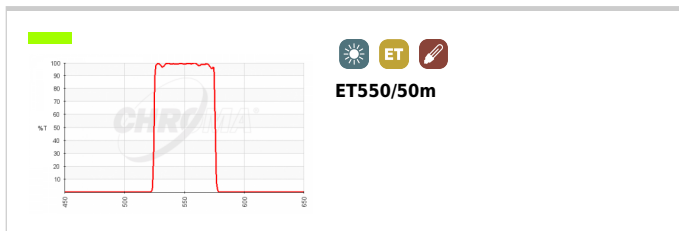
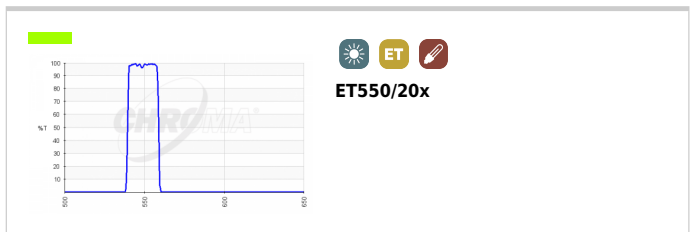
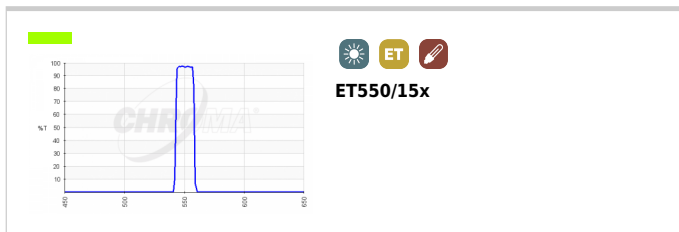
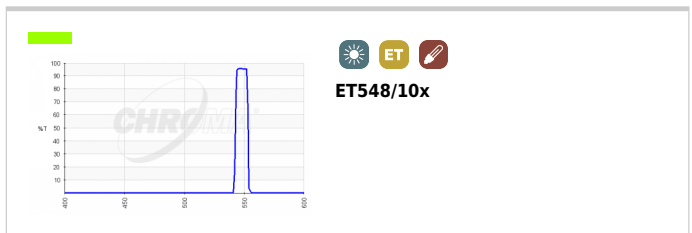
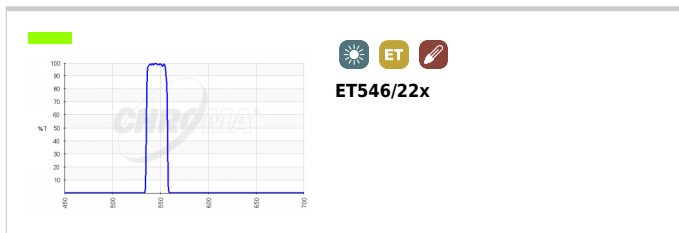
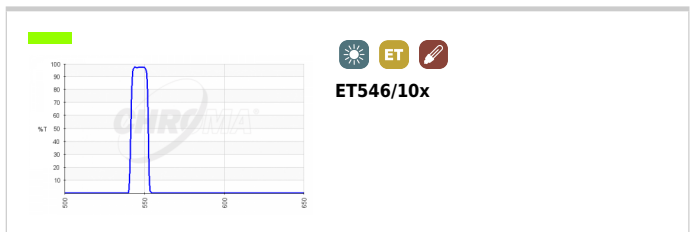
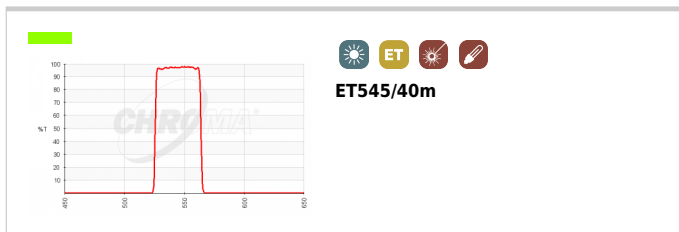
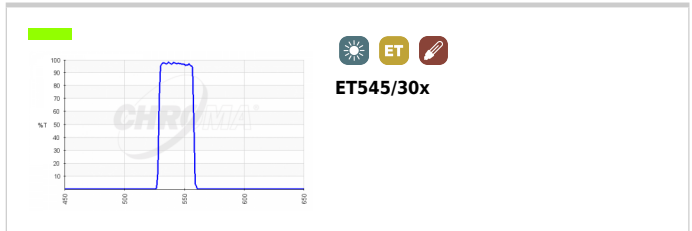
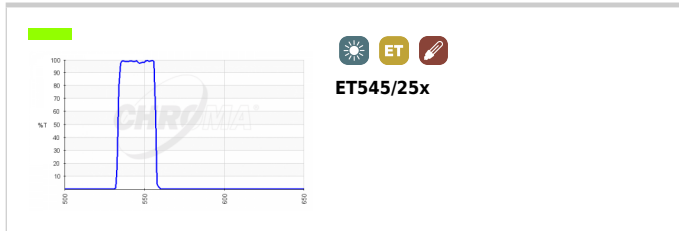
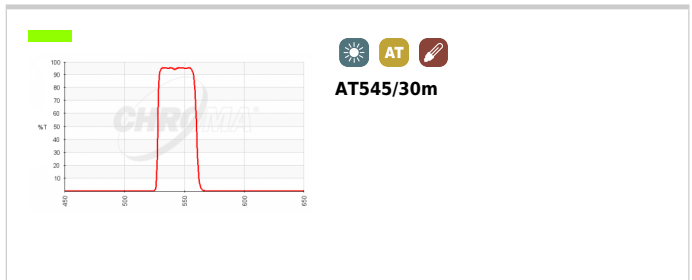
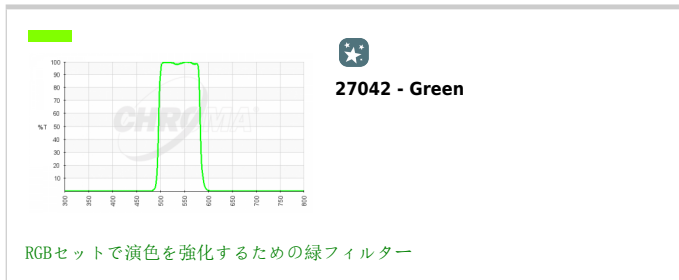
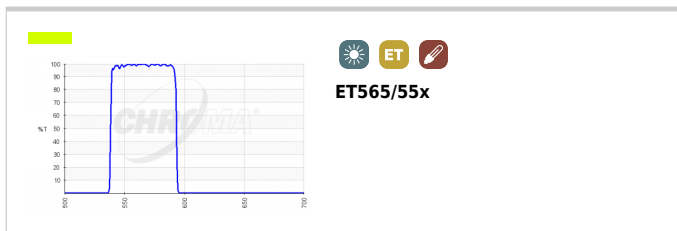
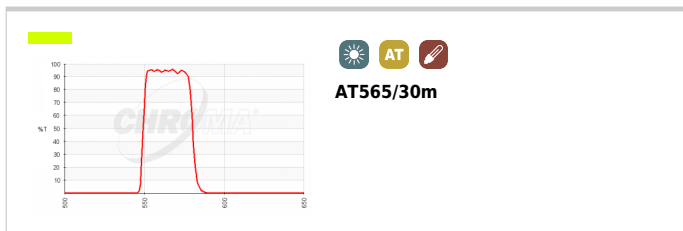
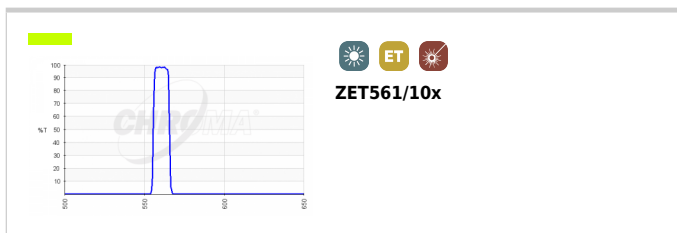
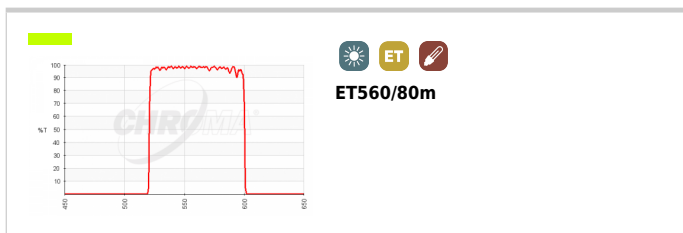
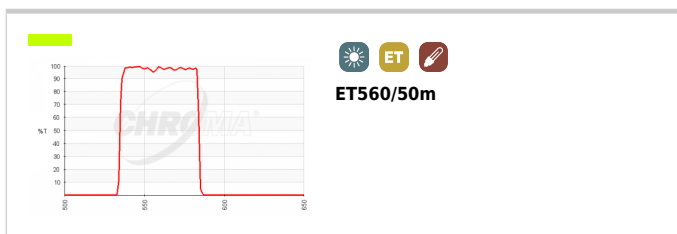
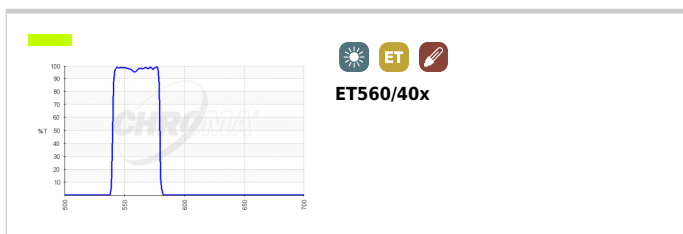
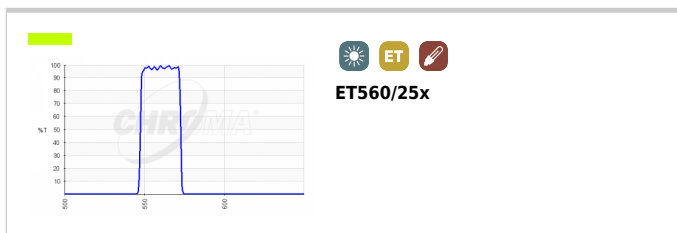
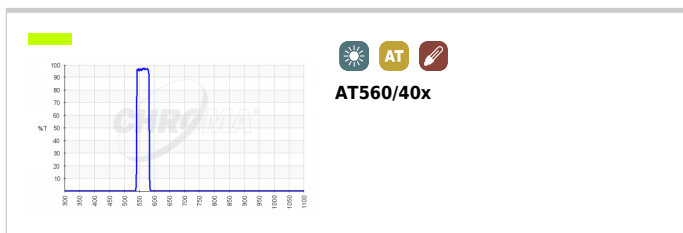
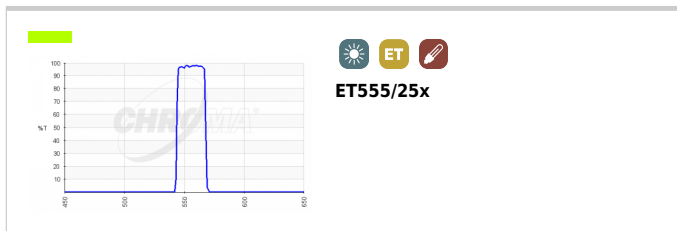
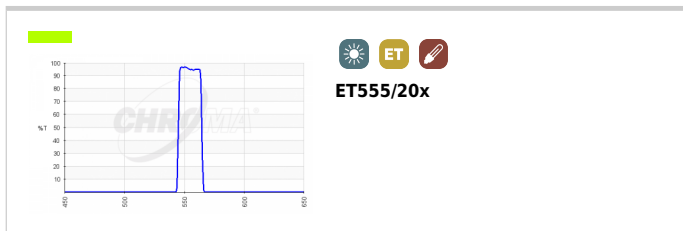
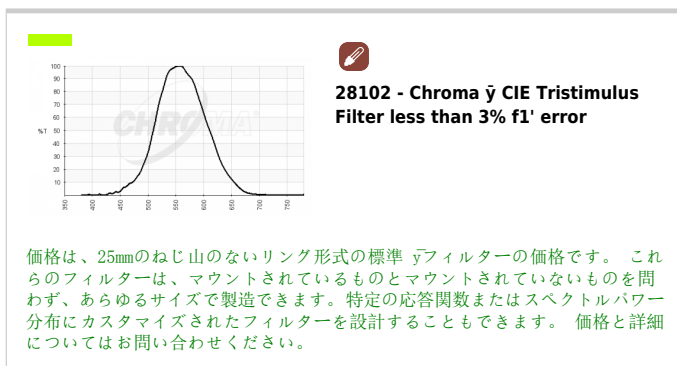
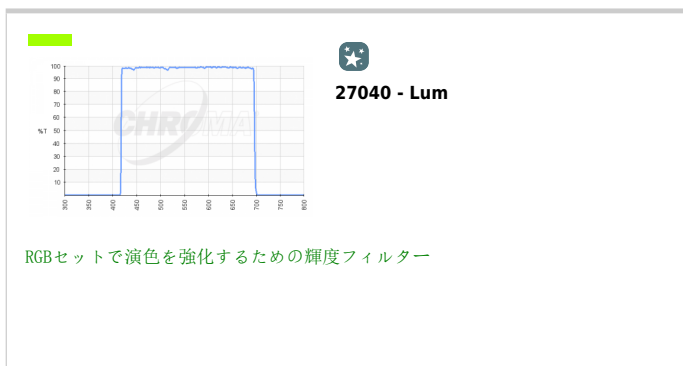
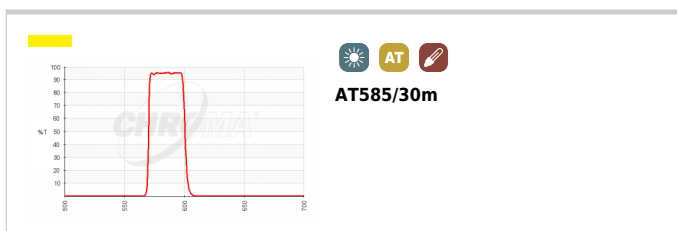
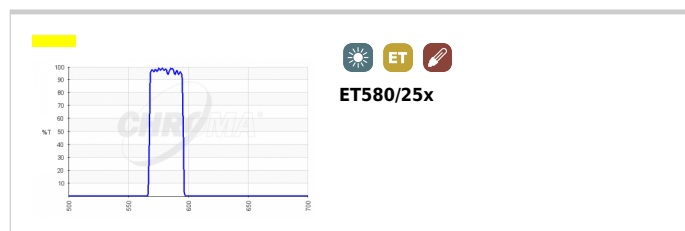
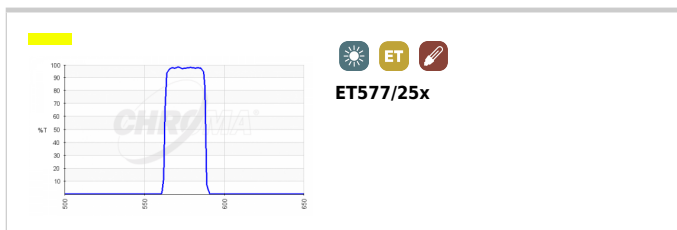
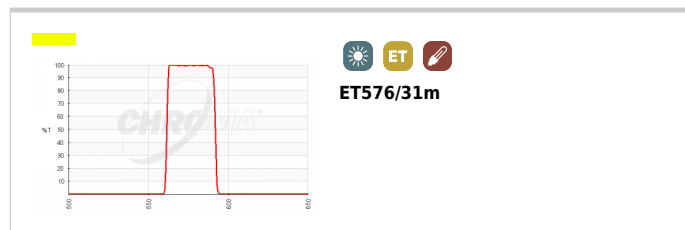
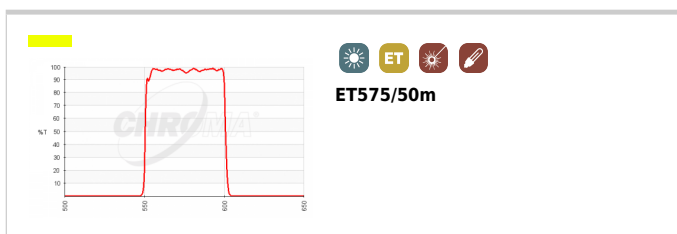
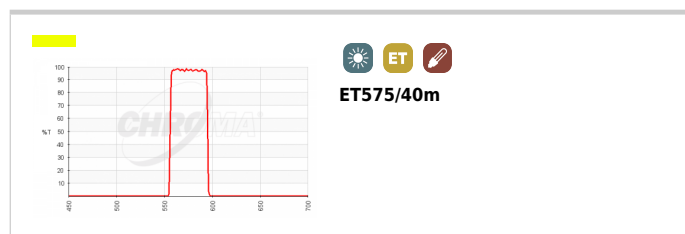
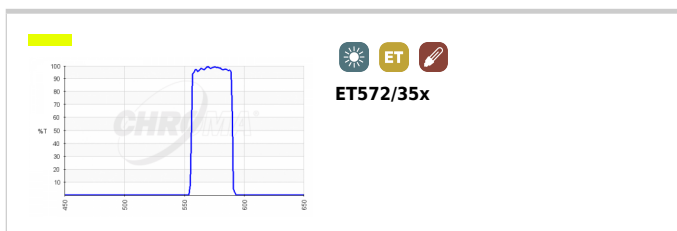
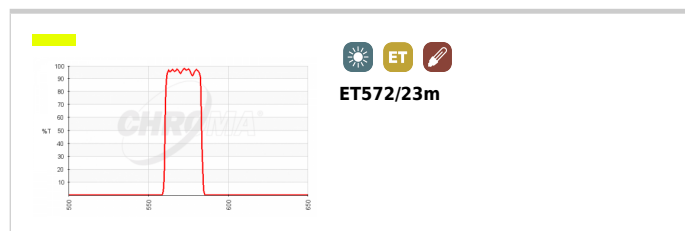
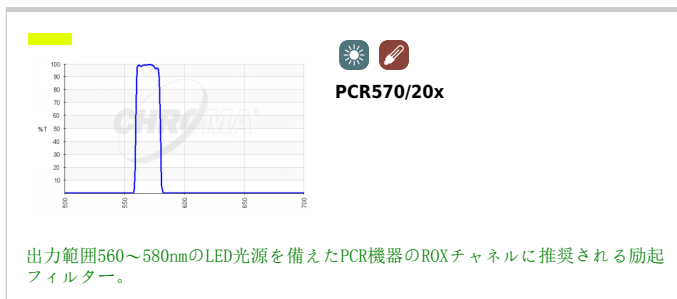
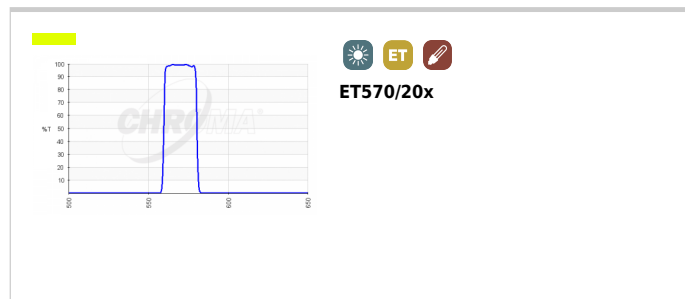
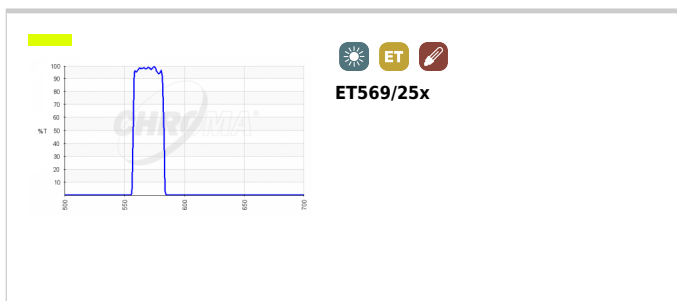
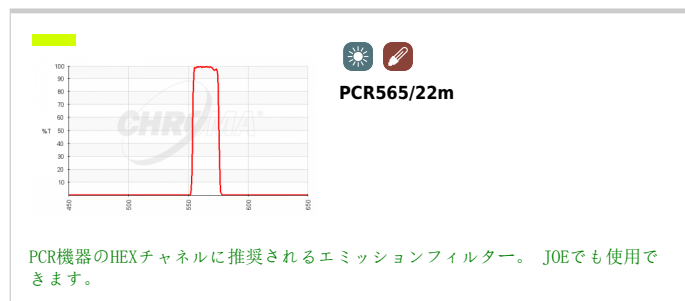
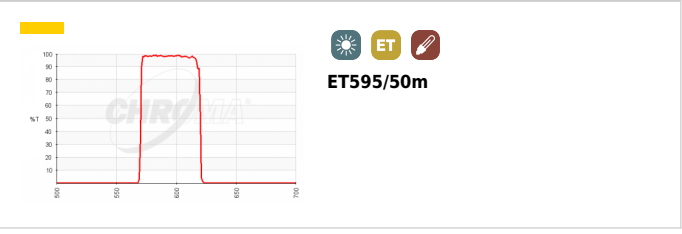
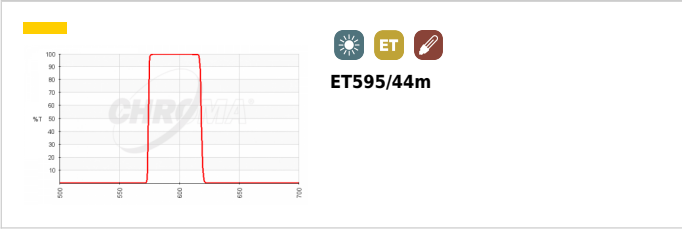
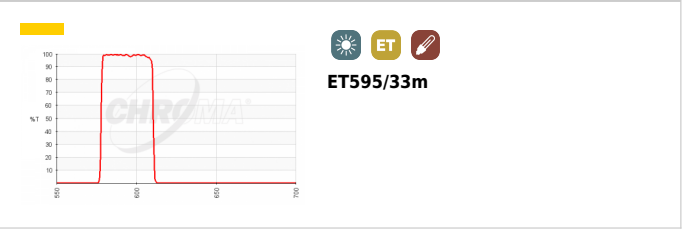
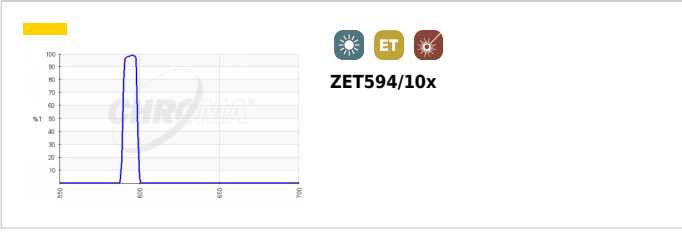
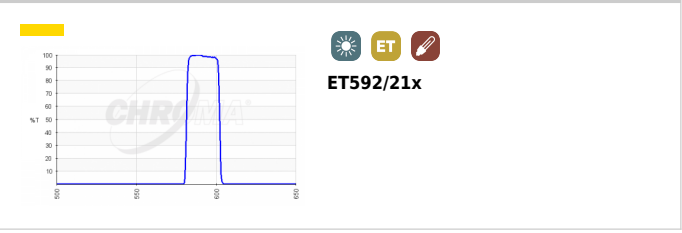
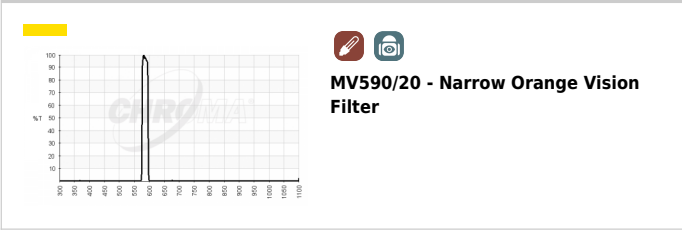
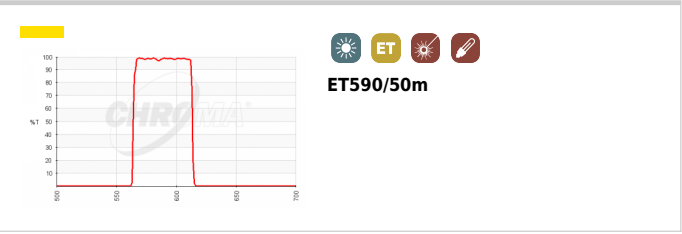
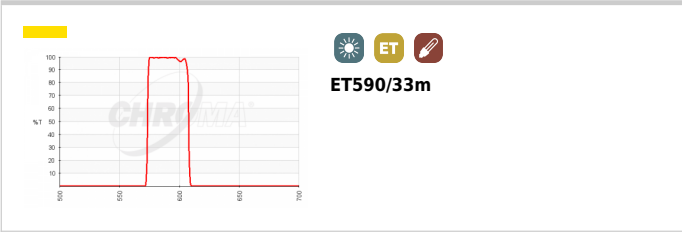
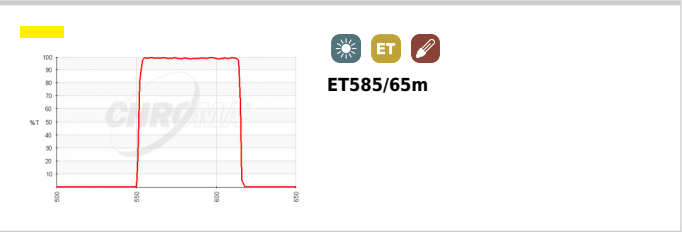
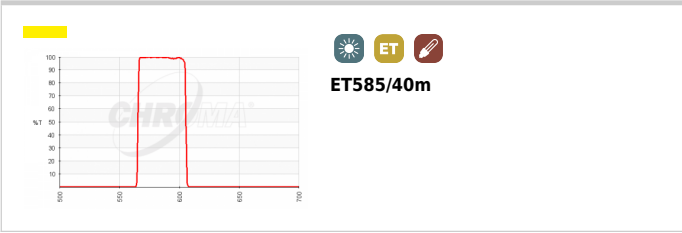
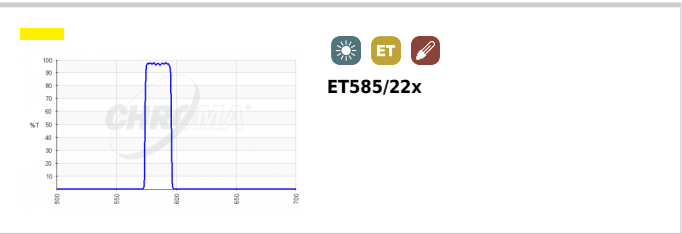
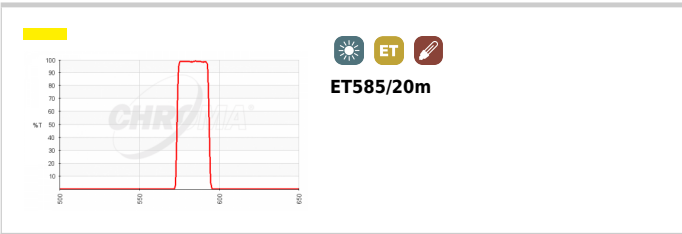
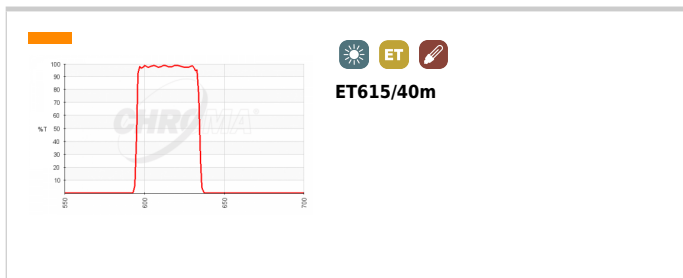
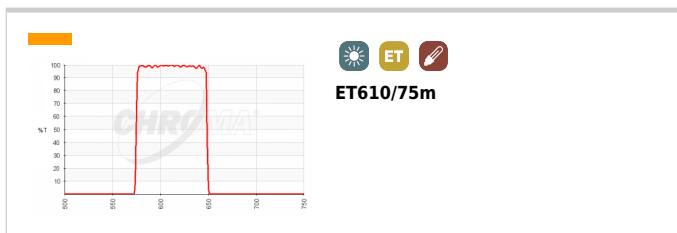
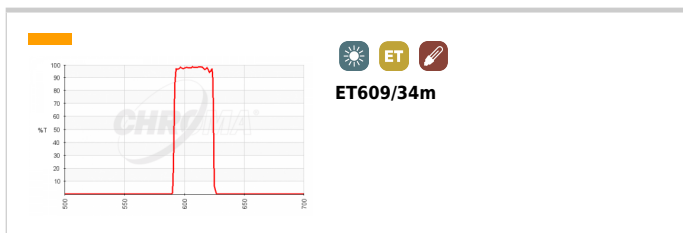
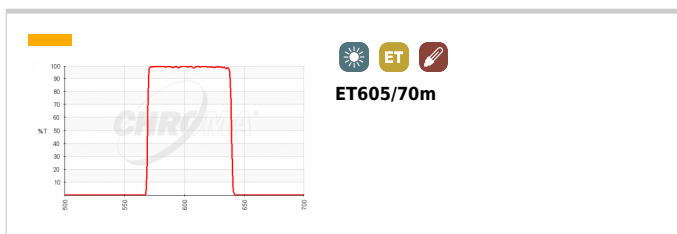
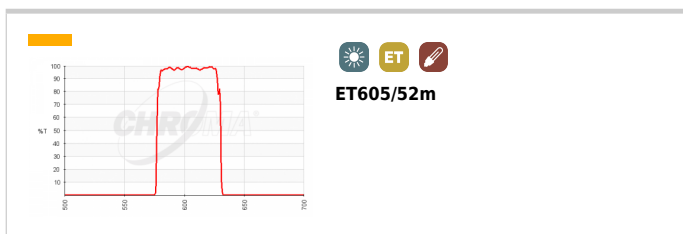
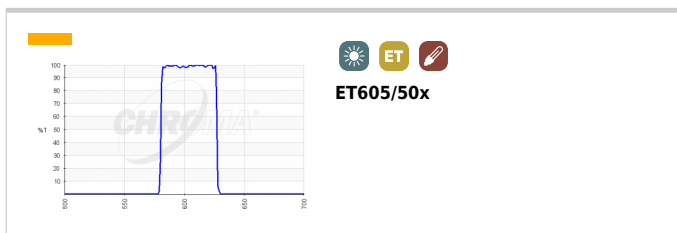
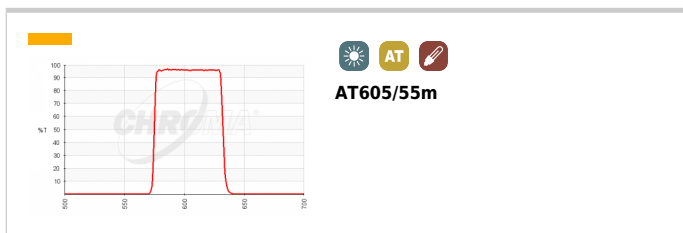
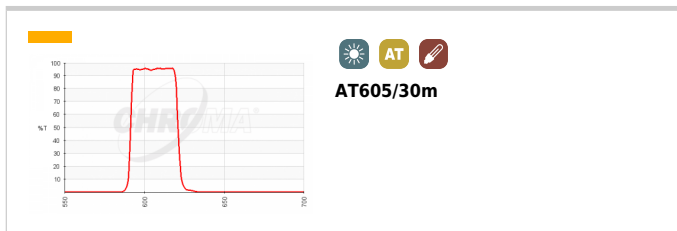
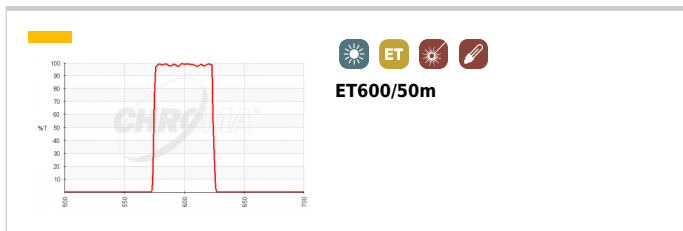
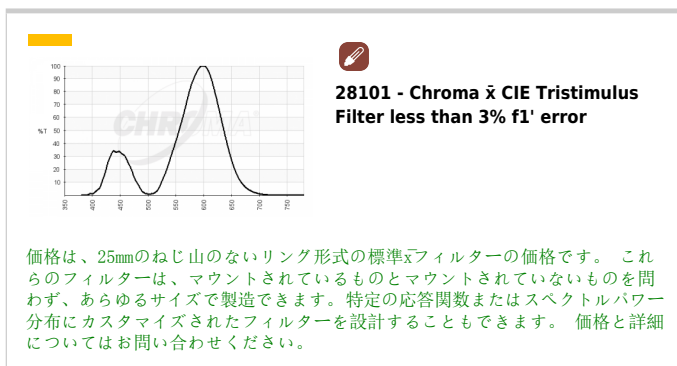
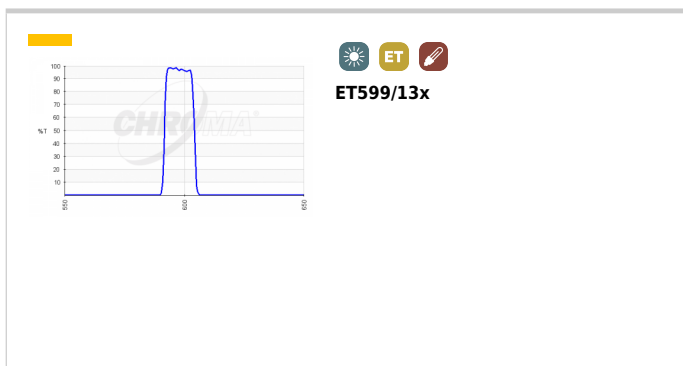
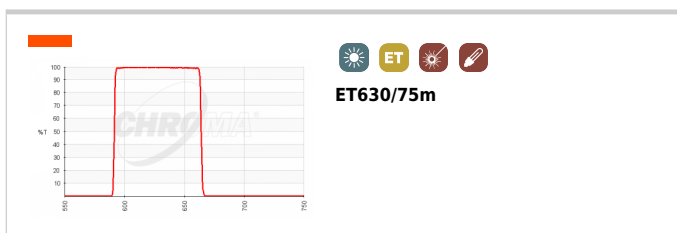
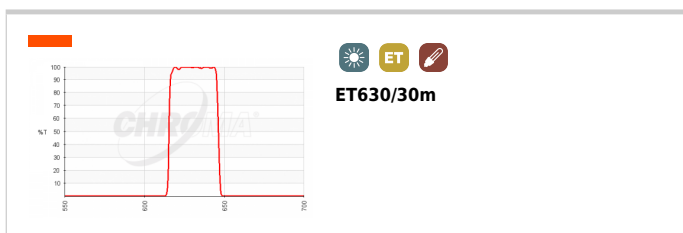
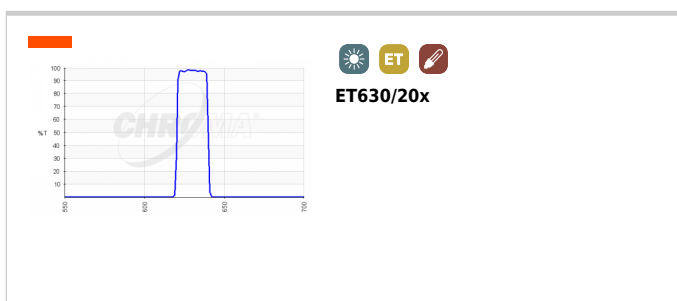
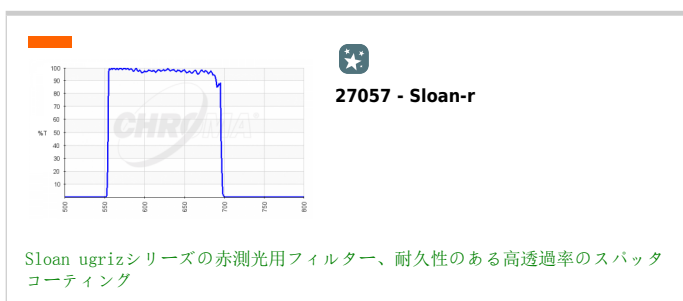
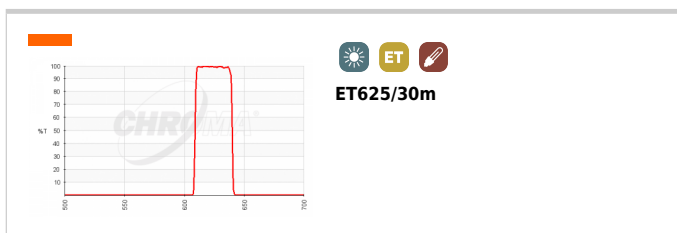
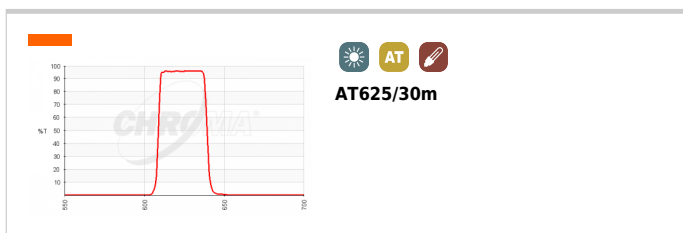
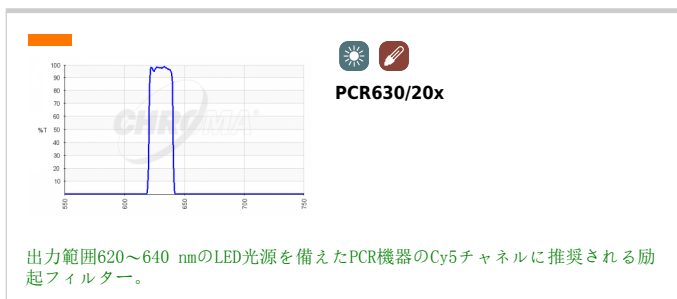
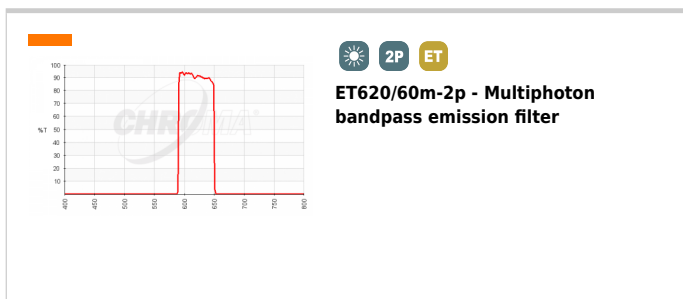
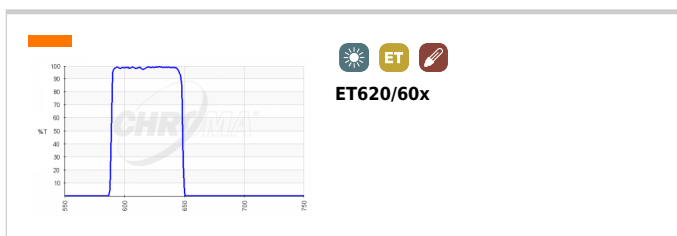
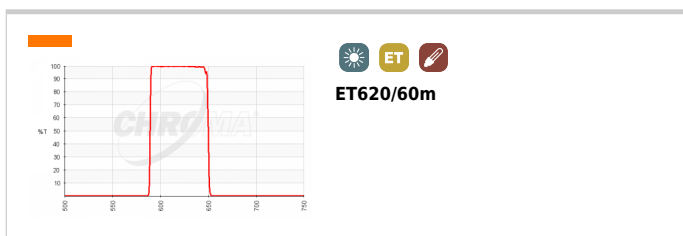
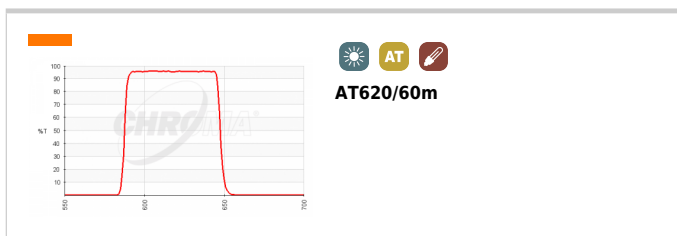
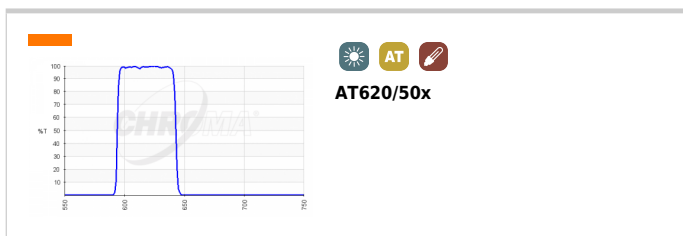
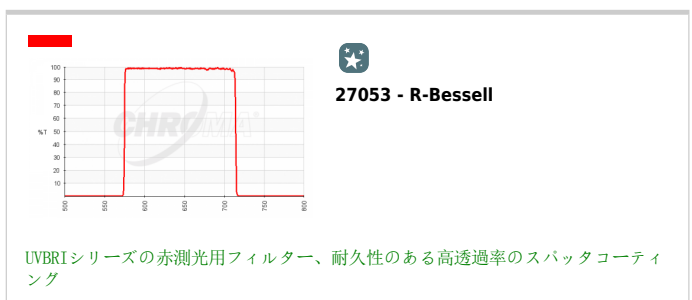
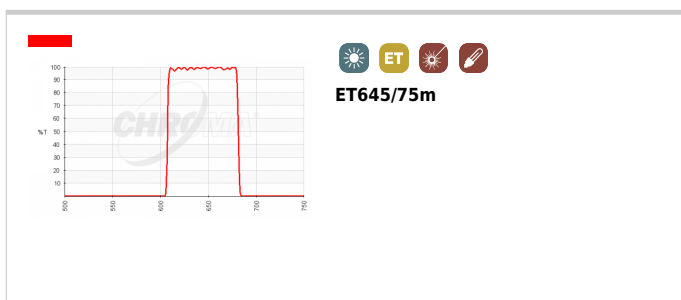
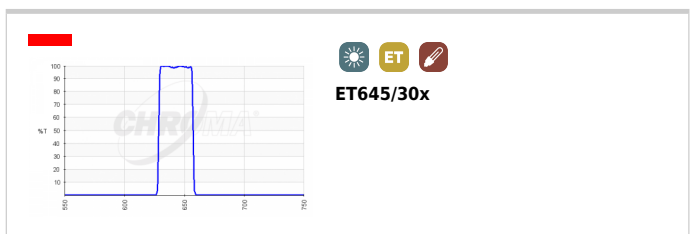
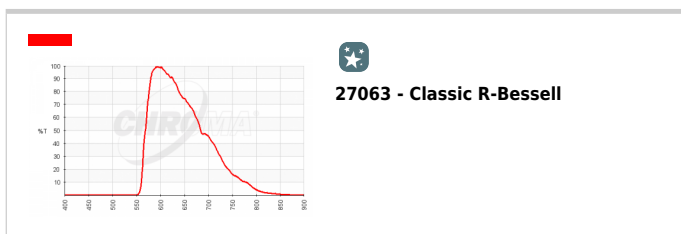
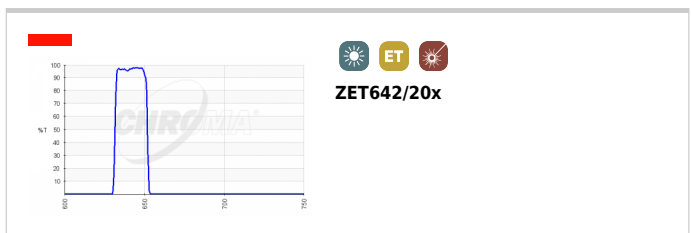
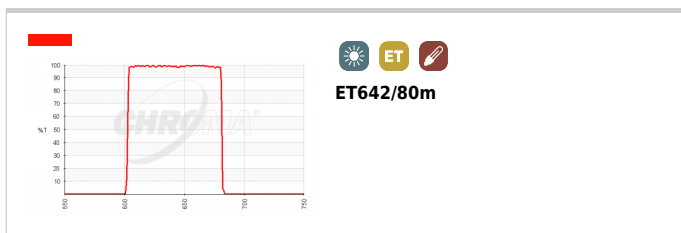
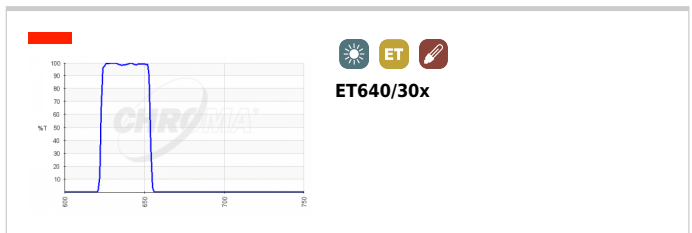
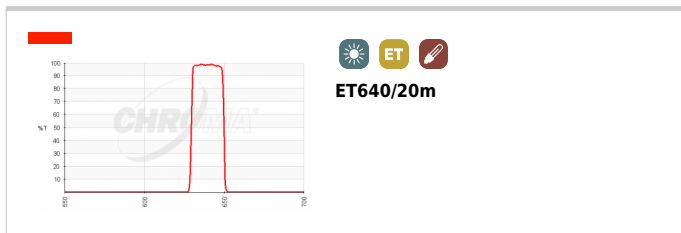
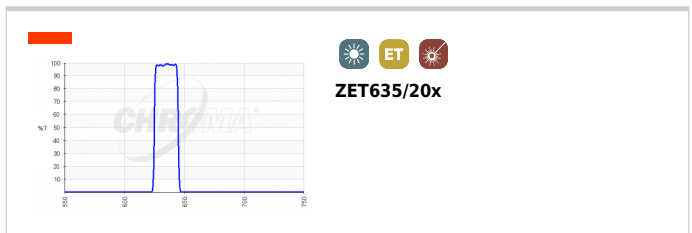
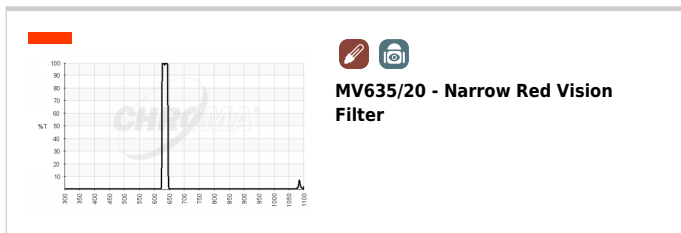
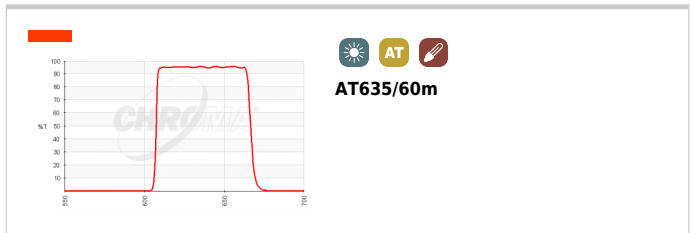
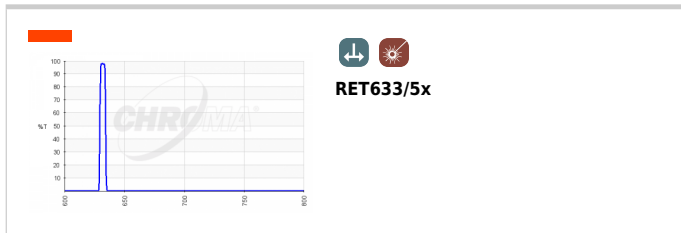
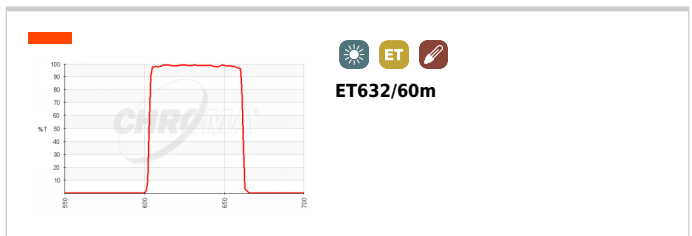
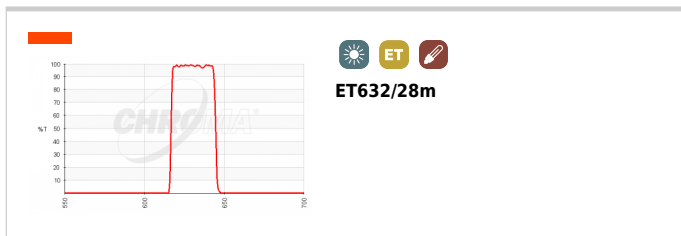
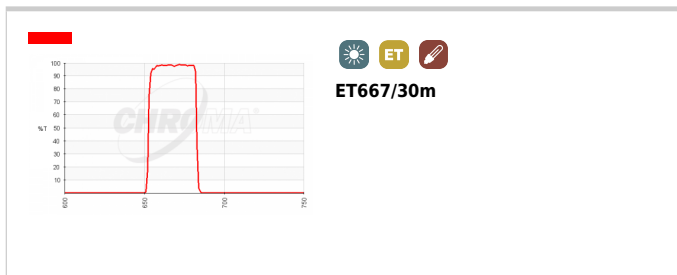
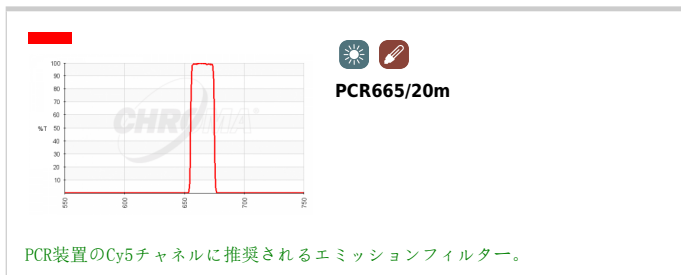
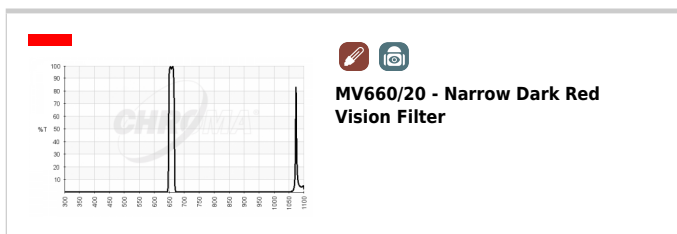
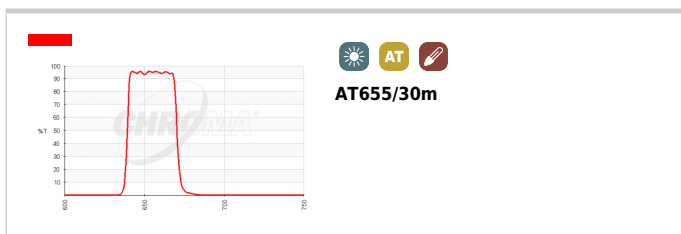
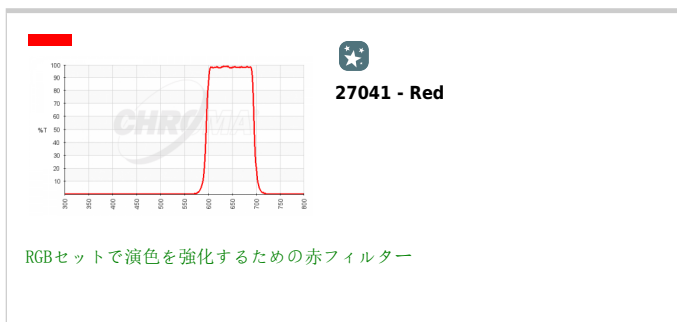
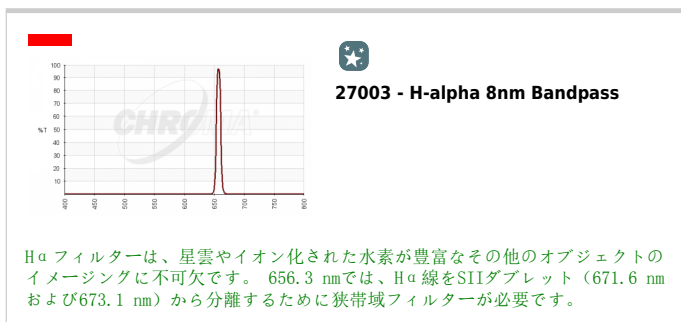
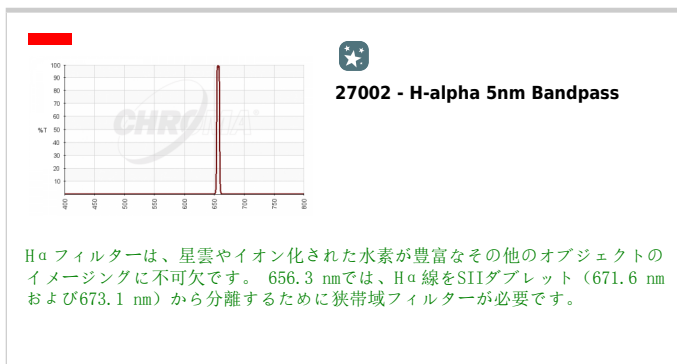
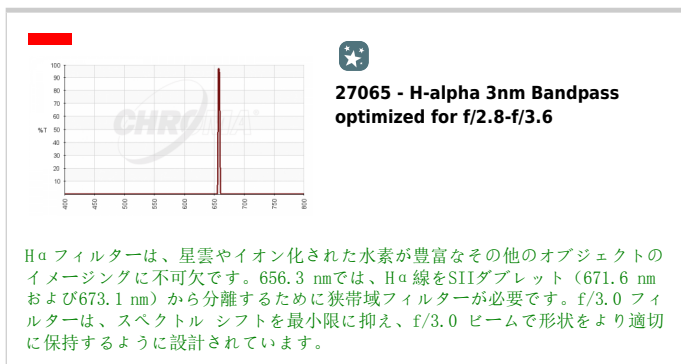
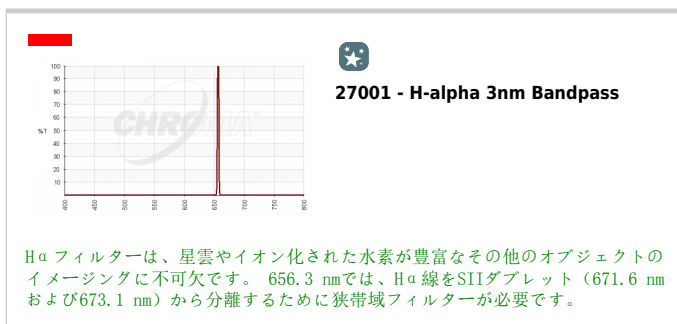
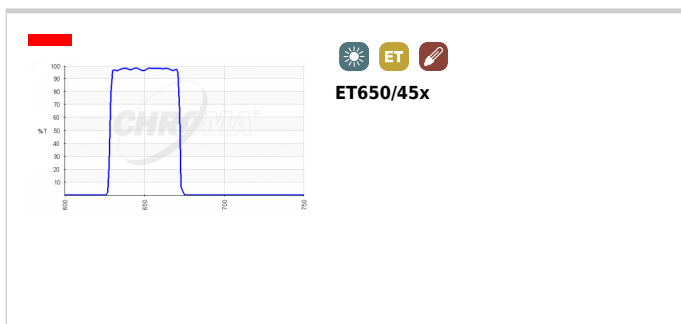
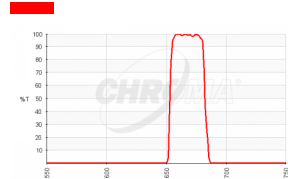
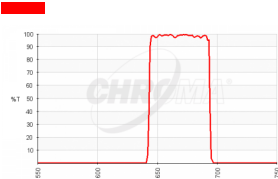
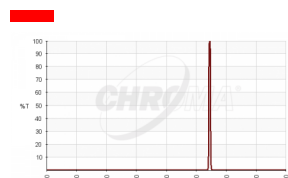
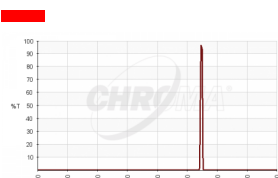
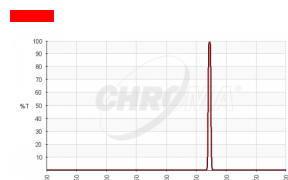
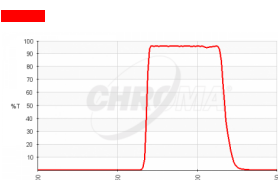
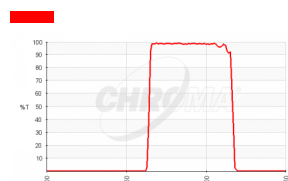
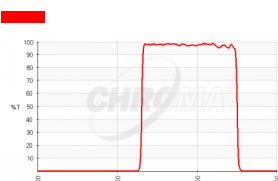
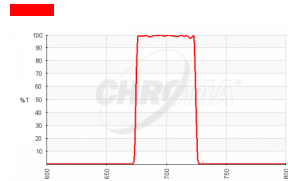
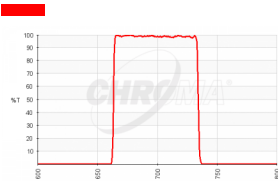
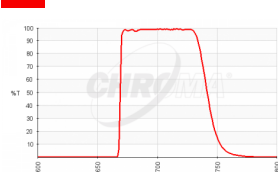
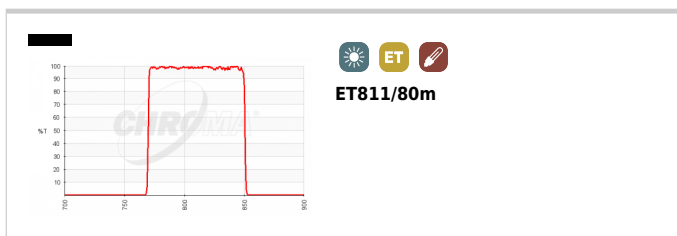
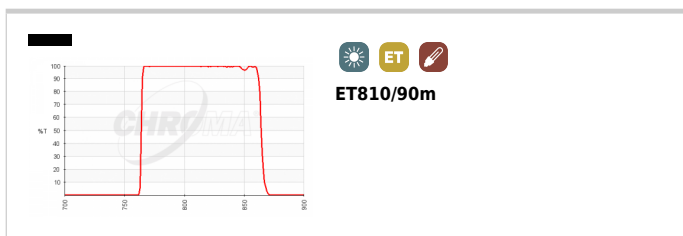
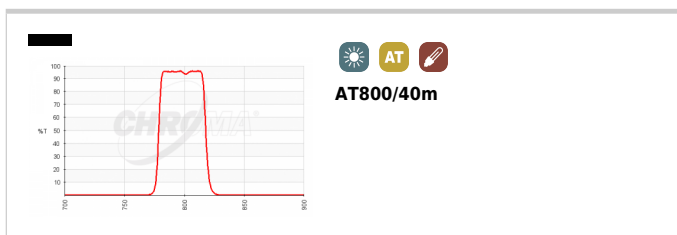
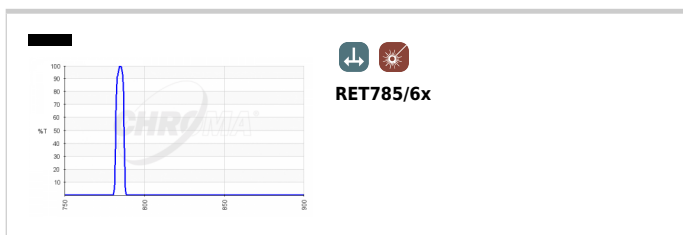
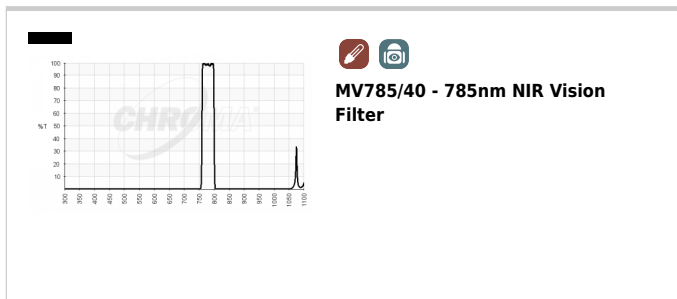
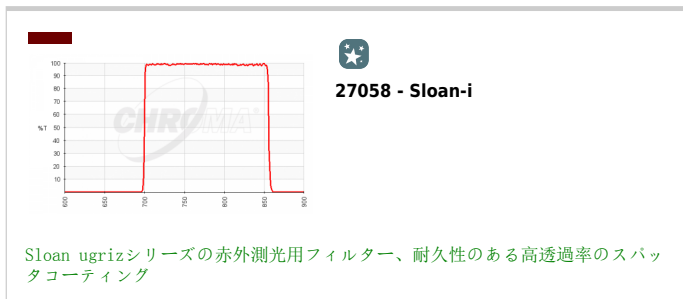
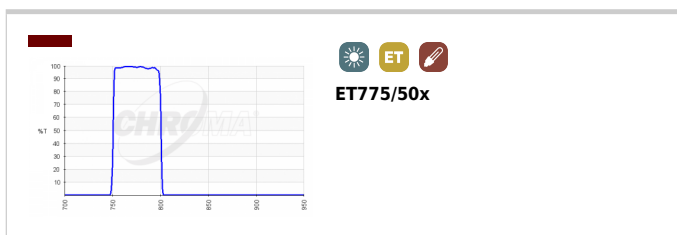
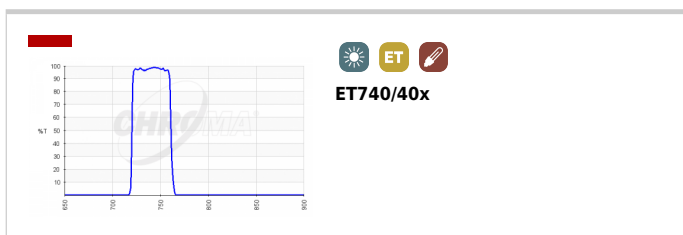
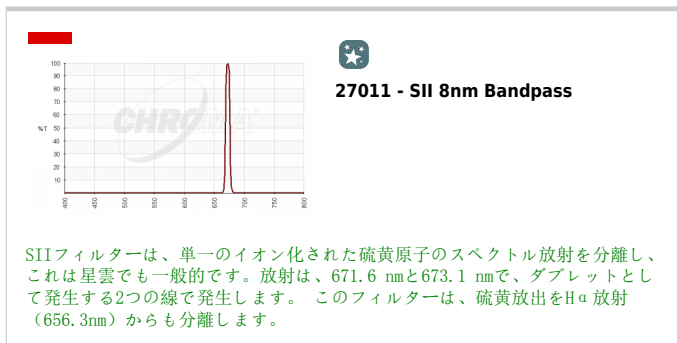
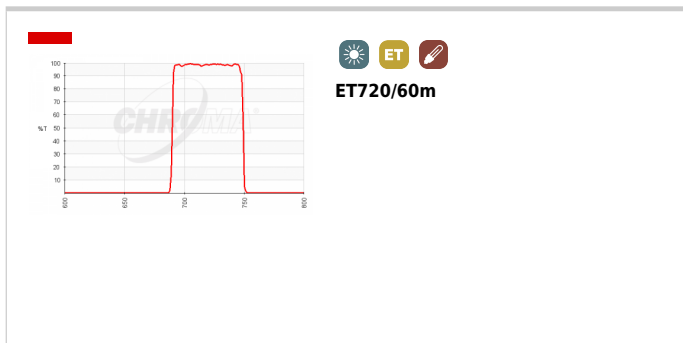
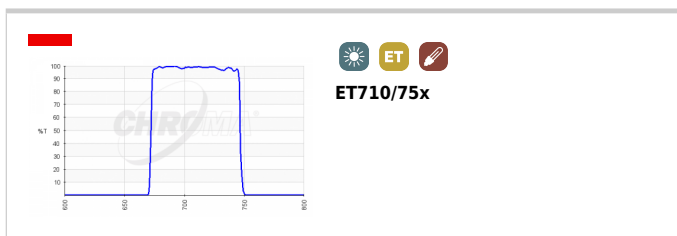
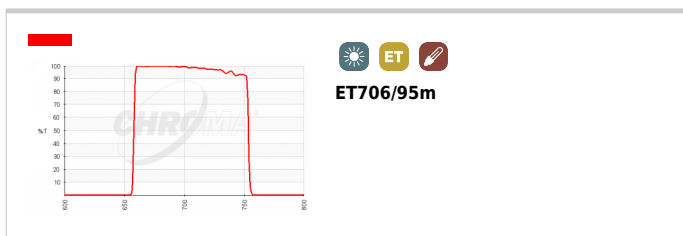
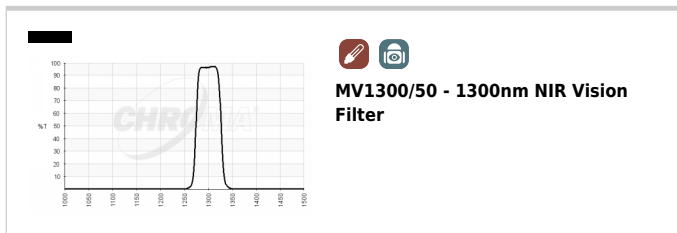
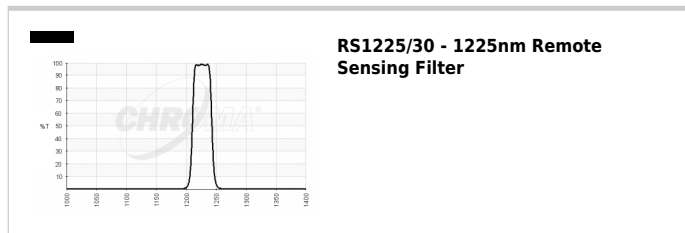
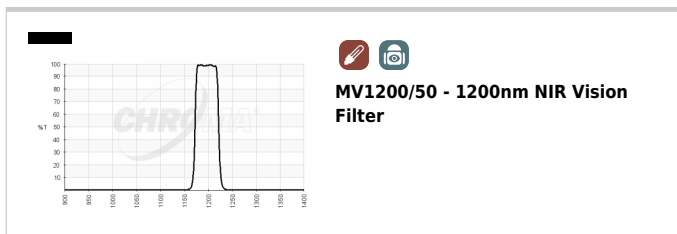
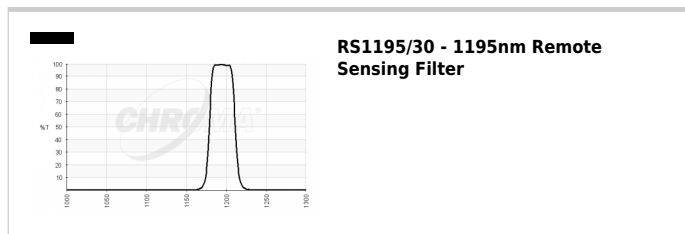
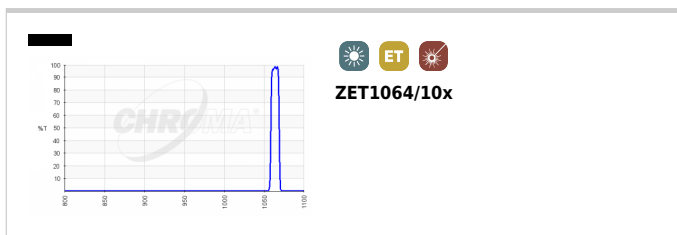
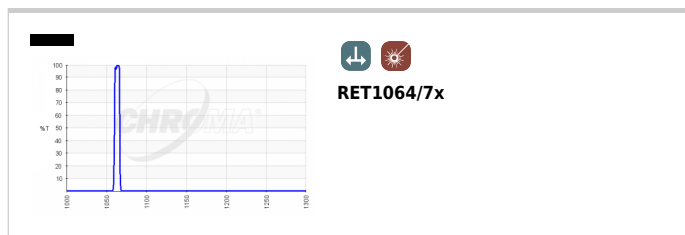
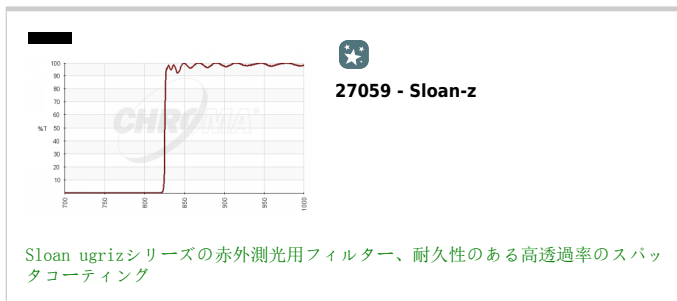
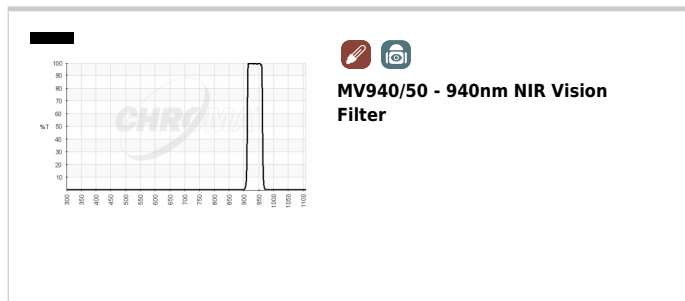
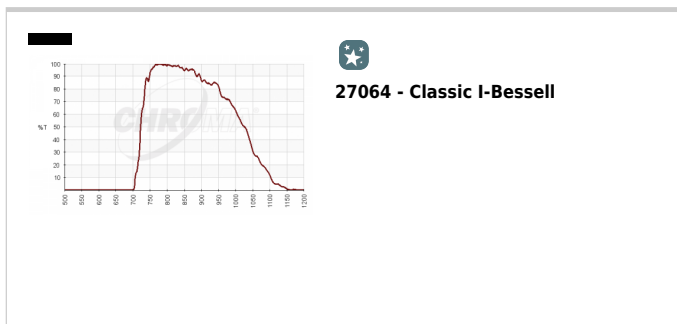
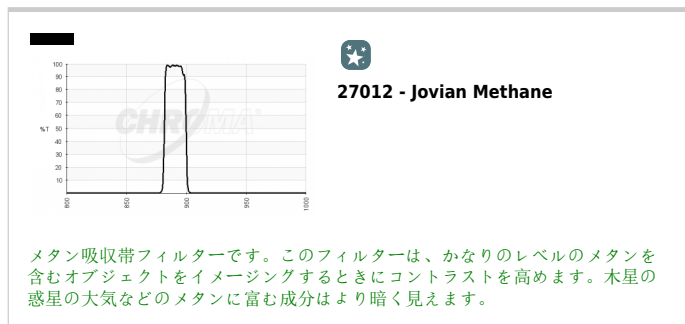
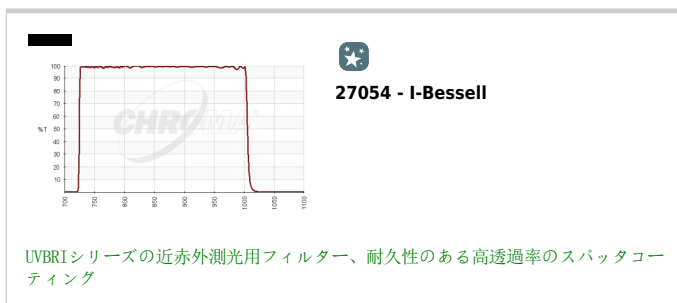
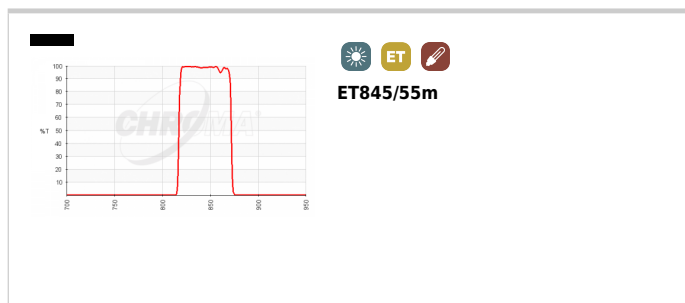
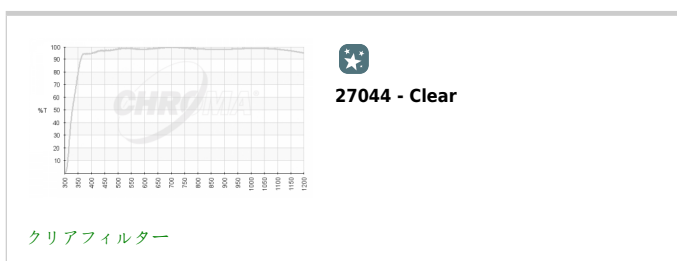
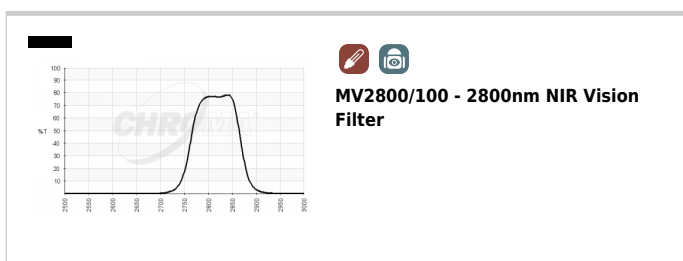
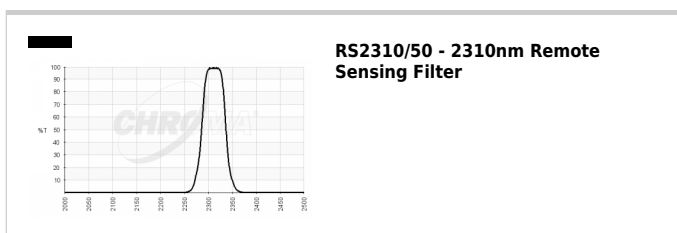
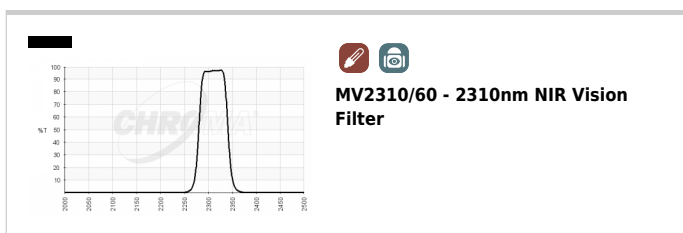
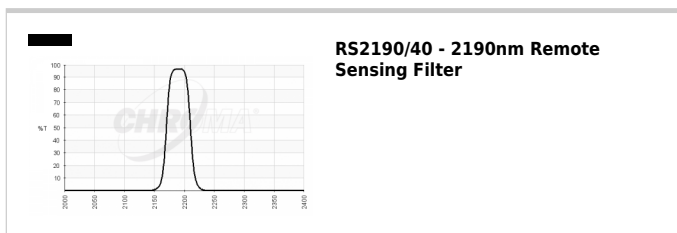
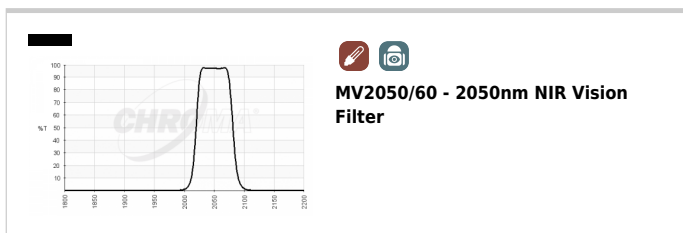
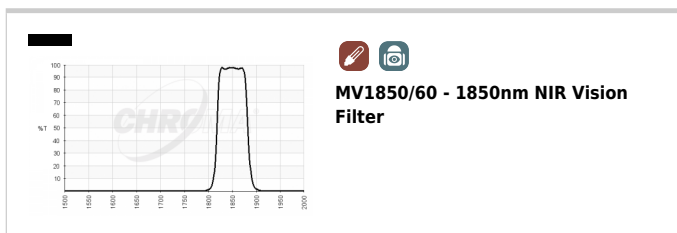
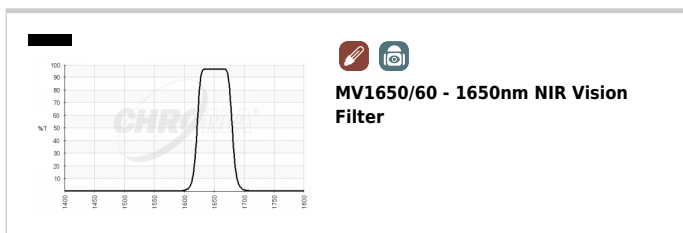
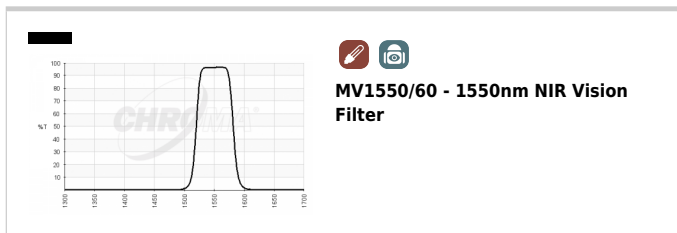
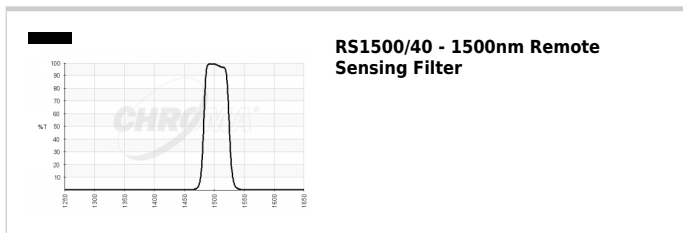
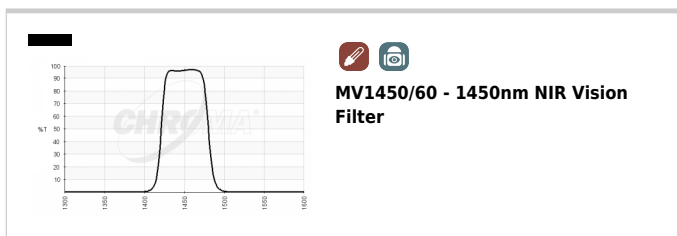
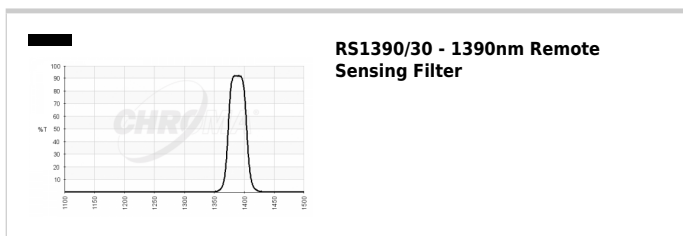
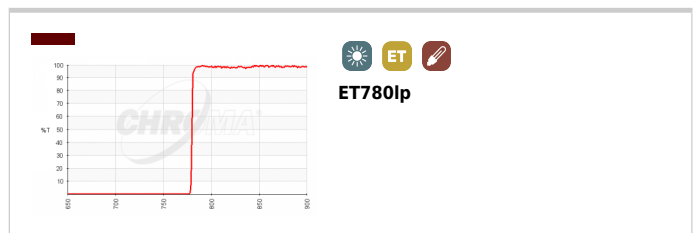
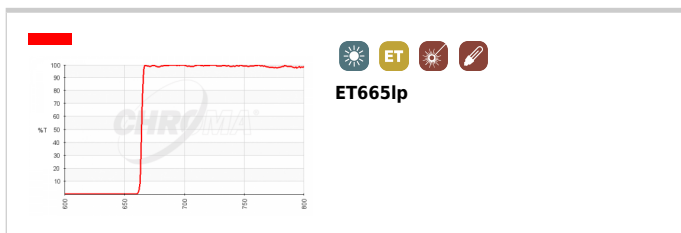
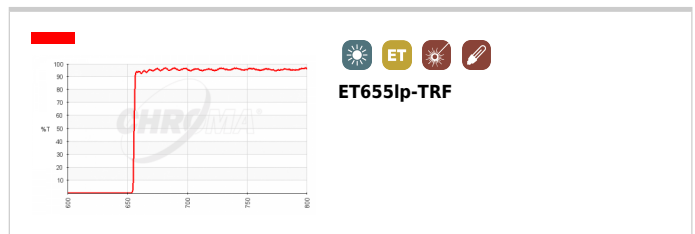
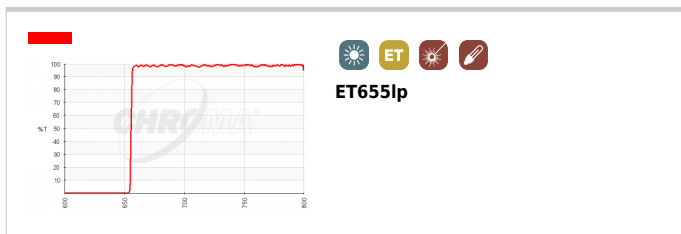
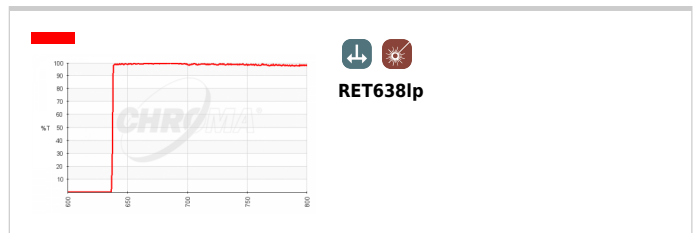
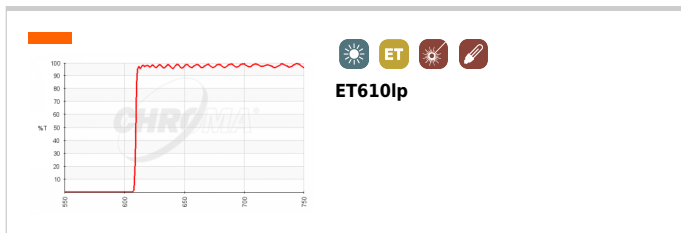
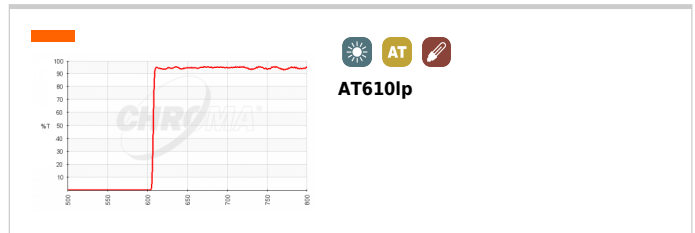
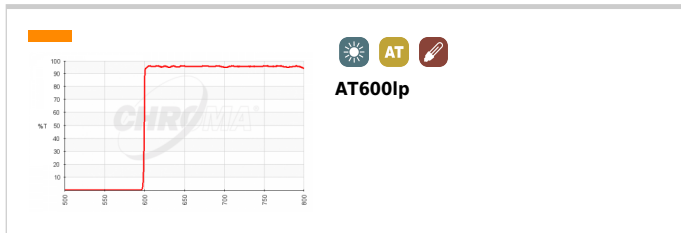
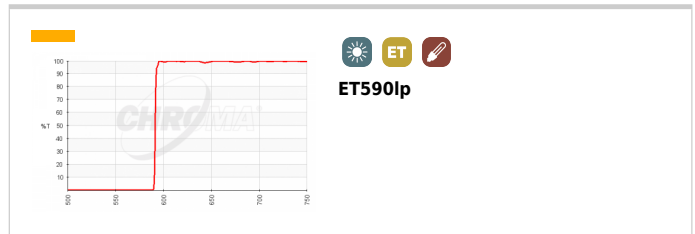
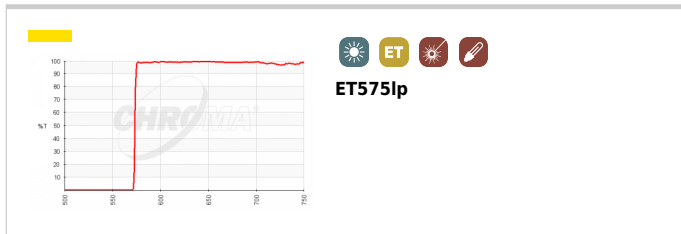
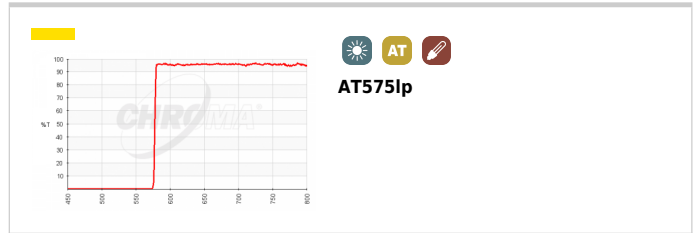
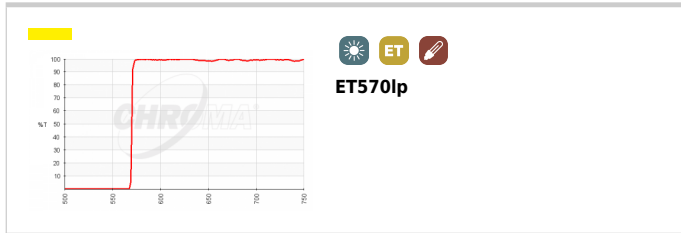
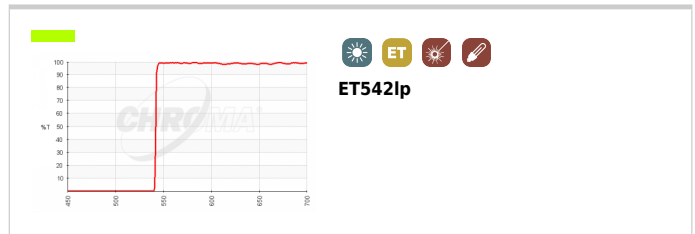
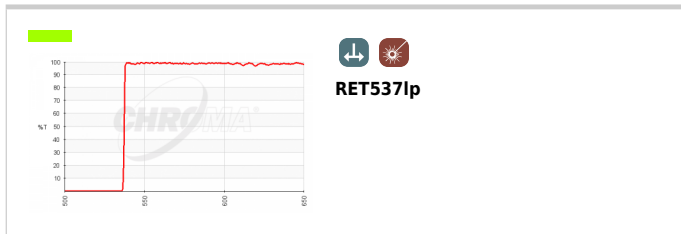
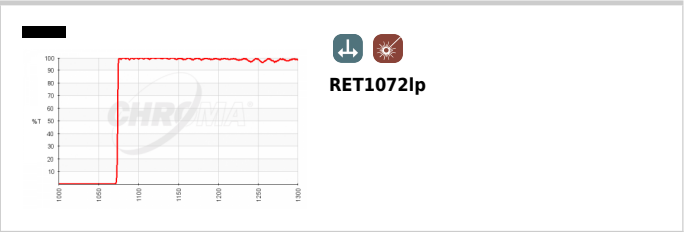
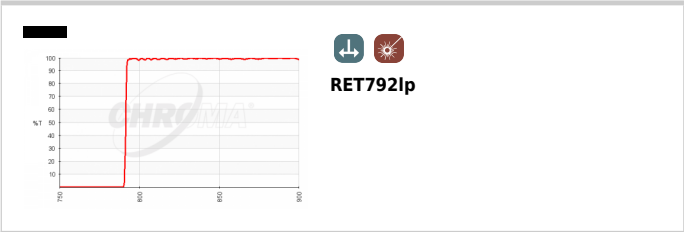
Shortpass Filters

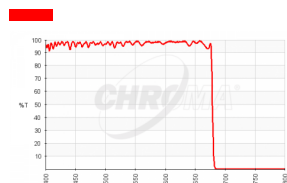**ET680sp-2p8 - Multiphoton shortpass emission filter**

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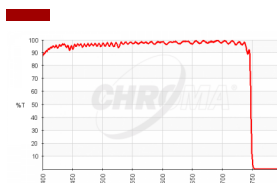
**ET750sp-2p8 - Multiphoton shortpass emission filter**

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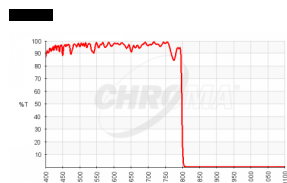**ET800sp-2p/1350 - Multiphoton shortpass emission filter**

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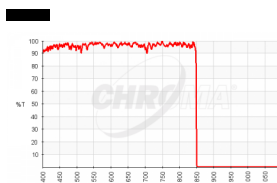
**ET850sp-2p - Multiphoton shortpass emission filter**

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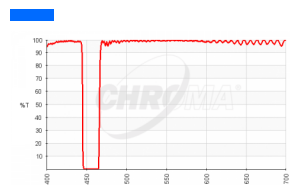
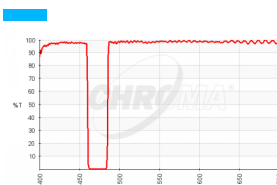
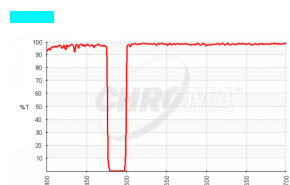
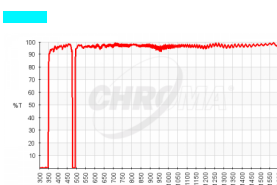
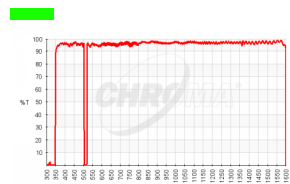
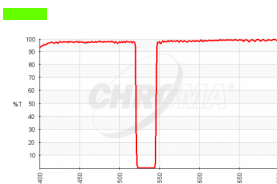
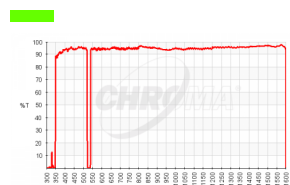
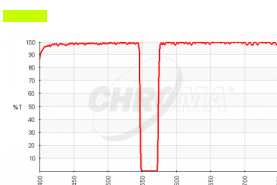
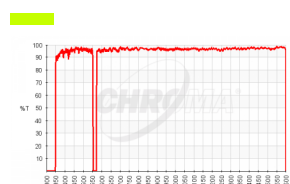
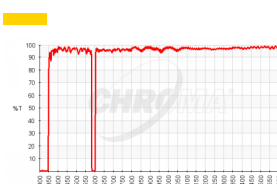
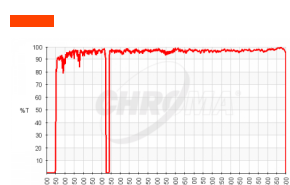
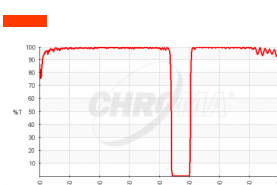
**ZET473NF - 473nm Laser Notch Filter****ZET488NF - 488nm Laser Notch Filter****ZET488TopNotch™ - High Performance Narrow Band 488nm Laser Notch Filter****ZET514TopNotch™ - High Performance Narrow Band 514nm Laser Notch Filter****ZET532NF - 532nm Laser Notch Filter****ZET532TopNotch™ - High Performance Narrow Band 532nm Laser Notch Filter****ZET561NF - 561nm Laser Notch Filter****ZET561TopNotch™ - High Performance Narrow Band 561nm Laser Notch Filter****ZET594TopNotch™ - High Performance Narrow Band 594nm Laser Notch Filter****ZET633TopNotch™ - High Performance Narrow Band 633nm Laser Notch Filter****ZET635NF - 635nm Laser Notch Filter**

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
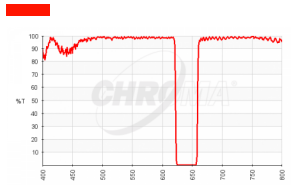
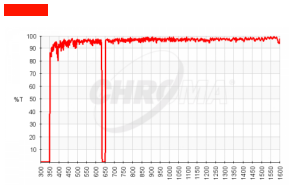
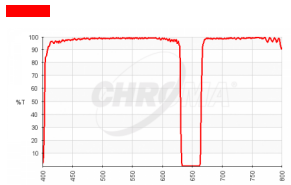
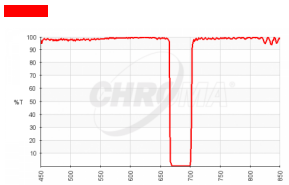
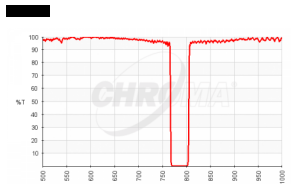
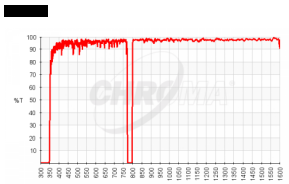
Beamsplitters (50/50, etc.)

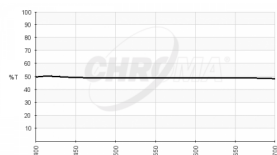
21014 - Silver sputtered non-polarizing 50/50 Vis Beamsplitter

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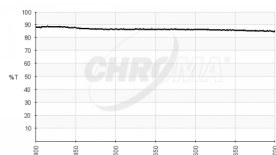
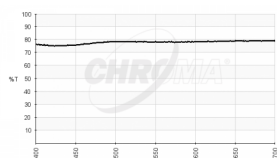
21017 - Sputter-coated, non-polarizing 20/80 Beamsplitter

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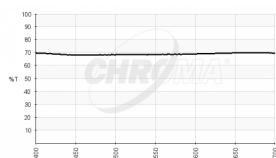
21018 - Sputter-coated, non-polarizing 30/70 Beamsplitter

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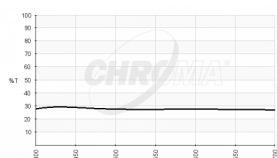
21019 - Sputter-coated, non-polarizing 70/30 Beamsplitter

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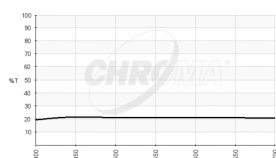21020 - Sputter-coated, non-polarizing 80/20 Beamsplitter

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
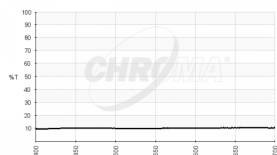
21021 - Sputter-coated, non-polarizing 90/10 Beamsplitter

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Single Edge and Single Band Dichroic Beamsplitters

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








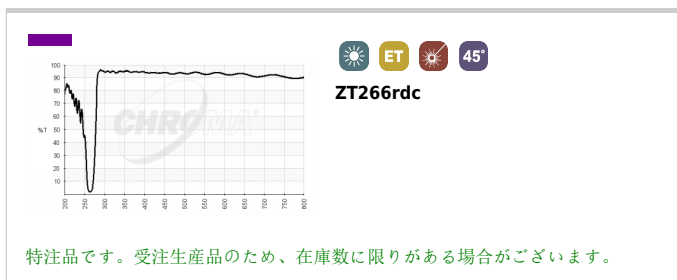
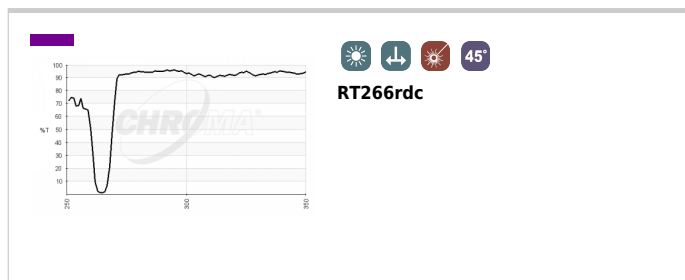
Fluorescence Filter Sets

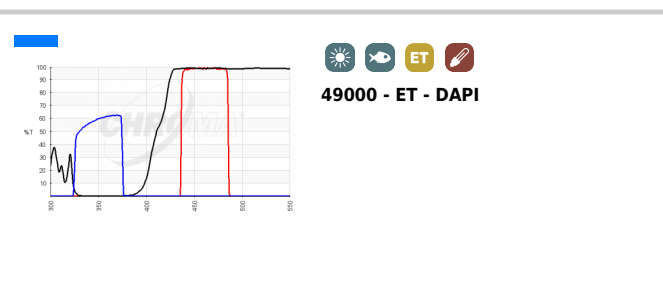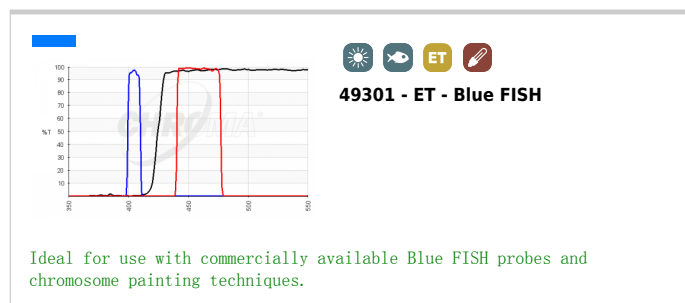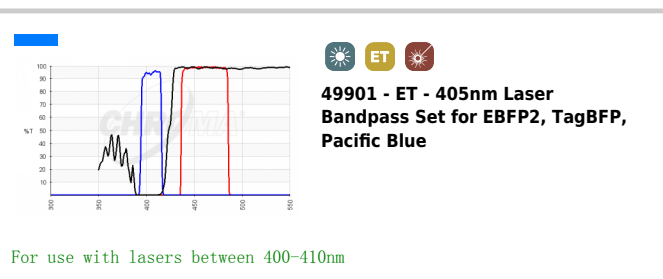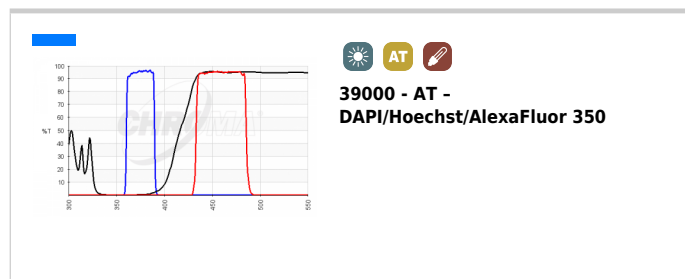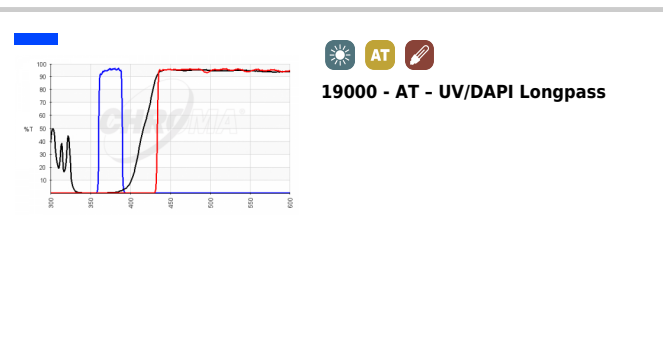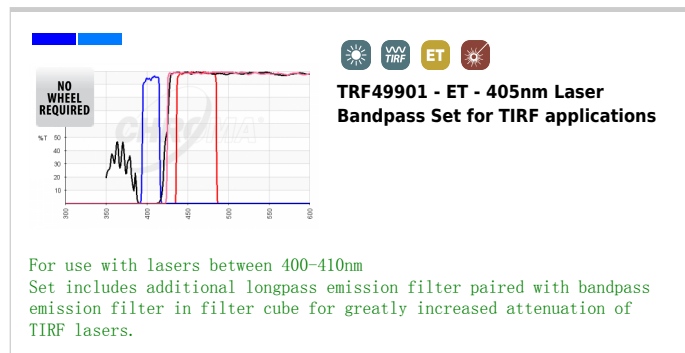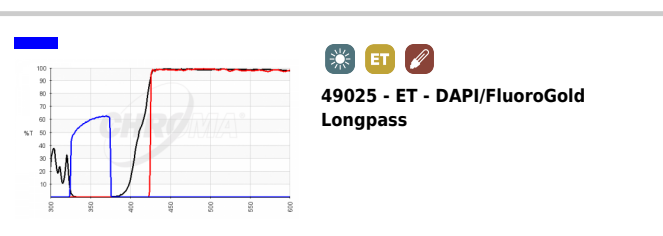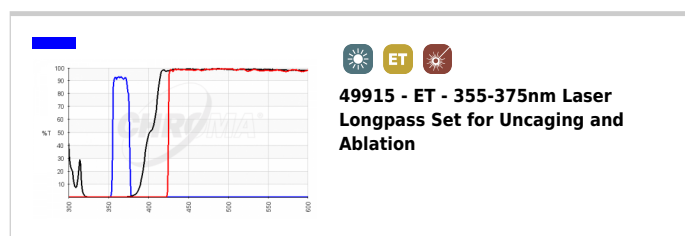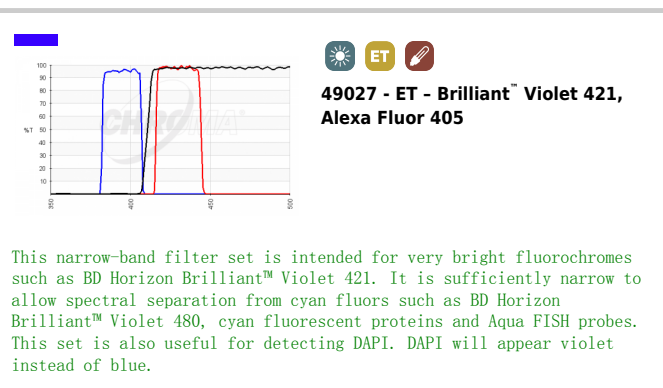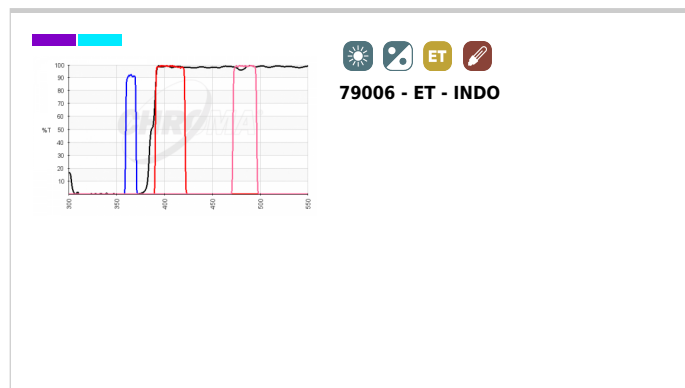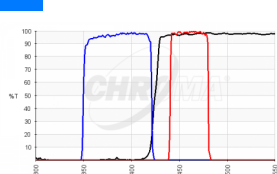
49053 - ET - DAPI for 355-405nm LEDs


Optimized for use with CoolLED light engines.

This set features a wide band excitation filter which is compatible with a broad range of LEDs with peaks between 355-405nm, intended for fluorophores such as DAPI. This set is intended for use with LEDs. If instead, it is used with a broad white light source such as a mercury or metal halide lamp, the excitation filter will transmit very high levels of light to the sample.

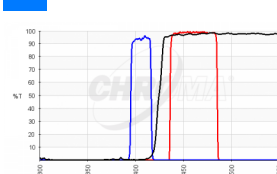
49021 - ET - EBFP2/Coumarin/Attenuated DAPI


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

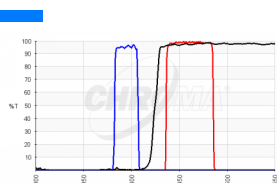49028 - ET - DAPI for 395nm light sources

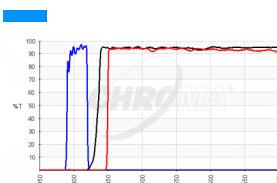
19011 - AT - Coumarin/Pacific Blue Longpass

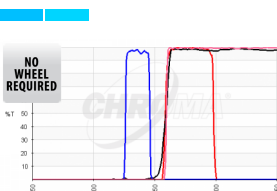TRF49902 - ET - 440nm Laser Bandpass Set for TIRF applications

For use with lasers between 430-445nm


Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

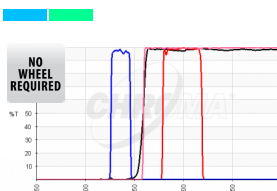
TRF49903 - ET - 440nm Laser Bandpass Set for TIRF applications

For use with lasers between 430-445nm


Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

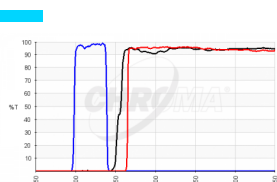19001 - AT - Aqua Longpass

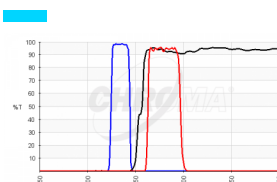
39001 - AT - ECFP/Cerulean

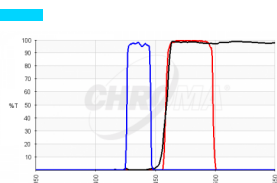
49902 - ET - 440nm Laser Bandpass Set for ECFP, Cerulean

For use with lasers between 430-445nm

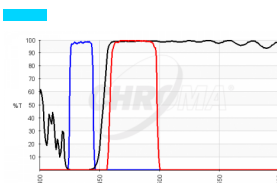
49001 - ET - ECFP


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

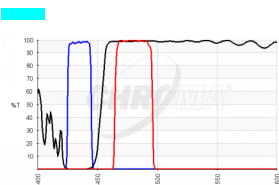**49032 - ET - Brilliant™ Violet 480**





This narrow-band filter set is intended for very bright fluorochromes such as BD Horizon Brilliant™ Violet 480 (BV 480). When used with BV480, it is sufficiently narrow to allow spectral separation from green fluors such as Alexa fluor 488, FITC, GFP and green FISH probes, as well as separation from violet fluors such as Brilliant Violet™ 421.

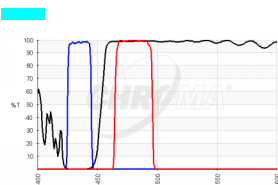**49302 - ET- Aqua FISH**


Ideal for use with commercially available Aqua FISH probes and chromosome painting techniques. Also ideal for use with Brilliant Violet 480 when multiplexing 4 or more colors.

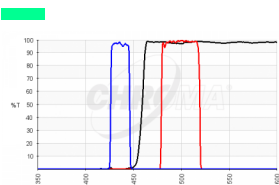**49903 - ET - 440nm Laser Bandpass Set for TFP (Teal FP)**




For use with lasers between 430-445nm

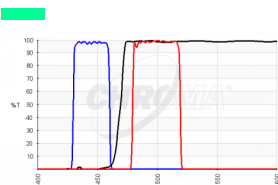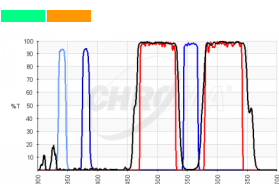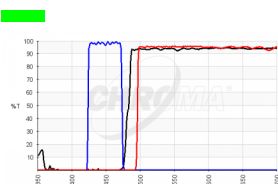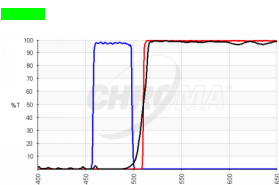**49013 - ET - TFP (Teal FP)****79003 - ET - Fura 2/TRITC****19008 - AT - Auramine Longpass****49012 - ET - FITC/EGFP Longpass**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

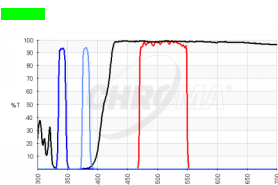**79001 - ET - Fura 2**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

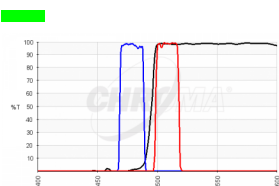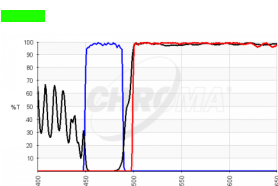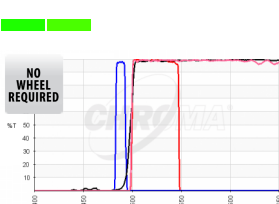**49020 - ET - Narrow Band EGFP to Minimize Autofluorescence****49018 - ET - EGFP Longpass****TRF49904 - ET - 488nm Laser Bandpass Set for TIRF applications**

For use with lasers between 485-491nm
Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

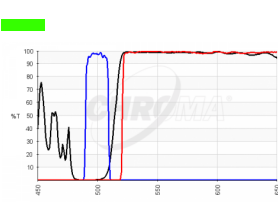
**49029 - ET - EYFP Longpass**

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

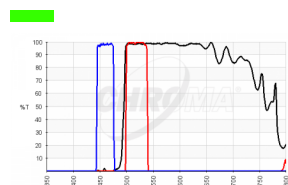
49054 - ET - GFP for 455-465nm LEDs

Optimized for use with CoolLED light engines. This set features an excitation filter which is shifted to slightly shorter wavelengths to be most compatible with blue LEDs with peaks between 455-465nm. GFP will be brighter compared with more standard sets that transmit 470-490nm LEDs.

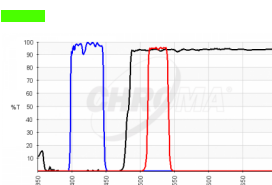
39101 - AT - Qdot 525

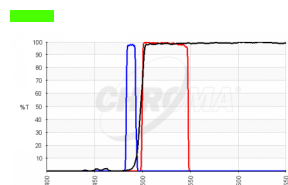
49904 - ET - 488nm Laser Bandpass Set for EGFP, AlexaFluor488, FITC, Fluo3

For use with lasers between 485-491nm

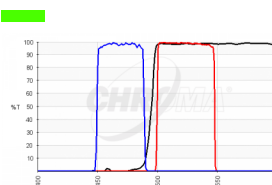
49002 - ET - EGFP (FITC/Cy2)

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

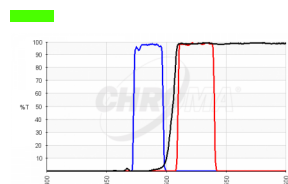
49308 - ET - Green#2 FISH

Ideal for use with commercially available Green FISH probes and chromosome painting techniques. Intended to avoid detection of Gold probes in problematic specimens. If Aqua signals are detected with this set and must be avoided, then please refer to set 49303. Recommended to detect CL-490 in Cell IDx Four Color assay.

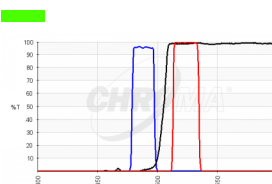
49312 - ET - Green#3 Narrow band FISH


4-color FISH assays utilizing Aqua, Green, Gold & Red are becoming increasingly popular. This FISH filter set is intended for detecting Green probes and completely avoiding both Aqua and Gold probes. It is a narrow band set and Green signals will be slightly less bright, but it will isolate Green signals very well.

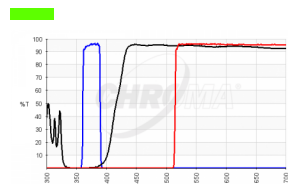
19012 - AT - FluoroGold/Calcofluor White Longpass


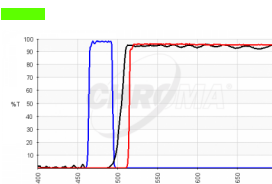
19002 - AT - GFP/FITC Longpass

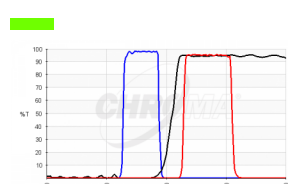
39002 - AT - EGFP/FITC/Cy2/AlexaFluor 488


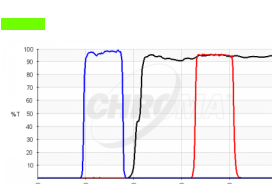
39008 - AT - Lucifer Yellow/Tetracycline/AlexaFluor 430

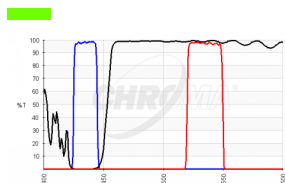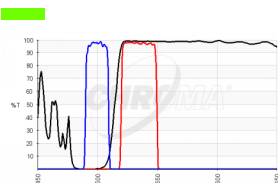**49052 - ET - ECFP/EYFP****49003 - ET - EYFP**




Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

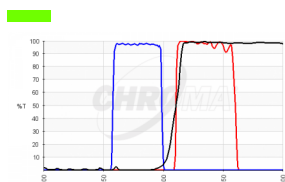**49011 - ET - FITC/Alexa Fluor 488/Fluo3/Oregon Green**



Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

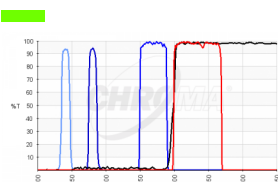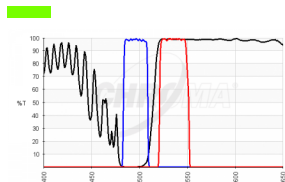**79004 - ET - Fura 2/GFP****49303 - ET - Green#1 FISH**



Ideal for use with commercially available Green FISH probes and chromosome painting techniques. Intended to avoid detection of Aqua probes in problematic specimens. If Gold signals are also detected with this set and must be avoided, then please refer to set 49308.

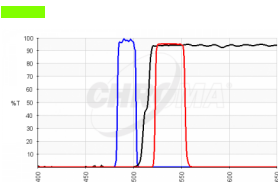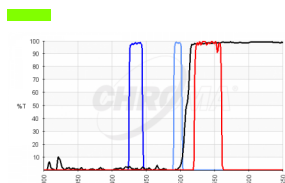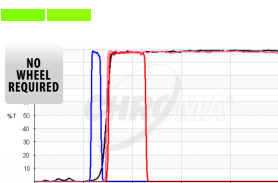**39003 - AT - EYFP/Venus/Citrine****79005 - ET - BCECF****TRF49905 - ET - 514nm Laser Bandpass Set for TIRF applications**

For use with lasers between 511-517nm
Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

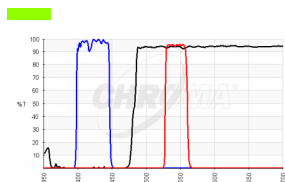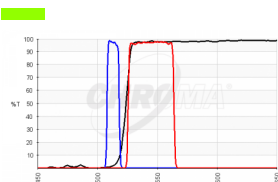
**39102 - AT - Qdot 545****49905 - ET - 514nm Laser Bandpass Set for EYFP, Venus, Citrine, Fluo3**

For use with lasers between 511-517nm

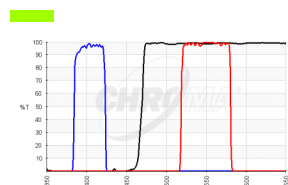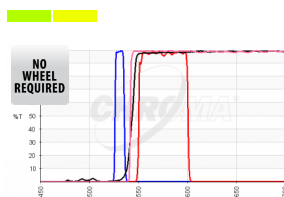**49026 - ET - Tetracycline****TRF49907 - ET - 532nm Laser
Bandpass Set for TIRF applications**





For use with lasers between 529-535nm

Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

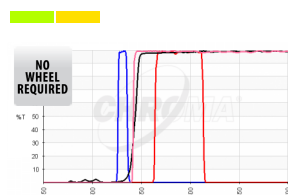
**TRF49908 - ET - 532nm Laser
Bandpass Set for TIRF applications**

For use with lasers between 529-535nm


Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

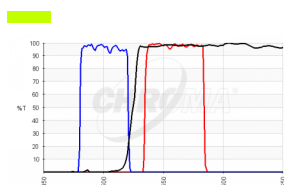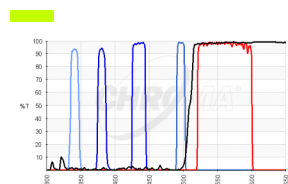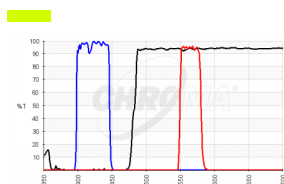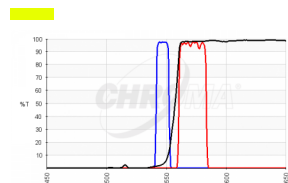**49023 - ET - AlexaFluor
514/ZsYellow/Eosin****79002 - ET - Fura 2/BCECF****39103 - AT - Qdot 565****49304 - ET - Gold FISH**


Ideal for use with commercially available Gold FISH probes and chromosome painting techniques. Will also detect Orange probes.

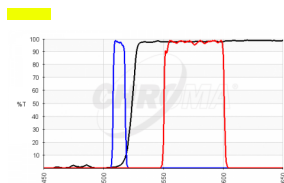
**49906 - ET - 514nm Laser
Bandpass Set for mKO, mOrange2**

For use with lasers between 511-517nm



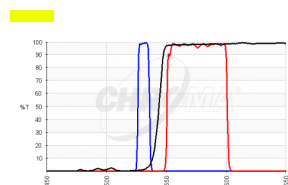
**49907 - ET - 532nm Laser
Bandpass Set for mKO, mOrange2**

For use with lasers between 529-535nm

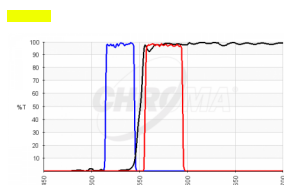**49014 - ET - mKO/mOrange**

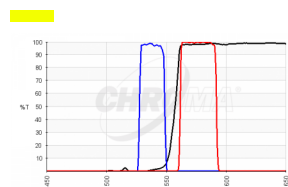
**49309 - ET - Orange#2 FISH**

Ideal for use with commercially available Orange FISH probes and chromosome painting techniques. Provides signals intended to be less bright than those obtained with the 49305 set. Intended to balance signal intensities with Green in problematic specimens which are too bright. This set can be used equally well with both Gold and Orange probes to eliminate the need for 2 different filter sets. Gold and Orange cannot be separated in a microscope.

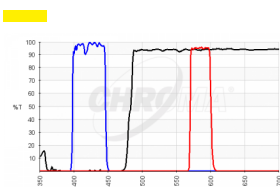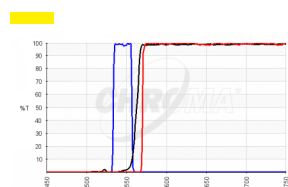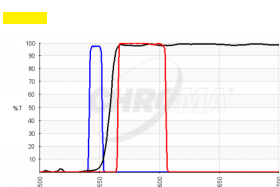
**39104 - AT - Qdot 585****49016 - ET - Cy3/TRITC Longpass****49010 - ET - R&B
Phycoerythrin/mOrange/mKO**

Recommended for use with mercury lamps or metal-halide lamps only.



Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

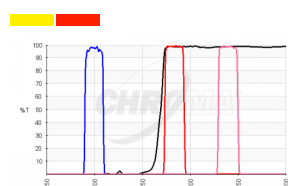**79009 - ET - SNARF with 500nm
excitation**


This filter set is intended for dual-emission ratiometric measurements of the SNARF family of pH indicators when excited with light in the 490-510nm wavelength range.

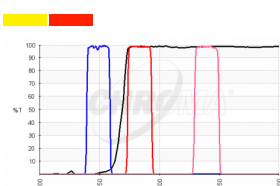**79010 - ET - SNARF with 550nm
excitation**


This filter set is intended for dual-emission ratiometric measurements of the SNARF family of pH indicators when excited with light in the 540-560nm wavelength range.

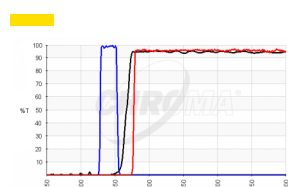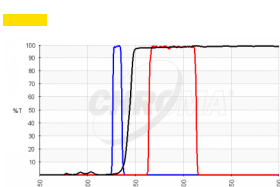**19004 - AT - TRITC/Cy3 Longpass****49908 - ET - 532nm Laser
Bandpass Set for TagRFP,
tdTomato, TRITC, Cy3, AlexaFluor
546**


For use with lasers between 529-535nm

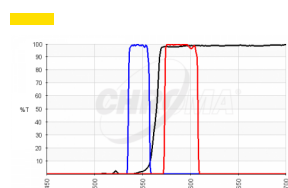**49305 - ET - Orange#1 FISH**

Ideal for use with commercially available Orange FISH probes and chromosome painting techniques. Provides very bright Orange signal. Also detects Gold probes and will detect Red probes. Please refer to set 49309 to avoid detection of Red probes or to provide Orange signals which are less bright. Recommended to detect CL-550 in Cell IDx Four Color assay.

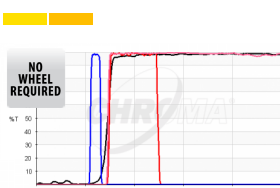**TRF49909 - ET - 561nm Laser
Bandpass Set for TIRF applications**


For use with lasers between 557-563nm
Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

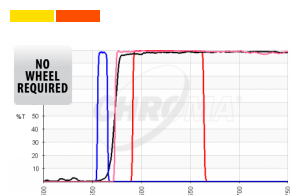**TRF49910 - ET - 561nm Laser
Bandpass Set for TIRF applications**


For use with lasers between 557-563nm

Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

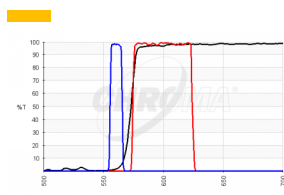
**49909 - ET - 561nm Laser
Bandpass Set for TagRFP,
tdTomato, AlexaFluor 546, TRITC**


For use with lasers between 557-563nm

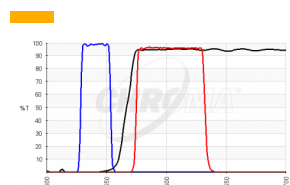**39004 - AT -
TRITC/Cy3/TagRFP/AlexaFluor 546**

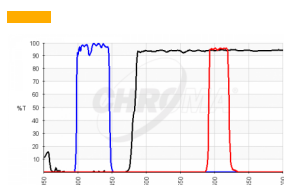
39105 - AT - Qdot 605

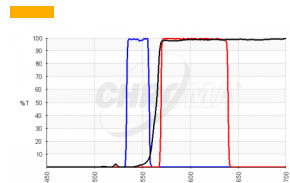
49004 - ET - CY3/TRITC


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

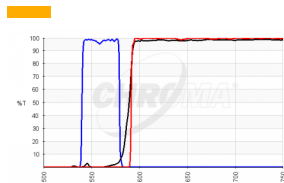**49017 - ET - mCherry/Texas Red
Longpass**

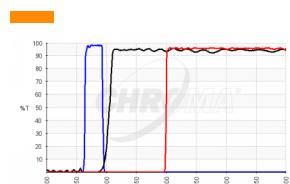
**19010 - AT - FM 1-43/Chlorophyll
Longpass**

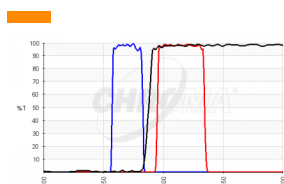
49031 - ET - Alexa Fluor 568, Cy3.5


This is a narrow-band filter set intended for applications where multiple fluorochromes are present in a sample and spectral separation from neighboring fluorochromes is essential. The 49008 set will provide a brighter signal, but may provide less effective spectral separation.

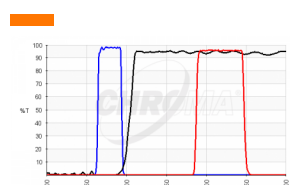**39009 - AT - Acridine Orange/Di-8-
ANEPPS**

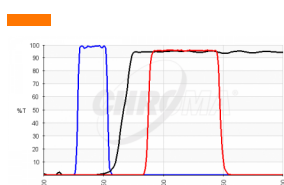
39005 - AT - TRITC-red-shifted


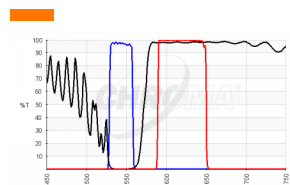
49005 - ET - DSRed (TRITC/Cy3)


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

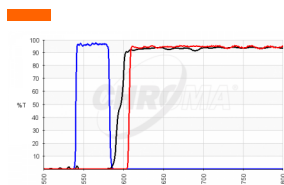
19006 - AT - Texas Red Longpass

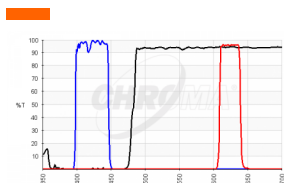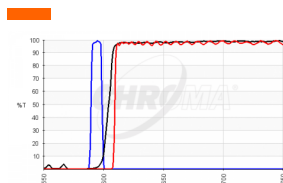**39106 - AT- Qdot 625****49912 - ET - 594nm Laser
Longpass Set for AlexaFluor 594,
mKate2, TurboFP650, NirFP**


For use with lasers between 591-596nm

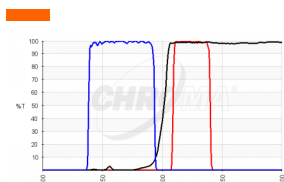**49056 - ET - Narrow
mCherry/Texas Red for 540-580nm
LEDs**


Optimized for use with CoolLED light engines. This set features a wide band excitation filter which transmits a broad range of the relatively weaker LEDs with peaks between 540-580nm. Useful for fluorophores related to rhodamine, Cy3 and fluorescent proteins with similar emission spectra. This set is specifically intended to avoid detection of CY5/AF647 and similar fluors if present in the same sample.

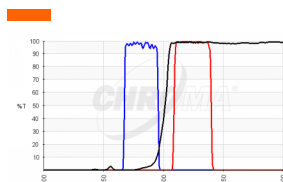**49306 - ET - Red#1 FISH**



Ideal for use with commercially available Red FISH probes and chromosome painting techniques. Provides very bright Red signals Will also detect Orange, and to a lesser extent will also detect Gold signals. Please refer to set 49310 to avoid Gold and Orange.

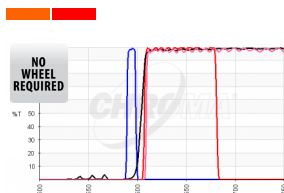**TRF49911 - ET - 594nm Laser
Bandpass Set for TIRF applications**


For use with lasers between 591-597nm
Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

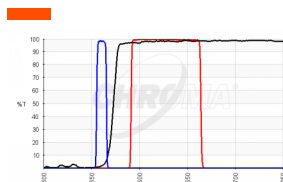**49910 - ET - 561nm Laser
Bandpass Set for AlexaFluor 568,
AlexaFluor 594, Texas Red®,
mCherry**


For use with lasers between 557-563nm

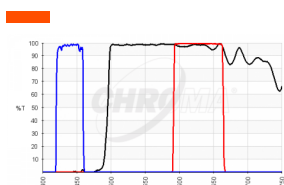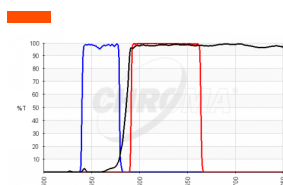**49024 - ET - Keima Red/DiA****49008 - ET - mCherry, Texas Red®**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

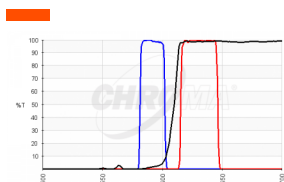**49310 - ET - Red#2 FISH**


Ideal for use with commercially available Red FISH probes and chromosome painting techniques. Red signals will be less bright than those obtained with 49306. Intended to avoid detection of Gold. Use 49311 to avoid detection of Orange in samples which also contain Red.

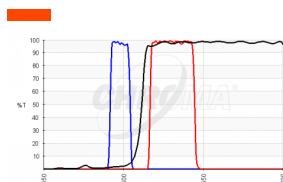**49311 - ET - Red#3 Narrow band
FISH**

Ideal for use with commercially available Red FISH probes and chromosome painting techniques, including mFISH. Red signals will be less bright than those obtained with 49306 or 49310. Intended to completely avoid detection of Orange. This set allows for more combinations in one hybridization, allowing Red and Orange to be run together when using filter set 49309 for Orange.

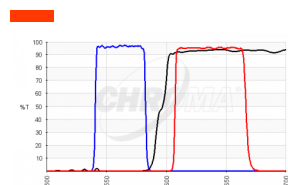
**39010 - AT - Texas
Red/mCherry/AlexaFluor 594**

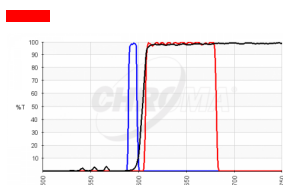
**49911 - ET - 594nm Laser
Bandpass Set for mCherry,
mKate2, AlexaFluor 594, Texas
Red®**

For use with lasers between 591-596nm

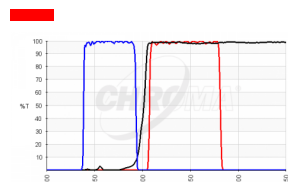
**49055 - ET - Wide mCherry/Texas
Red for 540-580nm LEDs**

Optimized for use with Cooled light engines.
This set features a wide band excitation filter which transmits a broad range of the relatively weaker LEDs with peaks between 540-580nm, resulting in bright fluorescence signal when used with the wide emission filter in this set. Useful for fluorophores related to rhodamine, Cy3 and fluorescent proteins with similar emission spectra. This filter set may also detect Cy5/AF 647 if present.

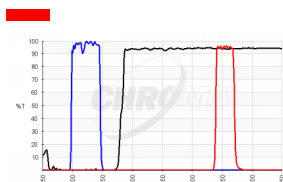
39107 - AT - Qdot 655

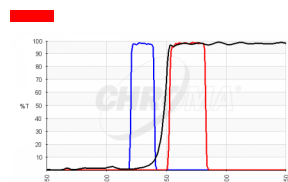
49307 - ET - Far Red FISH

Ideal for use with commercially available Far Red FISH probes and chromosome painting techniques.

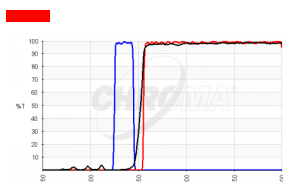
**49913 - ET - 633-640nm Laser
Longpass Set for AlexaFluor 647,
DyLight 649, Atto 647**

For use with lasers between 627-643nm

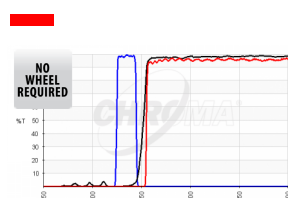
**TRF49913 - ET - 633-640nm Laser
Longpass Set for TIRF applications**


For use with lasers between 627-643nm
This set uses a longpass emission filter with greatly increased blocking for greater attenuation of TIRF lasers without the use of an additional emission filter.

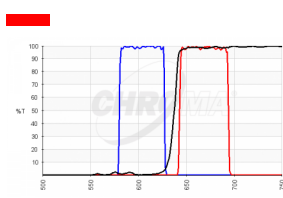49015 - ET - AlexaFluor 633











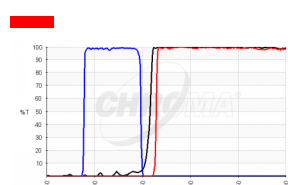
49019 - ET - Cy5 Longpass

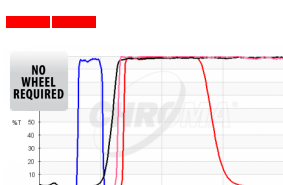**TRF49914 - ET - 640-647nm Laser
Longpass Set for TIRF applications**


For use with lasers between 633-650nm
Set includes additional longpass emission filter paired with bandpass emission filter in filter cube for greatly increased attenuation of TIRF lasers.

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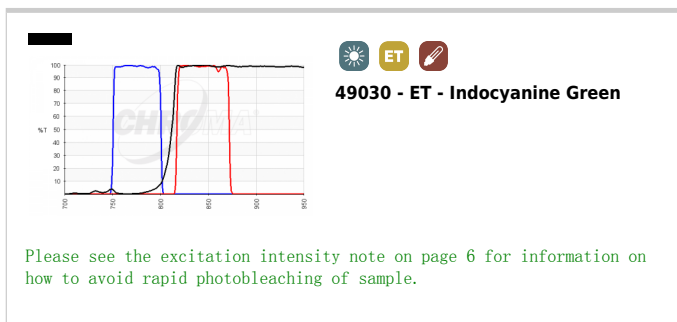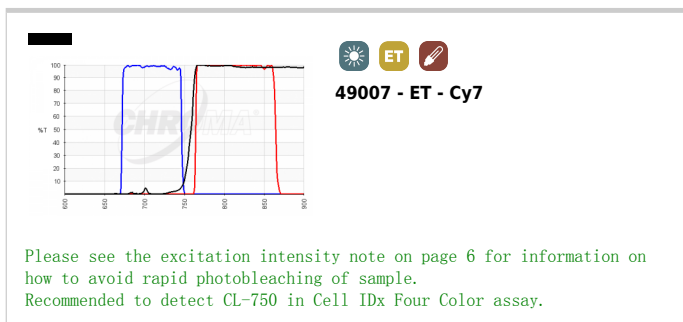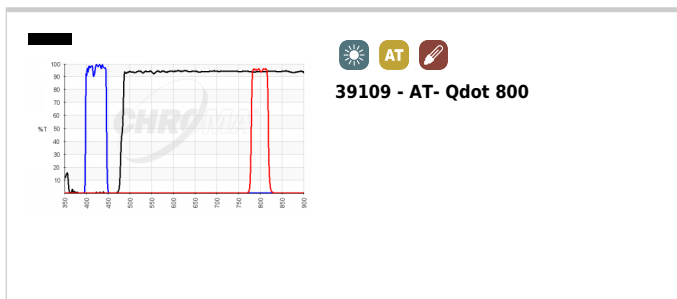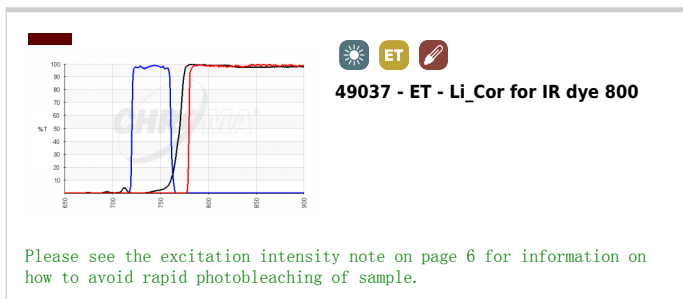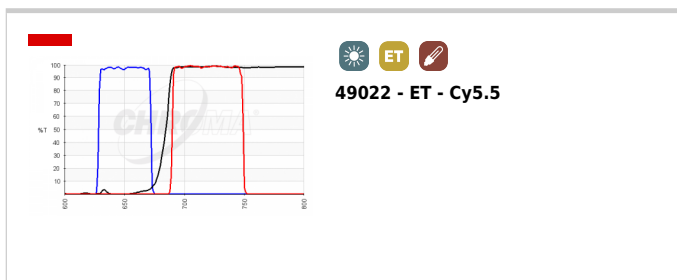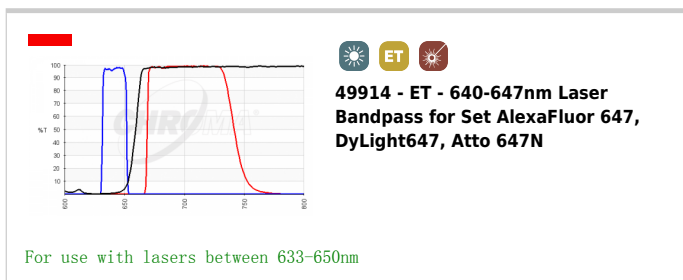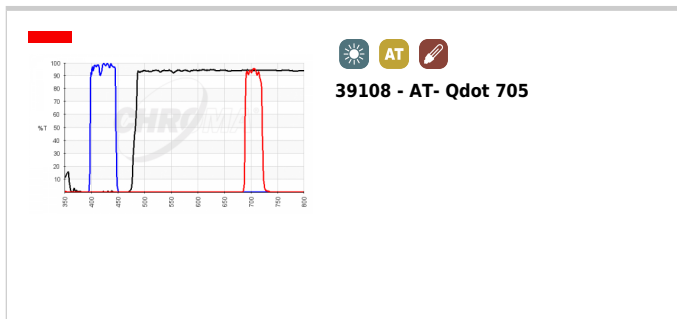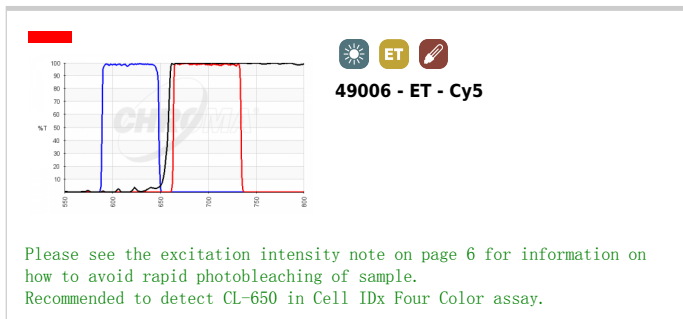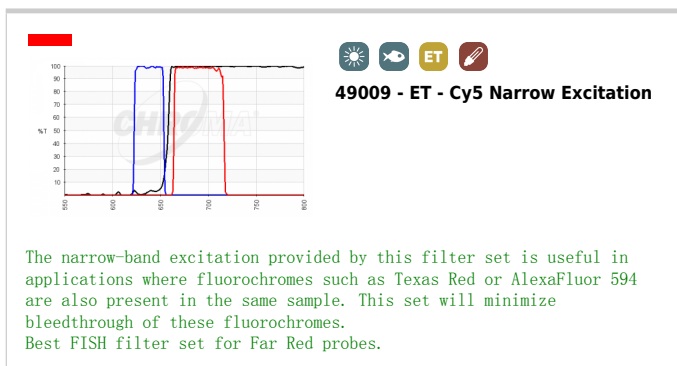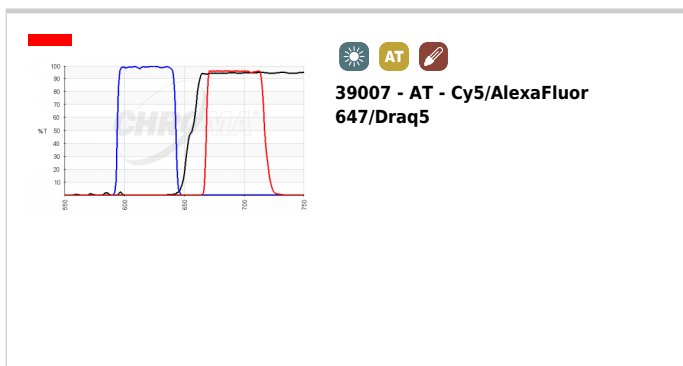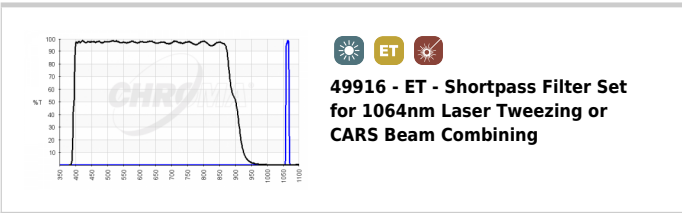
Raman Filter Sets

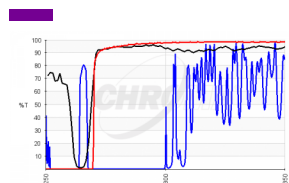
49953 - RT - Raman 266nm Laser Longpass Set

Dedicated Filter set for Raman spectroscopy. The ultra-steep 266nm laser dichroic mirror perfectly matches the ultra-steep longpass emission filter which blocks the 266nm laser. Set includes the narrow-band 266nm laser clean-up filter.

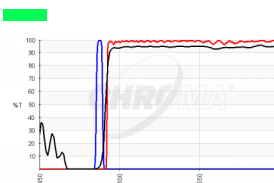
49951 - RT - Raman 488nm Laser Longpass Set

Dedicated Filter set for Raman spectroscopy. The ultra-steep 488nm laser dichroic mirror perfectly matches the ultra-steep longpass emission filter which blocks the 488nm laser. Set includes the narrow-band 488nm laser clean-up filter.

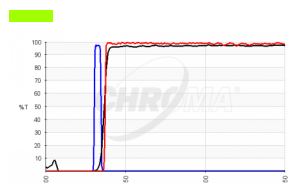
49952 - RT - Raman 532nm Laser Longpass Set

Dedicated Filter set for Raman spectroscopy. The ultra-steep 532nm laser dichroic mirror perfectly matches the ultra-steep longpass emission filter which blocks the 532nm laser. Set includes the narrow-band 532nm laser clean-up filter.

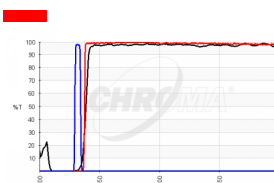
49954 - RT - Raman 633nm Laser Longpass Set

Dedicated Filter set for Raman spectroscopy. The ultra-steep 633nm laser dichroic mirror perfectly matches the ultra-steep longpass emission filter which blocks the 633nm laser. Set includes the narrow-band 633nm laser clean-up filter.

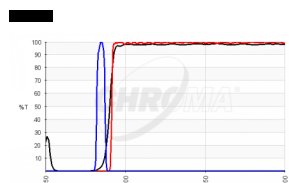49950 - RT - Raman 785nm Laser Longpass Set

Dedicated Filter set for Raman spectroscopy. The ultra-steep 785nm laser dichroic mirror perfectly matches the ultra-steep longpass emission filter which blocks the 785nm laser. Set includes the narrow-band 785nm laser clean-up filter.

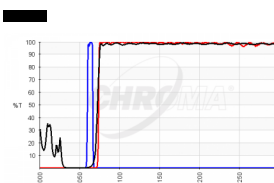
49955 - RT - Raman 1064nm Laser Longpass Set

Dedicated Filter set for Raman spectroscopy. The ultra-steep 1064nm laser dichroic mirror perfectly matches the ultra-steep longpass emission filter which blocks the 1064nm laser. Set includes the narrow-band 1064nm laser clean-up filter.

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
Astronomy Filter Sets

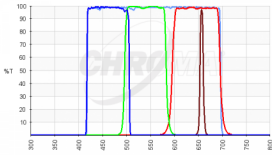**27102 - HaLRGB**

H-alpha, Luminance, Red, Green, and Blue astronomy filter set

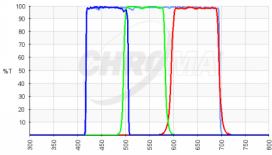
**27101 - LRGB**

Luminance, Red, Green, and Blue astronomy filter set

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
Multi Band Exciter/Multi Band Emitter

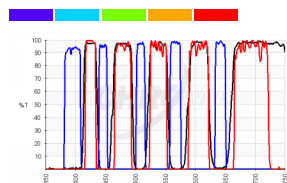
**89904 - ET -
405/445/514/561/640nm Laser
Quinta Band Set**

For use with BD Horizon Brilliant™ Violet 421 and Brilliant™ Violet 480 in 5-color fluorescence imaging scheme. Intended for use with laser excitation, however this set is also compatible with triggerable LED excitation if the appropriate LEDs are used.

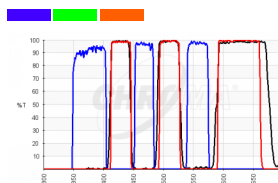
**69401 - ET - 380/55-470/30-557/35
Multi LED set**

Optimized for use with CoolLED light engines. This set provides a wide band for UV excitation, compatible with UV LEDs from 355-405nm. Also compatible with shorter-wavelength blue LEDs from 453-483nm and a wide band for green LEDs from 540-575nm.

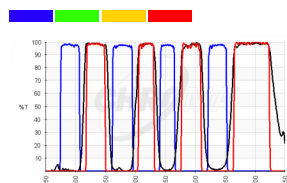
**89402 - ET -
391-32/479-33/554-24/638-31 Multi
LED set**

Quad band catalog set for Multi LED light engines like those manufactured by Lumencor, CoolLED, Excelitas and Cairn Research. Excitation passbands are matched to discrete, triggerable LEDs appropriate for excitation of DAPI/FITC/TRITC/Cy5 and similar fluors. Not intended for use with white light sources.

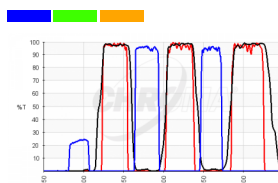
**69013v2 - ET - DAPI/Green/Orange
#1 FISH**

Triple band filter set ideal for simultaneous viewing of DAPI and commercially available Green and Orange FISH probes

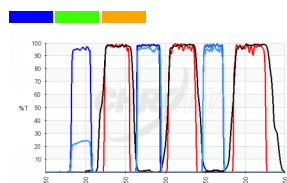
**89084v2 - ET - DAPI/Green/Orange
Multi-Exciter FISH**

This filter set replaces the original 82000 set and the recent 89084 set. It is a tripleband filter set with passbands equivalent to the 69013v2 set. In addition, This set also contains a Green/Orange dualband excitation filter and 3 individual single band excitation filters. This allows the user to detect DAPI/Green/Orange simultaneously with the tripleband exciter, or only Green/Orange with the dualband exciter, or single colors for imaging to a camera. The tripleband excitation filter features an attenuated level of DAPI excitation.

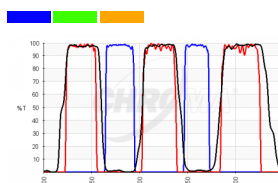
59011 - ET - Green/Orange #1 FISH

Dualband filter set ideal for simultaneous viewing of commercially available Green and Orange FISH probes. Intended to balance the two fluorescence signals by attenuating Orange signals relative to Green when used with mercury lamps and metal halide light sources.

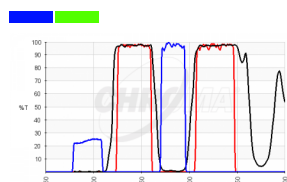
59001v2 - ET - DAPI/Green FISH

Dual band filter set for simultaneous viewing of Green FISH probes and a DAPI counterstain. Features an attenuated level of DAPI excitation to balance the intensity of DAPI and Green fluorescence signals. Replaces the original 59001.

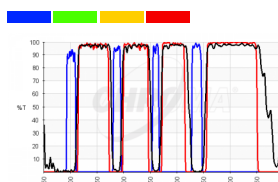
**89901v2 - ET - 405/488/561/640nm
Laser Quad Band Set**

For use with lasers between 400-410nm, 485-492nm, 558-563nm and 632-643nm

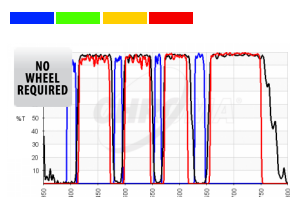
**TRF89901v2 - ET -
405/488/561/640nm Laser Quad
Band Set for TIRF applications**

Set includes no additional emission filters; instead, the quad band emission filter has increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

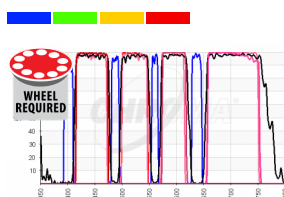
**TRF89901-EMv2 - ET -
405/488/561/640nm Laser Quad
Band Set for TIRF applications**

Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

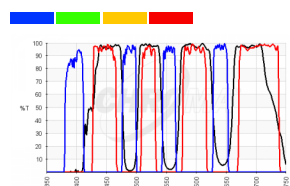
**89401 - ET - DAPI/FITC/TRITC/CYS
Quad**

A version of the 89000 Quad set with a quad band excitation filter, dichroic and emission filter. Useful for simultaneous live viewing of multiple channels or image-splitting applications. May also be paired with individual emission filters from 89000 set for highest SNR image-splitting when acquiring multiple channels simultaneously. May be used with color cameras for standard image acquisition. Use of monochrome cameras requires an image-splitting device.

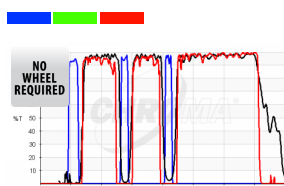
**TRF69901 - ET - 405/488/561nm
Laser Triple Band Set for TIRF
applications**

For use with lasers between 400-410nm, 486-492nm and 558-563nm
Set includes no additional emission filters; instead, the triple band emission filter has greatly increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

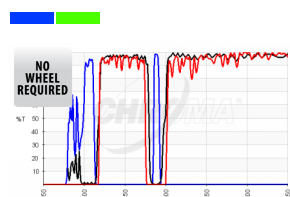
**TRF59901 - ET - 405/488nm Laser
Dual Band Set for TIRF
applications**

For use with lasers between 400-408nm and 485-491nm
Set includes no additional emission filters; instead, the dual band emission filter has increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

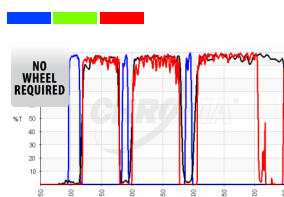
**TRF69902 - ET - 405/488/594nm
Laser Triple Band Set for TIRF
applications**

For use with lasers between 400-412nm, 486-492nm and 591-596nm
Set includes no additional emission filters; instead, the triple band emission filter has greatly increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

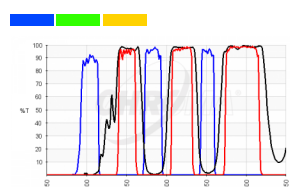
69000 - ET - DAPI/FITC/TRITC

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

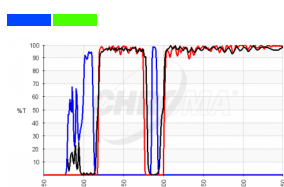
**59901 - ET - 405/488nm Laser Dual
Band Set**

For use with lasers between 400-408nm and 485-491nm

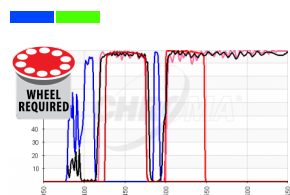
**TRF59901-EM - ET - 405/488nm
Laser Dual Band Set for TIRF
applications**

For use with lasers between 400-408nm and 485-491nm
Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

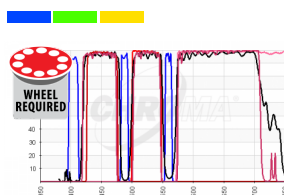
**TRF69901-EM - ET -
405/488/561nm Laser Triple Band
Set for TIRF applications**

For use with lasers between 400-410nm, 486-492nm and 558-563nm
Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

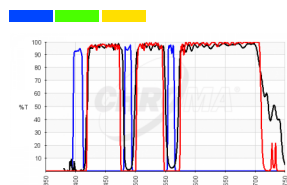69901 - ET - 405/488/561nm Laser Triple Band set

For use with lasers between 400-410nm, 486-492nm and 558-563nm

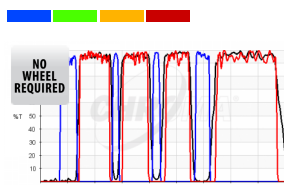
TRF89902 - ET - 405/488/561/647nm Laser Quad Band Set for TIRF applications

Set includes no additional emission filters; instead, the quad band emission filter has increased blocking for greater attenuation of TIRF lasers without the use of additional filters.


89902 - ET - 405/488/561/647nm Laser Quad Band Set

For use with lasers between 400-410nm, 486-492nm, 559-563nm and 640-658nm

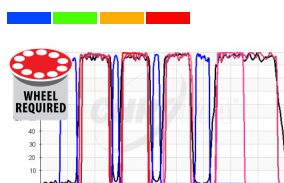
TRF89902-EM - ET - 405/488/561/647nm Laser Quad Band Set for TIRF applications

Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

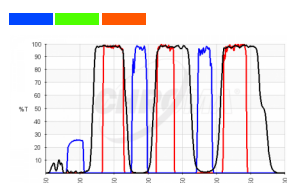
69015 - ET - DAPI/Green/Red FISH

Triple band filter set ideal for simultaneous viewing of DAPI and commercially available Green and Red FISH probes. Attenuated level of DAPI excitation provides better balance of colors so that DAPI signal is not overwhelming.

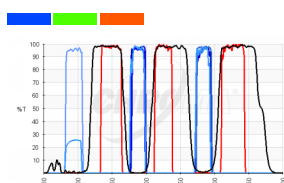
89085 - ET - DAPI/Green/Red Multi-Exciter FISH

This filter set replaces the original 83000 set and the recent 89083 set. Tripleband filter set with passbands equivalent to the 69015 set. In addition, This set also contains a Green/Red dualband excitation filter and 3 individual single band excitaiton filters. This allows the user to detect DAPI/Green/Red simultaneously with the tripleband exciter, or only Green/Red with the dualband exciter, or single colors for imaging to a camera.

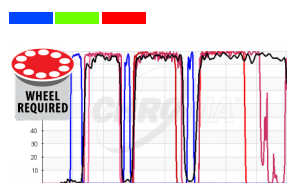
TRF69902-EM - ET - 405/488/594nm Laser Triple Band Set for TIRF applications

For use with lasers between 400-412nm, 486-492nm and 591-596nm
Set includes additional emission filters for use with external emission filter wheel for signal selection and greater attenuation of TIRF lasers.

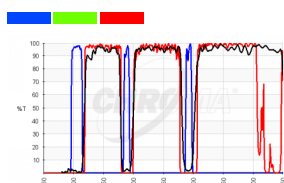
69902 - ET - 405/488/594nm Laser Triple Band set


For use with lasers between 400-412nm, 486-492nm and 591-596nm

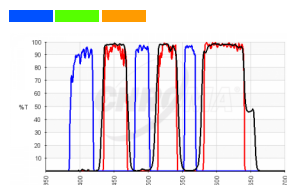**69010 - ET - DAPI/FITC/Cy3**


This triple band set is intended to provide an attenuated Cy3/TRITC/AF546 signal in applications where these orange channel signals overwhelm the green channel signals, and compensation is required.

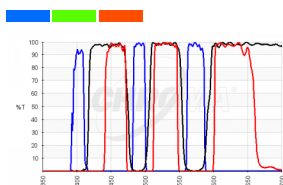**69002 - ET - DAPI/FITC/Texas Red®**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

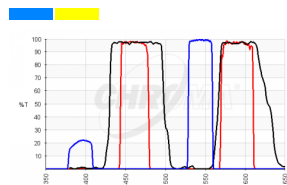**59002v2 - ET - DAPI/Orange FISH**


Dual band filter set for simultaneous viewing of OrangeFISH probes and a DAPI counterstain. Features an attenuated level of DAPI excitation to balance the intensity of DAPI and Orange fluorescence signals. Replaces the original 59002.

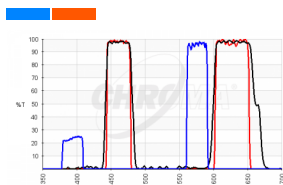**59003v2 - ET - DAPI/Red FISH**


Dual band filter set for simultaneous viewing of Red FISH probes and a DAPI counterstain. Features an attenuated level of DAPI excitation to balance the intensity of DAPI and Red fluorescence signals. Replaces the original 59003.

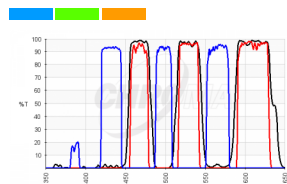**69014 - ET - DAPI/Aqua/Green/Orange FISH**


Similar to 69011 filter set, except for excitation filter. Excitation filter is a quad band, with an attenuated level of excitation for DAPI. This allows for a dim DAPI counterstain that does not overwhelm the Aqua signals. DAPI and Aqua use the same emission band for fluorescence. For simultaneous viewing, not for use with camera.

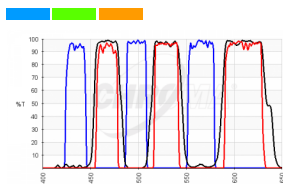**69011 - ET-Aqua/Green/Orange FISH**


Triple band filter set ideal for simultaneous viewing of commercially available Aqua, Green and Orange FISH probes.

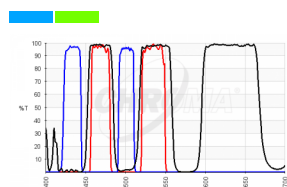**59017 - ET - ECFP/EYFP**

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

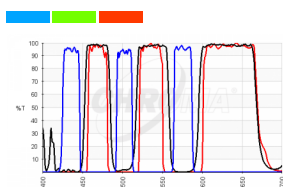
**69008 - ET - ECFP/EYFP/mCherry**

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample. This set is also ideal as a tripleband filter set in FISH applications for simultaneous viewing of commercially available Aqua, Green and Red FISH probes.

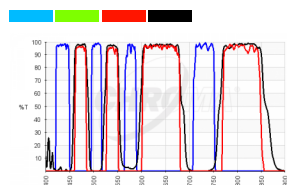
**89403 - ET -
436-28/506-21/578-24/730-40 Multi
LED set**

Quad band catalog set for Multi LED light engines like those manufactured by Lumencor, CoolLED, Excelitas and Cairn Research. Excitation passbands are matched to discrete, triggerable LEDs appropriate for excitation of CFP/YFP/RFP/Cy7 and similar fluors. Not intended for use with white light sources.

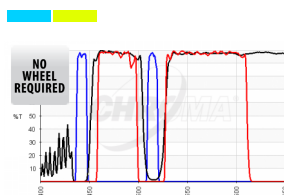
**TRF59902 - ET - 442/514nm Laser
Dual Band Set for TIRF
applications**

For use with lasers between 439-445nm and 512-517nm. Set includes no additional emission filters; instead, the dual band emission filter has increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

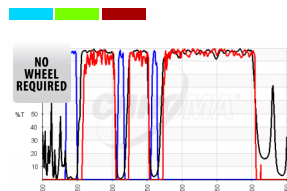
**TRF69904 - ET - 442/514/561nm
Laser Triple Band Set for TIRF
applications**

For use with lasers between 437-448nm, 511-516nm and 559-564nm. Set includes no additional emission filters; instead, the triple band emission filter has greatly increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

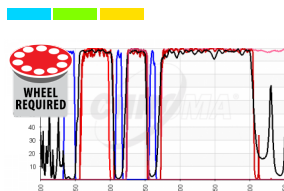
**TRF69904-EM - ET -
442/514/561nm Laser Triple Band
Set for TIRF applications**

For use with lasers between 437-448nm, 511-516nm and 559-564nm. Set includes additional emission filters for use with external emission filter wheel for signal selection and greater attenuation of TIRF lasers.

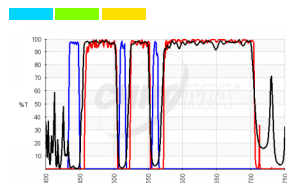
**69904 - ET - 442/514/561nm Laser
Triple Band set**

For use with lasers between 437-448nm, 511-516nm and 559-564nm.

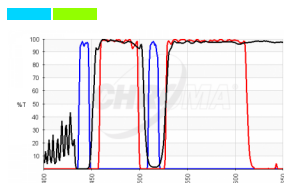
**59902 - ET - 442/514nm Laser Dual
Band Set**

For use with lasers between 439-445nm and 512-517nm.

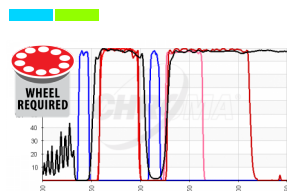
**TRF59902-EM - ET - 442/514nm
Laser Dual Band Set for TIRF
applications**

For use with lasers between 439-445nm and 512-517nm. Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

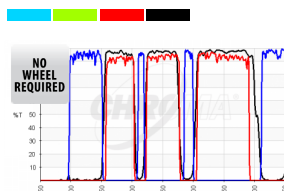
**TRF69905 - ET - 442/514/594nm
Laser Triple Band Set for TIRF
applications**

For use with lasers between 400-449nm, 511-518nm and 588-597nm. Set includes no additional emission filters; instead, the triple band emission filter has greatly increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

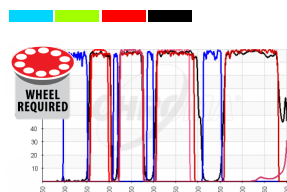
**TRF69905-EM - ET -
442/514/594nm Laser Triple Band
Set for TIRF applications**

For use with lasers between 400-449nm, 511-518nm and 588-597 nm. Set includes additional emission filters for use with external emission filter wheel for signal selection and greater attenuation of TIRF lasers.

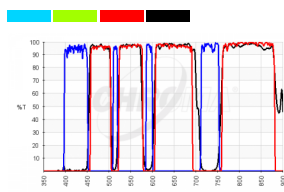
**69905 - ET - 445/514/594nm Laser
Triple Band set**


For use with lasers between 405-448nm, 511-516nm and 589-596nm.

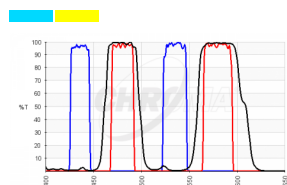**59033 - ET - Aqua/Gold FISH**


Dualband filter set ideal for simultaneous viewing of commercially available Aqua and Gold FISH probes.

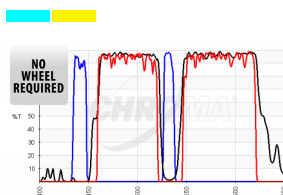**TRF59903 - ET - 442/532nm Laser Dual Band Set for TIRF applications**


For use with lasers between 438-445nm and 530-535nm
Set includes no additional emission filters; instead, the dual band emission filter has increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

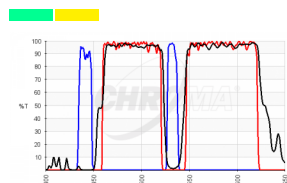**59903 - ET - 442/532nm Laser Dual Band Set**


For use with lasers between 438-445nm and 530-535nm

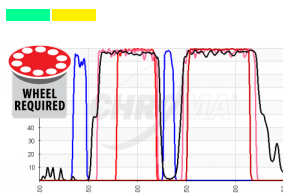**TRF59903-EM - ET - 442/532nm Laser Dual Band Set for TIRF applications**


For use with lasers between 438-445nm and 530-535nm
Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

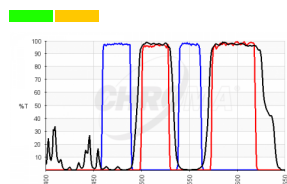**59012 - ET - Green/Orange #2 FISH**


Dualband filter set ideal for simultaneous viewing of commercially available Green and Orange FISH probes. Intended to balance the two fluorescence signals by increasing Orange signal intensity relative to Green when used with LED sources which tend to be weaker in 540-560nm wavelength region.

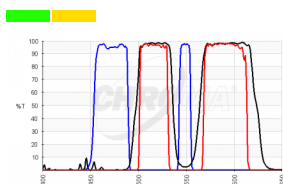**59004 - ET - FITC/TRITC**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

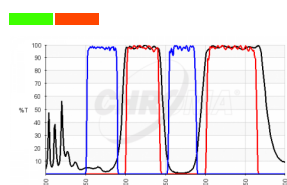**59022 - ET - EGFP/mCherry (or FITC/TxRed)**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

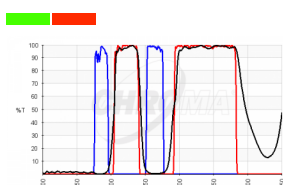**59009 - ET - FITC/CY3**


This dual band set is intended to provide an attenuated Cy3/TRITC/AF546 signal in applications where these orange channel signals overwhelm the green channel signals, and compensation is required.

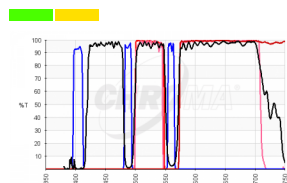**59908 - ET - 488/561 PALM Dual Band Emission Set with 405nm activation**

For use with imaging lasers between 486-492nm and 558-563nm and photoactivation laser 400-410nm

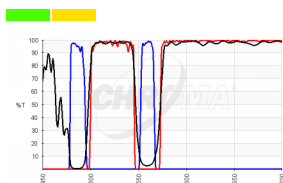
**59904 - ET - 488/561nm Laser Dual Band Set**

For use with lasers between 486-492nm and 556-563nm

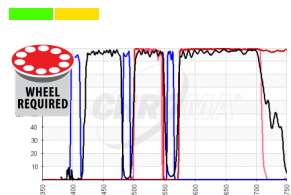
TRF59908-EM - ET - PALM Dual Band Emission Set for 488/561nm TIRF applications with 405nm excitation

For use with imaging lasers between 486-492nm and 558-563nm and photoactivation laser 400-410nm
Set includes the same filter combination as the 59908 PALM set, but with the 2mm TIRF dichroic.

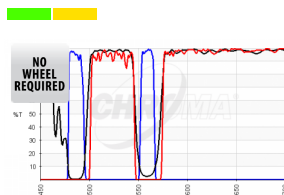
TRF59904 - ET - 488/561nm Laser Dual Band Set for TIRF applications

For use with lasers between 486-492nm and 556-563nm
Set includes no additional emission filters; instead, the dual band emission filter has increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

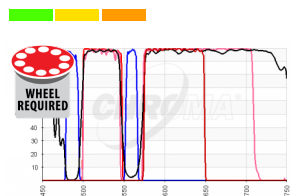
TRF59904-EM - ET - 488/561nm Laser Dual Band Set for TIRF applications

For use with lasers between 486-492nm and 556-563nm
Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

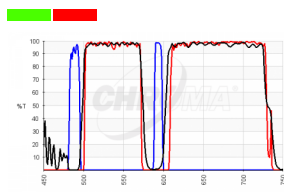
59905 - ET - 488/594nm Laser Dual Band Set

For use with lasers between 486-492nm and 591-596nm

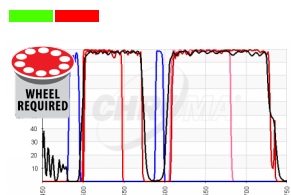
TRF59905-EM - ET - 488/594nm Laser Dual Band Set for TIRF applications

For use with lasers between 486-492nm and 591-596nm
Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

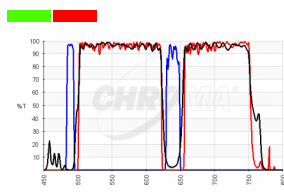
59906 - ET - 488/640nm Laser Dual Band Set

For use with lasers between 486-492nm and 633-647nm

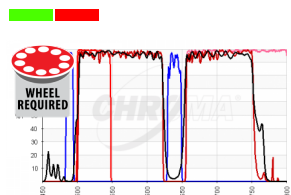
TRF59906-EM - ET - 488/640nm Laser Dual Band Set for TIRF applications

For use with lasers between 486-492nm and 633-647nm
Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

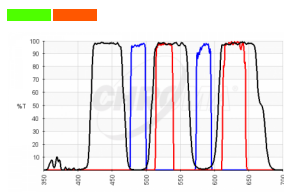
59010 - ET - Green/Red #1 FISH

This dual band filter set is ideal for simultaneous viewing of commercially available Green and Red FISH probes. Intended to balance the two fluorescence signals and also avoid detection of Gold probes if also present.

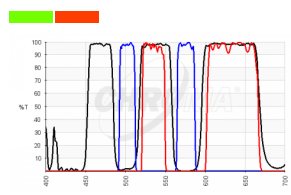
59026 - ET - EYFP/mCherry

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

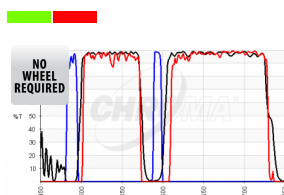
TRF59905 - ET - 488/594nm Laser Dual Band Set for TIRF applications

For use with lasers between 486-492nm and 591-596nm
Set includes no additional emission filters; instead, the dual band emission filter has increased blocking for greater attenuation of TIRF lasers without the use of additional filters.

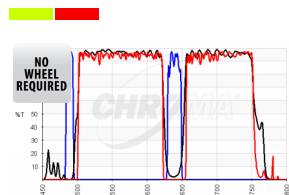
**TRF59906 - ET - 488/640nm Laser
Dual Band Set for TIRF
applications**

For use with lasers between 486-492nm and 633-647nm
Set includes no additional emission filters; instead, the dual band
emission filter has increased blocking for greater attenuation of TIRF
lasers without the use of additional filters.

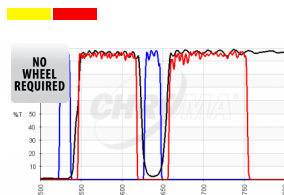**TRF59907 - ET - 532/640nm Laser
Dual Band Set for TIRF
applications**

For use with lasers between 528-535nm and 633-647nm
Set includes no additional emission filters; instead, the dual band
emission filter has increased blocking for greater attenuation of TIRF
lasers without the use of additional filters.

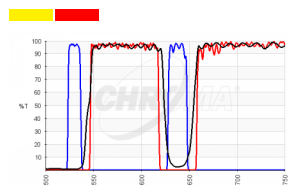**59907 - ET - 532/640nm Laser Dual
Band Set**

For use with lasers between 528-535nm and 633-647nm


**TRF59907-EM - ET - 532/640nm
Laser Dual Band Set for TIRF
applications**

For use with lasers between 528-535nm and 633-647nm
Set includes additional emission filters for use with external
emission filter wheel for signal selection and greatest attenuation of
TIRF lasers.

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
Single Band Exciters/Multi Band Emitter

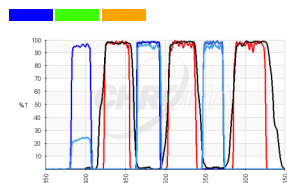
**89084v2 - ET - DAPI/Green/Orange
Multi-Exciter FISH**

This filter set replaces the original 82000 set and the recent 89084 set. It is a tripleband filter set with passbands equivalent to the 69013v2 set. In addition, This set also contains a Green/Orange dualband excitation filter and 3 individual single band excitaiton filters. This allows the user to detect DAPI/Green/Orange simultaneously with the tripleband exciter, or only Green/Orange with the dualband exciter, or single colors for imaging to a camera. The tripleband excitation filter features an attenuated level of DAPI excitation.

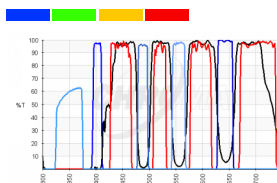
**89400 - ET - DAPI/FITC/TRITC/Cy5
Quad with Single Band Exciters**

A version of the 89000 Quad set with a quad band emission filter and four individual excitation filters for use in an external filter housing such as a filter wheel, similar to Chroma's original Pinkel filter set. Allows for selective excitation of individual fluors in sample containing multiple fluors. Appropriate for live viewing or image acquisition with monochrome or color cameras.

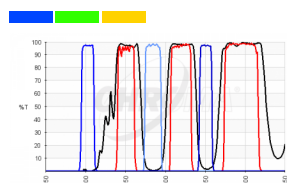
**69300 - ET - DAPI/FITC/TRITC with
single band exciters**

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

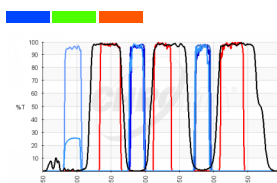
**89085 - ET - DAPI/Green/Red Multi-
Exciter FISH**

This filter set replaces the original 83000 set and the recent 89083 set. Tripleband filter set with passbands equivalent to the 69015 set. In addition, This set also contains a Green/Red dualband excitation filter and 3 individual single band excitaiton filters. This allows the user to detect DAPI/Green/Red simultaneously with the tripleband exciter, or only Green/Red with the dualband exciter, or single colors for imaging to a camera.

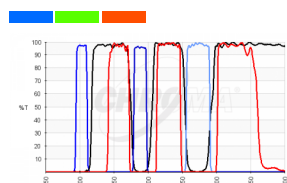
**69302 - ET - DAPI/FITC/Texas Red®
with single band exciters**

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

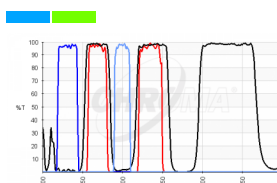
**59217 - ET - ECFP/EYFP with single
band exciters**

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

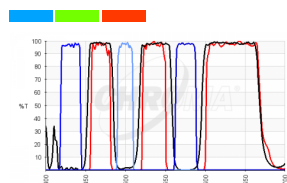
**69308 - ET - ECFP/EYFP/mCherry
with single band exciters**

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

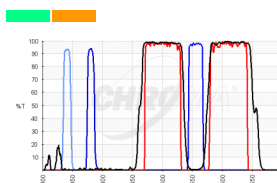
79003 - ET - Fura 2/TRITC

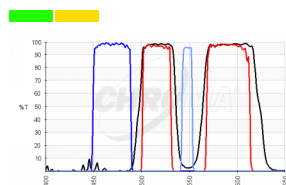
59204 - ET - FITC/TRITC with single band exciters

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

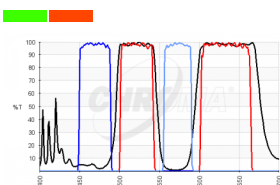59222 - ET - EGFP/mCherry (or FITC/TxRed) with single band exciters


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

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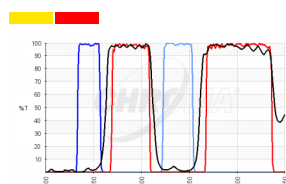
Multi-Band Exciters/Single Band Emitters

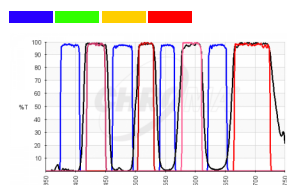
89404 - ET -
391-32/479-33/554-24/638-31 Multi
LED set with single band emitters

Single emission filter version of 89402 set.
 Quad band catalog set for Multi LED light engines like those manufactured by Lumencor, CoolLED, Excelitas and Cairn Research. Excitation passbands are matched to discrete, triggerable LEDs appropriate for excitation of DAPI/FITC/TRITC/Cy5 and similar fluorors. Not intended for use with white light sources.

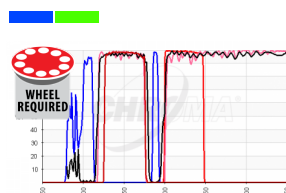
TRF59901-EM - ET - 405/488nm
Laser Dual Band Set for TIRF
applications

For use with lasers between 400-408nm and 485-491nm
 Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

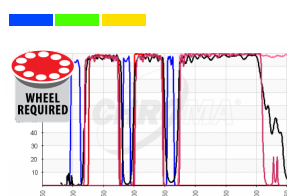
TRF69901-EM - ET -
405/488/561nm Laser Triple Band
Set for TIRF applications

For use with lasers between 400-410nm, 486-492nm and 558-563nm
 Set includes additional emission filters for use with external emission filter wheel for signal selection and greater attenuation of TIRF lasers.

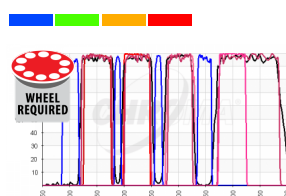
TRF89902-EM - ET -
405/488/561/647nm Laser Quad
Band Set for TIRF applications

Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.

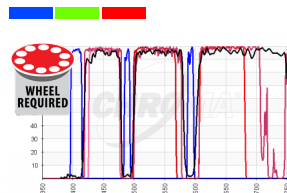
TRF69902-EM - ET -
405/488/594nm Laser Triple Band
Set for TIRF applications

For use with lasers between 400-412nm, 486-492nm and 591-596nm
 Set includes additional emission filters for use with external emission filter wheel for signal selection and greater attenuation of TIRF lasers.

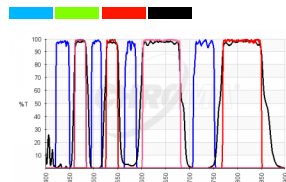
89405 - ET -
436-28/506-21/578-24/730-40 Multi
LED set with single band emitters

Single emission filter version of 89403 set.
 Quad band catalog set for Multi LED light engines like those manufactured by Lumencor, CoolLED, Excelitas and Cairn Research. Excitation passbands are matched to discrete, triggerable LEDs appropriate for excitation of CFP/YFP/RFP/Cy7 and similar fluorors. Not intended for use with white light sources.

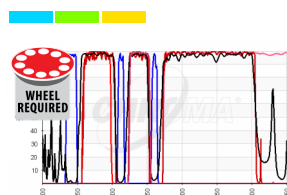TRF69904-EM - ET -
442/514/561nm Laser Triple Band
Set for TIRF applications

For use with lasers between 437-448nm, 511-516nm and 559-564nm
 Set includes additional emission filters for use with external emission filter wheel for signal selection and greater attenuation of TIRF lasers.

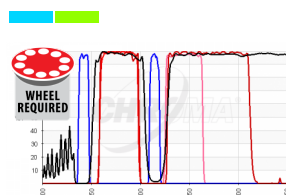
TRF59902-EM - ET - 442/514nm
Laser Dual Band Set for TIRF
applications

For use with lasers between 439-445nm and 512-517nm
 Set includes additional emission filters for use with external emission filter wheel for signal selection and greatest attenuation of TIRF lasers.


**TRF69905-EM - ET -
442/514/594nm Laser Triple Band
Set for TIRF applications**

For use with lasers between 400-449nm, 511-518nm and 588-597 nm
Set includes additional emission filters for use with external
emission filter wheel for signal selection and greater attenuation of
TIRF lasers.

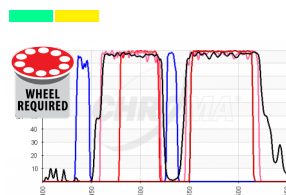
**TRF59903-EM - ET - 442/532nm
Laser Dual Band Set for TIRF
applications**

For use with lasers between 438-445nm and 530-535nm
Set includes additional emission filters for use with external
emission filter wheel for signal selection and greatest attenuation of
TIRF lasers.

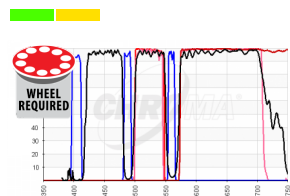**TRF59908-EM - ET - PALM Dual
Band Emission Set for 488/561nm
TIRF applications with 405nm
excitation**

For use with imaging lasers between 486-492nm and 558-563nm and
photoactivation laser 400-410nm
Set includes the same filter combination as the 59908 PALM set, but
with the 2nm TIRF dichroic.

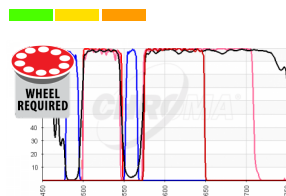
**TRF59904-EM - ET - 488/561nm
Laser Dual Band Set for TIRF
applications**

For use with lasers between 486-492nm and 556-563nm
Set includes additional emission filters for use with external
emission filter wheel for signal selection and greatest attenuation of
TIRF lasers.


**TRF59905-EM - ET - 488/594nm
Laser Dual Band Set for TIRF
applications**

For use with lasers between 486-492nm and 591-596nm
Set includes additional emission filters for use with external
emission filter wheel for signal selection and greatest attenuation of
TIRF lasers.

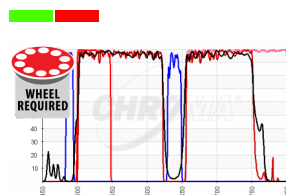**TRF59906-EM - ET - 488/640nm
Laser Dual Band Set for TIRF
applications**

For use with lasers between 486-492nm and 633-647nm
Set includes additional emission filters for use with external
emission filter wheel for signal selection and greatest attenuation of
TIRF lasers.

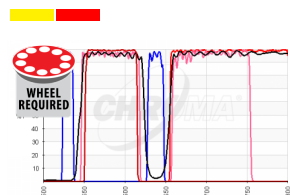
**TRF59907-EM - ET - 532/640nm
Laser Dual Band Set for TIRF
applications**


For use with lasers between 528-535nm and 633-647nm
Set includes additional emission filters for use with external
emission filter wheel for signal selection and greatest attenuation of
TIRF lasers.

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Single Band Exciters/Single Band Emitters

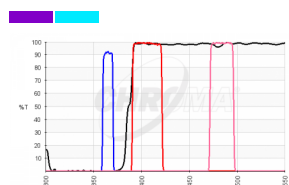
79006 - ET - INDO

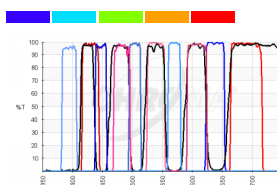
89903 - ET -
BV421/BV480/AF488/AF568/AF647
Quinta Band set

For use with BD Horizon Brilliant™ Violet 421 and Brilliant™ Violet 480 in 5-color fluorescence imaging scheme. Separate filter housings (such as external filter wheels) required for individual excitation and emission filters.

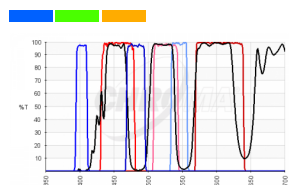
89013 - ET - DAPI/FITC/TRITC

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

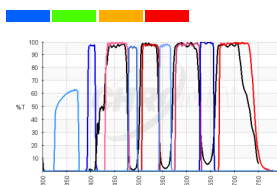
89000 - ET - Sedat Quad

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

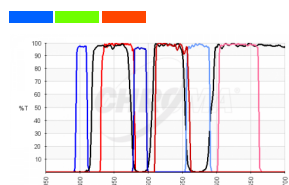
89014 - ET - DAPI/FITC/Texas Red®

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

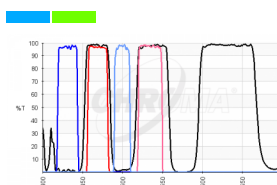
89002 - ET - ECFP/EYFP

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

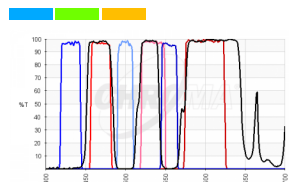
89007 - ET - ECFP/EYFP/DsRed

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

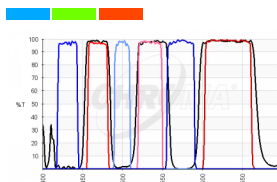
89006 - ET - ECFP/EYFP/mCherry

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

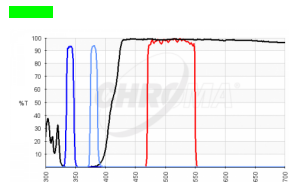
79001 - ET - Fura 2

Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

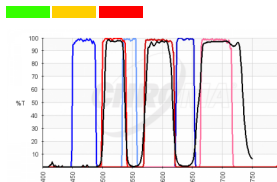
89016 - ET - FITC/Cy3/Cy5



Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

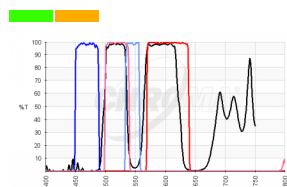**89019 - ET - FITC/TRITC**




Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

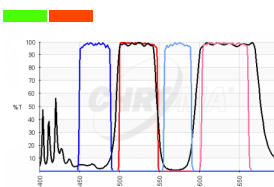**89021 - ET - EGFP/mCherry**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

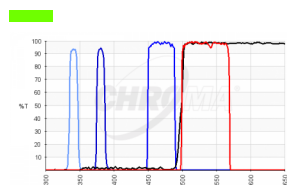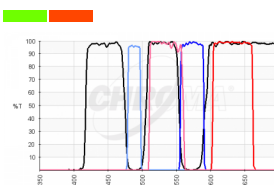**79004 - ET - Fura 2/GFP****89017 - ET - FITC/Texas Red®**


Please see the excitation intensity note on page 6 for information on how to avoid rapid photobleaching of sample.

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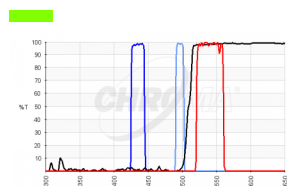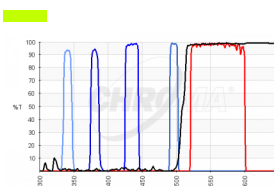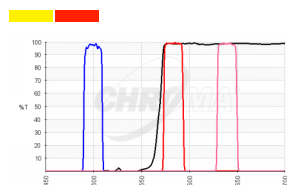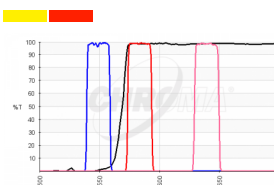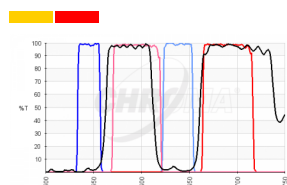Cubes, Sliders and Rings

**Product Number 91000 - 18mm
Nikon Microscopes (Diaphot, TMD)**

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**Product Number 91003 - Olympus
IMT2 Filter Cube**

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**Product Number 91005 - Leica
Aristoplan Filter Cube**

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**Product Number 91006 - Zeiss
Axioline 3FL Slider**

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**Product Number 91007 - Leica DM
Filter Cube**

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**Product Number91009 - Leica
DMIRB Filter Cube**

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**Product Number91011 - Zeiss
Axioline 4FL Slider**

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**Product Number91012 - Zeiss
Axioline 5FL Slider**

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**Product Number91013 - Zeiss
Axioline 6FL Slider**

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**Product Number91015 - Zeiss Axio
screws into turret**

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**Product Number91016 - Zeiss
Axioplan with non-tilted emission
filter (Push & Click)**

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**Product Number91018 - BX2 filter
cube (U-MF2) For Olympus BX2
and IX2 models**

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Product Number91026 - Nikon prism-to-dichroic conversion mount

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Product Number91027 - Zeiss Axio inverted 2F

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Product Number91029 - Zeiss Axio Imager with tilted emission filter (Push & Click)

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Product Number91032 - Laser TIRF for Nikon TE2000/Ti

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Product Number91033 - Nikon Excitation Filter Ring 18mm (for 91000 cube)

Product Number 91034 - Nikon Emission Filter Ring 18mm (for 91000 cube)

Product Number 91035 - Nikon Excitation Filter Ring 25mm (for 91001 cube)

Product Number 91036 - Nikon Emission Filter Ring 25mm (for 91001 cube)

Product Number 91037 - Zeiss Axio Filter Retaining Ring 25mm (for all cubes and sliders)

Product Number 91038 - BX3/IX3 filter cube (U-FF, U-M710) For Olympus BX3 and IX3 models for 25mm filters

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Product Number 91040 - Olympus U-M1000 for CX41

Exciter (x) = 20mm, Emitter (m) = 20mm, Beamsplitter = 21x29x1mm.

Product Number 91041 - Laser TIRF for Olympus BX2 models

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Product Number 91042 - Laser TIRF for Zeiss Axio

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Product Number 91043 - IX3-FFXL filter cube (U-M101) for Olympus IX3 models for 32mm filters

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Product Number 91044 - Laser TIRF for Olympus BX3/IX3 models, for 25mm filters

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Product Number 91046 - Nikon P2-EFLC filter cube for SMZ18 & SMZ25 Stereo Microscope with 25mm filters

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Product Number 91049 - Nikon P-EFLC filter cube for SMZ800n & SMZ1270 Stereo Microscope with 18mm filters

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Product Number 91050 - Leica DMi8 standard filter cube, for 25mm filters, fits 1mm-thick dichroics

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Product Number 91051 - Leica DMi8 filter cube, for 25mm filters, fits dichroics up to 3mm thick

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Product Number 91054 - Nikon C-FLL-C Large FOV cube

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Product Number 91056 - Keyence filter cube for BZX series of fluorescence microscopes

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Diagnostic Slides

CHROMA TECHNOLOGY CORP

An Employee-owned Company Producing the World's Finest Optical Filters and Filter Sets.

US Corporate Headquarters
10 Imtec Lane
Bellows Falls, Vermont 05101 USA
sales@chroma.com
Tel: +1-800-824-7662
+1-802-428-2500

Chroma Technology GmbH
Maximilianstrasse 33
D-82140 Olching, Germany
europe@chroma.com
Tel: +49-8142-2847525

Chroma China
Rm 803, Building No.16
Yuanbo Wuli,
Xingjin Road, Jimei district
Xiamen City, Fujian Province,
361021 China
china@cn.chroma.com
Tel: +86-0592-5062089

Chroma Technology Japan
8F Yokohama Onoecho Building
4-57 Onoecho Naka-ku
Yokohama 231-0015 Japan
japan@chroma.com
Tel: +81 (0) 45 285 1583